

Welcome to the SYNAPT G2-Si HDMS Online Information System



What do you want to do?

- Find out what's new
- Browse or search the Help system



Do you need additional help?

- Contact Waters Technical Support
- Search our knowledge base
- Visit a Waters User Forum
- Access additional product documentation

Copyright

Copyright © Waters Corporation 2011 – 2015 All Rights Reserved.

This document contains information that is proprietary to Waters Corporation. The original recipient of this document may duplicate this document in whole or in part for internal business purposes only, provided that this entire notice appears in all copies. In duplicating any part of this document, the recipient agrees to make every reasonable effort to prevent the unauthorized use and distribution of the proprietary information.

The information in this Help system is subject to change without notice and should not be construed as a commitment by Waters Corporation. Waters Corporation assumes no responsibility for any errors that may appear in this Help system. This Help system is believed to be complete and accurate at the time of publication. In no event shall Waters Corporation be liable for incidental or consequential damages in connection with, or arising from, the use of this document.

ACQUITY, ACQUITY UPLC, UPLC, Xevo, and Waters are registered trademarks of Waters Corporation. UNIFI, Empower, MassLynx, and "THE SCIENCE OF WHAT'S POSSIBLE." are trademarks of Waters Corporation.

Microsoft, Internet Explorer, Windows, and Windows 7 are registered trademarks of Microsoft Corporation in the United States and/or other countries.

Oracle is a registered trademark of Oracle Corporation and/or its affiliates.

Other trademarks or registered trademarks are the sole property of their respective owners.

Related Topics

- Waters contact information

Safety Advice

Waters instruments display hazard symbols designed to alert you to the hidden dangers of operating and maintaining the instruments. Their corresponding user guides include those symbols, accompanying them with text statements describing the hazards and telling you how to avoid them. The safety symbols and statements that apply to the entire line of Waters products appears below.

Warning symbols

Warning symbols alert you to the risk of death, injury, or seriously adverse physiological reactions associated with an instrument's use or misuse. Heed all warnings when you install, repair, and operate Waters instruments. Waters assumes no liability for the failure of those who install, repair, or operate its instruments to comply with any safety precaution.

Task-specific hazard warnings

The following warning symbols alert you to risks that can arise when you operate or maintain an instrument or instrument component. Such risks include burn injuries, electric shocks, ultraviolet radiation exposures, and others.

The text accompanying the following symbols appears in parenthesis because when the following symbols appear in a manual's narratives or procedures, their accompanying text identifies the specific risk and explains how to avoid it.



Warning: (General risk of danger. When this symbol appears on an instrument, consult the instrument's user documentation for important safety-related information before you use the instrument)



Warning: (Risk of burn injury from contacting hot surfaces)



Warning: (Risk of electric shock)



Warning: (Risk of fire)



Warning: (Risk of needle puncture)



Warning: (Risk of injury caused by moving machinery)



Warning: (Risk of exposure to ultraviolet radiation)



Warning: (Risk of contacting corrosive substances)



Warning: (Risk of personal contamination with a toxic substance)



Warning: (Risk of personal exposure to laser radiation)



Warning: (Risk of exposure to biological agents that can pose a serious health threat)



Warning: (Risk of crushing. Tipping hazard. Maximum weight as indicated.)

Warnings that apply to particular instruments, instrument components, and sample types

The following warnings can appear in the user manuals of particular instruments and on labels affixed to them or their component parts.

Burst warning

This warning applies to Waters instruments fitted with nonmetallic tubing.

Warning: Pressurized nonmetallic, or polymer, tubing can burst. Observe these precautions when working around such tubing:



- Wear eye protection.
- Extinguish all nearby flames.
- Do not use tubing that is, or has been, stressed or kinked.
- Do not expose nonmetallic tubing to incompatible compounds like tetrahydrofuran (THF) and nitric or sulfuric acids.
- Be aware that some compounds, like methylene chloride and dimethyl sulfoxide, can cause nonmetallic tubing to swell, which significantly reduces the pressure at which the tubing can rupture.

Mass spectrometer flammable solvents warning

This warning applies to instruments operated with flammable solvents.



Warning: Where significant quantities of flammable solvents are involved, a continuous flow of nitrogen into the ion source is required to prevent possible ignition in that enclosed space.

Ensure that the nitrogen supply pressure never falls below 400 kPa (4 bar, 58 psi) during an analysis in which flammable solvents are used. Also ensure a gas-fail connection is connected to the HPLC system so that the LC solvent flow stops if the nitrogen supply fails.

Mass spectrometer shock hazard

This warning applies to certain instruments when they are in Operate mode.

Biohazard warning

This warning applies to Waters instruments that can be used to process material that might contain biohazards: substances that contain biological agents capable of producing harmful effects in humans.



Warning: Waters instruments and software can be used to analyze or process potentially infectious human-sourced products, inactivated microorganisms, and other biological materials. To avoid infection with these agents, assume that all biological fluids are infectious, observe Good Laboratory Practices and, consult your organization's biohazard safety representative regarding their proper use and handling. Specific precautions appear in the latest edition of the US National Institutes of Health (NIH) publication, *Biosafety in Microbiological and Biomedical Laboratories*.

Chemical hazard warning

This warning applies to Waters instruments that can process corrosive, toxic, flammable, or other types of hazardous material.



Warning: Waters instruments can be used to analyze or process potentially hazardous substances. To avoid injury with any of these materials, familiarize yourself with the materials and their hazards, observe Good Laboratory Practices (GLP), and consult your organization's safety representative regarding proper use and handling. Guidelines are provided in the latest edition of the National Research Council's publication, *Prudent Practices in the Laboratory: Handling and Disposal of Chemicals*.

Warnings that apply to all Waters instruments



Warning: To avoid electric shock, when the safety protection of an instrument has been compromised in any way (which would include damage sustained to its protective covers), disconnect the instrument from all power sources, and secure it against unintended operation.



Warning: To avoid electric shock, disconnect the instrument from all power sources before servicing it.



Warning: To protect against fire, replace open or otherwise defective fuses with fuses of the same type and rating.



Warning: If an instrument or component is operated in a manner not specified by its manufacturer, safety protection can be impaired.



Warning: To avoid electric shock, when the safety protection of an instrument has been compromised in any way (which would include damage sustained to its protective covers), disconnect the instrument from all power sources, and secure it against unintended operation.



Warning: To avoid electric shock, disconnect the instrument from all power sources before servicing it.



Warning: To protect against fire, replace open or otherwise defective fuses with fuses of the same type and rating.


Caution symbol

The caution symbol signifies that an instrument's use or misuse can damage the instrument or compromise a sample's integrity. The following symbol and its associated statement are typical of the kind that alert you to the risk of damaging the instrument or sample.



Notice: To prevent static charges and resultant circuit damage, do not remove the covers that protect integrated circuit chips.

Introducing the (Undefined variable: MyVariables.Product_Name_i)

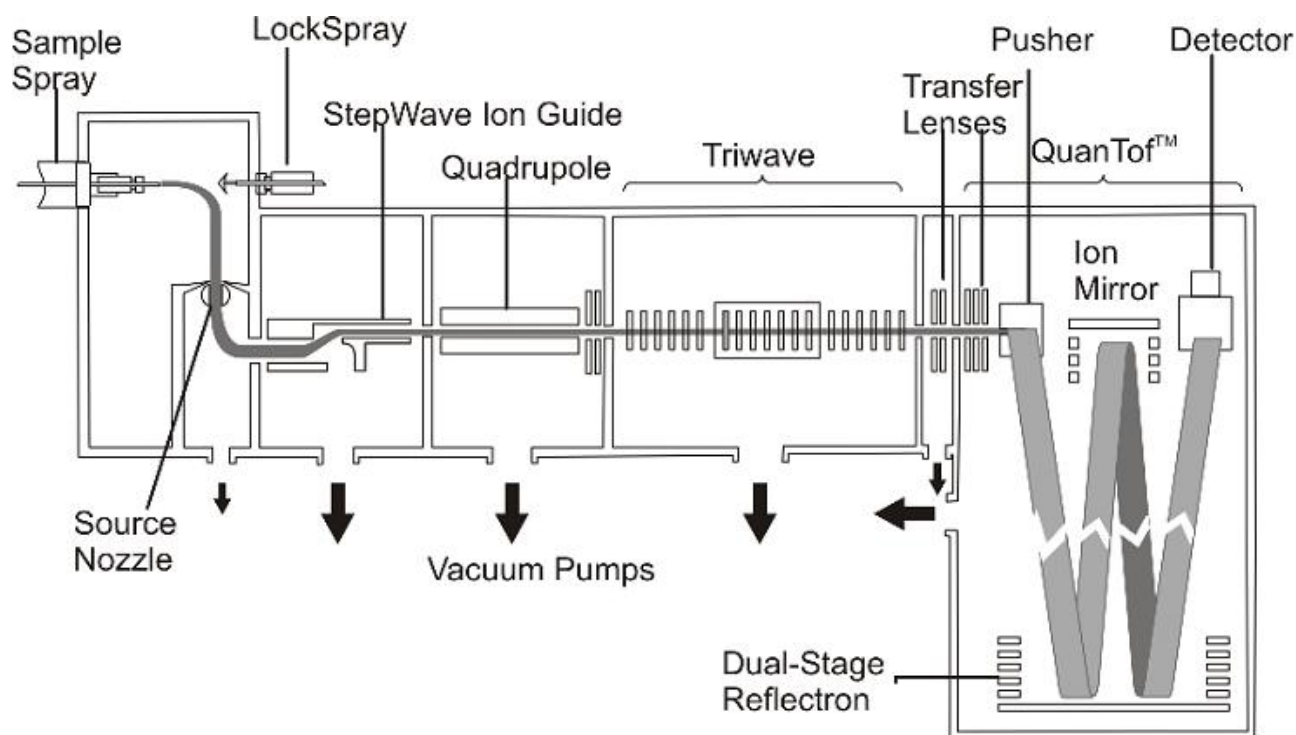

The (Undefined variable: MyVariables.Product_Name_i) is a hybrid quadrupole, time-of-flight (TOF), orthogonal acceleration, mass spectrometer that delivers routine high-resolution, exact mass measurement over a wide range of experimental options. SYNAPT G2-Si provides an operating mode in which ion mobility separations are either enabled (HDMS mode) or disabled (Tof mode).

The source regions consists of a dual-orthogonal 'Z-spray' ion extraction geometry and off-axis 'Stepwave' ion transfer optics to maximize the efficiency of ion extraction and transfer from the ionisation source to the quadrupole and Triwave ion optics.

Triwave is the enabling technology within the SYNAPT G2-Si HDMS System. Three T-Wave ion guides are employed, allowing ions to be trapped, separated based on their mobility, and transferred to the Quantitative Tof (QuanTof™) analyzer for high-resolution mass analysis. The innovative configuration of the Triwave ensures the introduction of IMS is not made at the expense of sensitivity. Enhanced ion mobility resolution has been achieved through increased pressure (with the use of a novel Helium filled entry cell in the IMS T-Wave), increased length of the IMS T-Wave and increased frequency/amplitude of the T-wave operation. The TRAP and TRANSFER T-Wave regions can also be used as collision cells, and in combination with ion mobility separations, can provide a unique route to more comprehensive structural characterization.

QuanTof is the enabling next- generation TOF analyzer of SYNAPT G2-Si. QuanTof's high field pusher and dual-stage reflectron, incorporating high transmission parallel wire grids, reduce ion turnaround times due to pre-push kinetic energy spread and improve focusing of high-energy ions, respectively. High resolution and exact mass measurement is delivered at up to 30 spectra/sec in Tof mode and up to 10 spectra/sec in HDMS mode. The ion detection system combines an ultrafast electron multiplier and hybrid analog-to-digital conversion (ADC) detector electronics to provide high resolution and low noise at low ion currents, high quantitative performance and exact mass at high ion currents, all at the elevated data acquisition rates of ion mobility analyses (the integrated IMS-ToF acquisition delivers the equivalent of 10 x 200 spectra/sec). The hybrid ADC system comprises an 8-bit ADC sampling at 3 GHz feeding to a field-programmable gate array (FPGA) for signal processing, and subsequently to a block of memory for accumulating the 200 sequential mass spectra that form the mobility arrival time spectrum.

SYNAPT G2-Si incorporates new or enhanced data acquisition techniques and processing workflows for untargeted or targeted applications - from discovery and quantitation, to screening, imaging and compound characterisation.

Ionization options

The instrument's sources, except for NanoFlow and ionKey sources, are designed for the high flow rates that UPLC use requires. The following table describes source options. Note that multi-mode ionization in a single run is possible using UNIFI software. The NanoFlow and ionKey sources are designed for the very low flow rates that a nanoACQUITY UPLC or ACQUITY UPLC M-Class system requires.

You can interchange sources without using tools. When you do so, the software automatically detects the new source.

Source options:	Source Description
Electrospray (ESI)	Ideal for ionic, polar, thermo labile compounds, or large biomolecules (m/z greater than 1000) such as proteins, peptides and oligonucleotides. Atmospheric pressure chemical ionization (APCI) produces singly charged protonated or deprotonated molecules for a broad range of nonvolatile samples.
nanoLockSpray	Allows electrospray ionization performed in the flow rate range of 5 to 1000 nL/min. For a given sample concentration, the ion currents for similar experiments approximate to those in normal flow rate electrospray. However, because sample consumption is greatly reduced, the sensitivity gains are significant when you adopt similar scan parameters
Combined electrospray and atmospheric pressure chemical ionization (ESCI)	The standard ESI probe is used in conjunction with a corona pin to allow alternating acquisition of ESI and APCI ionization data, facilitating high-throughput processing and wider compound coverage.
Atmospheric pressure chemical ionization (APCI)	For applications similar to those that perform with ESCI. APCI uses an optional, dedicated APCI probe for efficient droplet evaporation.
Atmospheric solids analysis probe (ASAP)	Facilitates rapid analysis of volatile and semivolatile compounds in solids, liquids, and polymers. It is particularly suited to analyzing low-polarity compounds. The ASAP directly replaces the electrospray or IonSABRE II probe in the instrument's source housing and has no external gas or electrical connections.
Atmospheric Pressure Photo Ionization (APPI)	Uses photons generated by a discharge UV lamp (~10.2 eV) to produce sample ions from vaporized LC eluent. Direct photoionization of the sample molecule occurs when the photon energy exceeds the ionization potential of the sample molecule. The optional dual-mode (APPI/APCI) ionization source incorporates an APPI source enclosure used in conjunction with a standard APCI probe. You can operate the source in APPI, APCI, or dual-mode, which switches rapidly between ionization modes, facilitating high-throughput analyses
Matrix-Assisted Laser Desorption Ionization (MALDI)	The optional (MALDI) interface enables rapid, tool-free switching between API and MALDI modes. A motorized stage moves the MALDI source into position.
ionKey M/S	The Waters ionKey/MS system integrates micro-scale UPLC separation into the source of the mass spectrometer, combining the ACQUITY UPLC M-Class and the SYNAPT G2-Si. Powered by the iKey, this technology simplifies the implementation of micro-scale chromatography and analysis of limited-volume samples, with significant increases in sensitivity compared to 2.1 mm ID chromatography.
Desorption Electro-Spray Ionisation (DESI)	The DESI ion source gives the full advantage of desorption electrospray ionisation technology. Equipped with a mechanism for automated sample positioning, it is designed for use with SYNAPT G2-Si mass spectrometers. The combination of the SYNAPT G2-Si and DESI ion source eliminates the need for complex and time consuming sample preparation.

Source options:	Source Description
Rapid Evaporative Ionisation Mass Spectrometry (REIMS)	<p>REIMS is a technique that allows for near real-time characterization of tissue. REIMS can be used to analyze the aerosol or "smoke" that is released during an electro-surgical dissection procedure.</p> <p>REIMS when coupled with electro-surgery for the diagnosis of tissue samples is known as intelligent knife (iKnife).</p>

What's new in this version

Enhancements introduced with this latest version of the SYNAPT G2-Si mainly concern introduction of a new ion source, thus extending the instrument's capabilities as follows:

Desorption Electrospray Ionisation (DESI) ion source

Desorption Electrospray ionization (DESI) is a pneumatically assisted electrospray technique in which a solvent mixture is directed toward a sample of interest and the resultant ions sampled directly, under ambient conditions, in the SYNAPT G2-Si.

The efficacy of DESI, as a direct-desorption technique, has been demonstrated in a wide range of applications, from explosives to proteomics.

Rapid Evaporative Ionisation Mass Spectrometry (REIMS) source

REIMS is a technique that allows for near real-time characterization of tissue. REIMS can be used to analyze the aerosol or "smoke" that is released during an electrosurgical dissection procedure.

REIMS when coupled with electro-surgery for the diagnosis of tissue samples is known as intelligent knife (iKnife).

Related Topics

- [Introduction](#)

Waters Corporation contact information

Technical support

If you have enhancement requests or technical questions regarding the use of any Waters Corporation product, contact us as follows:

Internet	<p>The Waters Web site includes contact information for Waters locations worldwide.</p> <p>Visit www.waters.com</p>
Telephone and fax	<p>From the USA or Canada, phone 800 252-HPLC, or fax 508 872 1990.</p> <p>For other locations worldwide, phone and fax numbers appear in the Waters Web site.</p>
Conventional mail	<p>Waters Corporation</p> <p>34 Maple Street</p> <p>Milford, MA 01757</p> <p>USA</p>

Your comments

We welcome your feedback on this online information system documentation or any other Waters Corporation product or document. We are also interested in your suggestions for any additional topics you want included in our product documentation. Please contact us at the following email address with your comments:

Tech_Comm@waters.com

Related Topics

- How to use this Online Information System
- Product documentation overview

Access spectrum tools

- Windows Desktop > Start > All Programs > MassLynx > Spectrum Tools

Two tools are provided with MassLynx but they are accessed separately as a Windows PC tool, not from the MassLynx application .

Mass Difference Calculator

To compare peak information obtained experimentally with standard reference peaks, use MassLynx's Mass Difference Calculator. The calculator calculates the mass difference between each reference peak and the nearest entry in the peak list, reporting the difference as the absolute error (mDa) or fractional error (ppm).

To start Mass Difference Calculator:

At the workstation Windows desktop, click Start > All Programs > MassLynx > Spectrum Tools > MassDiff.

Result: The calculator opens with a set of empty data grids.

To exit the calculator:

Click File > Exit.

ResCalc (Resolution Calculator)

You can examine spectral data in detail using ResCalc, the resolution calculator, which determines and reports the spectral resolution from spectra produced with MassLynx software.

ResCalc produces a report showing the resolution of the most intense peak displayed in a spectrum at FWHH (full-width-half-height) intensity.

To start ResCalc:

Click Start > All Programs > MassLynx > Spectrum Tools > ResCalc.

Result: ResCalc opens with an empty screen.

To exit from ResCalc:

Click File > Exit.

Before you begin using the instrument

Waters designed the SYNAPT G2-Si instrument to be used as a research tool to deliver authenticated, exact-mass measurement. It is not for use in diagnostic procedures.

Requirements: The SYNAPT G2-Si has been configured for your application by the Waters installer.

To prepare for data acquisition:

Before starting the instrument, you must complete the following tasks:

- Ensure that the instrument is pumped down (evacuated) and the ionization source for your application fitted. See the *Waters SYNAPT G2-Si Mass Spectrometry System Operator's Overview and Maintenance Guide* (part number 715004102) on the SYNAPT G2-Si Documentation CD (part number 715004101).
- Visually inspect the ion probe position.
- Ensure that the wash bottle contains a sufficient volume of a suitable wash solution.

Tip: The composition of the wash solution is sample-dependent. However, for most applications, Waters recommends an 80:20 methanol/water mixture.

- Ensure that your sample, lock-mass reference, and calibrant solutions are ready and mounted on the instrument.

Recommendation: Use reservoir bottles in this order: A for sample, B for Lockmass reference, C for calibrant.

Tip: Filter or centrifuge the sample solution to remove particulates before placing it in reservoir A, to avoid blocking the transfer line.

- Ensure the waste line is correctly plumbed.
- Ensure the LC system is configured and ready to use.
- If your SYNAPT G2-Si is the MALDI version, ensure the front support legs are correctly adjusted. See the *Waters MALDI SYNAPT G2-Si Operator's Overview and Maintenance Guide*.

Related Topics

- [Start the instrument](#)
- [Shut down the instrument](#)
- [Prepare for use](#)

Control fluidics and external devices

- [MassLynx > MS Method > Method Events](#)

Events can provide additional control over the onboard sample delivery (fluidics) system during a method. Events repeat for each sample and require a finite amount of time for completion.

Tip: To clarify the operation of an event, find it on the IntelliStart Fluidics console. With the exception of Event Out and Solvent Delay, the method events mirror the controls on the IntelliStart Fluidics console.

Requirement: To start using method events, select the "Enable method events" checkbox.

Events and their actions:



Event	Actions	Description
Event out 1	On Off Pulse Toggle	Controls the Event out 1 relay connection.
Event out 2	On Off Pulse Toggle	Controls the Event out 2 relay connection.
Infusion	Start Stop	Controls syringe infusion from on-board solvent delivery system
Injection	Injection	Selects inlet injection as sample source

Event	Actions	Description
Flow state	LC Combined Infusion Waste	Selects the position of the diverter valve.
Flow rate	Numeric	The desired rate in $\mu\text{L}/\text{min}$
Refill	Refill Always Refill Auto	The Refill Always function refills at the start of every sample run The Refill Auto function refills when the syringe empties. Refill Auto is the recommended setting.
Reservoir	A B C Wash	The instrument reservoir selected for the given system.
Fill Volume	Numeric	The syringe fill volume, in increments of 50 μL .
Solvent Delay	Time	A timed pause starting from the beginning of the acquisition to allow for sample retention times.

Example: Diverting solvent delay to waste

Diverting solvent delay to waste can reduce contamination from lighter molecules that can degrade mass spectrometer performance: for example, sodium salts in dirty analytes such as blood plasma, can circulate in the source and adversely affect results.

To divert solvent delay to waste:

- Click MS Method  , click Method Events  .
- Select the Enable checkbox.
- Add events to the events table, and then specify values for Time, Event, Action and System.

Example event selections:

Option	Description
Time	Select your expected solvent delay, in minutes
Event	Select Solvent Delay.
Action	Select Waste.
System	Select Sample.

- Click OK.
- To prevent unwanted heating while the solvent flow is directed waste, in the "Solvent Delay Options" box, set the API probe temperature to 20°C.
- Click OK.

Related Topics

- Prepare instrument for use
- Define Reference Compounds

Monitor and control solvent delivery

The IntelliStart Fluidics console consists of two areas described here:

- Monitoring and controlling sample flows
- Monitoring and controlling LockSpray flow

The buttons and selections described in this topic are also available from the Tune Window.




Warning: To avoid cuts from broken glass, injuries from falling objects, and exposure to toxic substances, never place solvent containers or other objects on the front covers, or on top of the SYNAPT G2-Si.

Monitoring and controlling sample flows

The IntelliStart Fluidics, solvent-delivery system consists of a multi-position selector valve with input connections from an infusion syringe and LC column. The system's output flows to the mass spectrometer's source (probe) and a waste line. You can monitor this system in real time on the IntelliStart Fluidics console.

The syringe infuses solutions at a controllable flow rate from one of three sample reservoirs A, B, C, or the wash reservoir. The syringe is a displacement type and therefore requires the wash solution for priming. The system primes the syringe automatically.

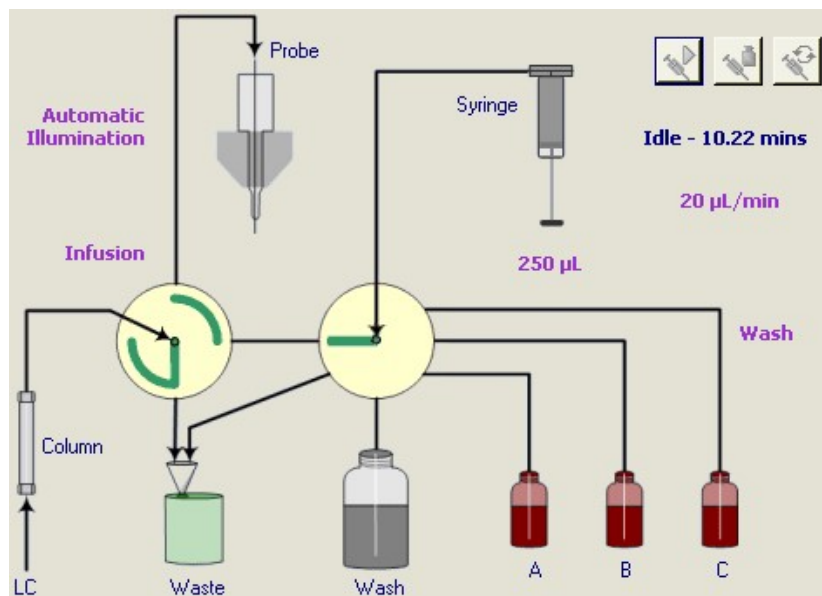

You can manually operate the solvent delivery system from the console, or automate it as part of an IntelliStart sequence or timed event in the MS instrument method. The console's Interactive Display tab provides a snapshot view of the system.

The recommended reservoir usage is as follows:

- A – Sample
- B – LockSpray reference
- C – Calibrant




Notice: To avoid damage, ensure that the wash reservoir contains a sufficient volume of a suitable wash solution.

Requirement: Ensure that reservoirs A, B, C, and the wash reservoir are filled with the appropriate volume of solution for your application.

Tip: The wash solution is sample-dependent. However, for most applications, Waters recommends 80:20 methanol/water.

To control the solvent delivery system from the console:

1. Click Control .
2. Click the current flow rate and in the Set Flow Rate dialog box, enter a new flow rate.
3. Click the currently selected reservoir A, B, C, or wash.
4. In the Select Reservoir dialog box click A, B, C, or wash.

Tip: The system purges if the reservoir changes.

5. Click the currently selected valve position label.
6. In the Select Flow State dialog box select a flow state.

The system's diverter valve operates in the following flow states:

Flow states:

Flow state	Description
LC	The flow goes from the LC to the mass spectrometer (that is, the syringe flow is diverted to the waste bottle).
Infusion	The flow goes directly from the syringe to the mass spectrometer's probe, and the LC flow goes to the waste bottle.
Combined	The syringe flow combines with the LC flow to the mass spectrometer's probe.
Waste	All the flow (LC and syringe) goes to the waste bottle


7. Click  to start.

Result: The status display indicates the amount of fluid in the syringe and the amount of time left before the syringe empties. When the syringe is empty, the system becomes idle.

8. Click  to refill the syringe or  to purge the system.


The Purge command flushes out the system with wash solution and then fills the syringe with the selected sample.

Tip: Purge the infusion syringe whenever you replace a solution bottle. For a sample volume of 250 μL , system purge uses a total volume of approximately 550 μL .

9. Click  to stop the current action.

Tip: Experienced users can also control the solvent delivery system from the Tune window.

To change the solvent delivery system illumination:

1. Click Control .
2. Click the current illumination.

- In the dialog box, select the illumination from the list.

Option	Description
Automatic illumination	The reservoir in use is lit.
All on	All reservoirs are lit.
All off	No reservoirs are lit.

[Back to top](#)

Monitoring and controlling LockSpray flow

You can monitor the status of the lock-mass spray system from the IntelliStart fluidics console. The lock-mass system consists of a multi-position selector valve with input connections from the sample reservoirs and wash solution. The system's output goes to the mass spectrometer's LockSpray reference probe and a waste line. The system infuses solutions at a controllable flow rate from one of three sample reservoirs (A, B, C) or the wash reservoir.

Lock-mass information on the IntelliStart fluidics panel:


Reference Fluidics:	
Selected Reservoir:	C Wait...
Flow State:	Infusion 10 µL/min
Ref Syringe State:	Empty

To control the solvent delivery system from the console:

- !** **Notice:** To avoid damage, ensure that the wash reservoir contains a sufficient volume of a suitable wash solution.

Requirement: Ensure that reservoirs A, B, C, and the wash reservoir are filled with the appropriate volume of solution for your application.

Tip: The wash solution is sample-dependent. However, for most applications, Waters recommends 80:20 methanol/water.

- Click Control .
- Click the flow rate.
- In the Set Flow Rate dialog box, enter a new flow rate.
- Click the currently selected reservoir: A, B, or C.
- In the Select Reservoir dialog box click A, B, or C.

Tip: If the selected reservoir differs from the current reservoir, the system purges as the reservoir changes.

- Click the currently selected valve position label.
- In the Select Flow State dialog box, select one of these flow states:

Selector valve flow state	Description
Infusion	The flow is direct from the reservoir to the LockSpray probe


Selector valve flow state	Description
Waste	All the flow goes to the waste bottle

8. Click  to start.

Result: The status display indicates the amount of fluid in the syringe and the amount of time remaining before the syringe empties. When the syringe is empty, the system becomes idle.

9. Click  to refill the syringe or  to purge the system.

Tip: The fill volume for the LockSpray is fixed at 500 µL.

10. Click  to stop the current action.


Prepare the instrument for use


▪ → [MassLynx](#) > [MS Console](#) > [Instruments: SYNAPT G2-Si](#)

Perform the following tasks to prepare the SYNAPT G2-Si for use.

Requirement: The instrument is configured for your application.



To prepare for data acquisition:

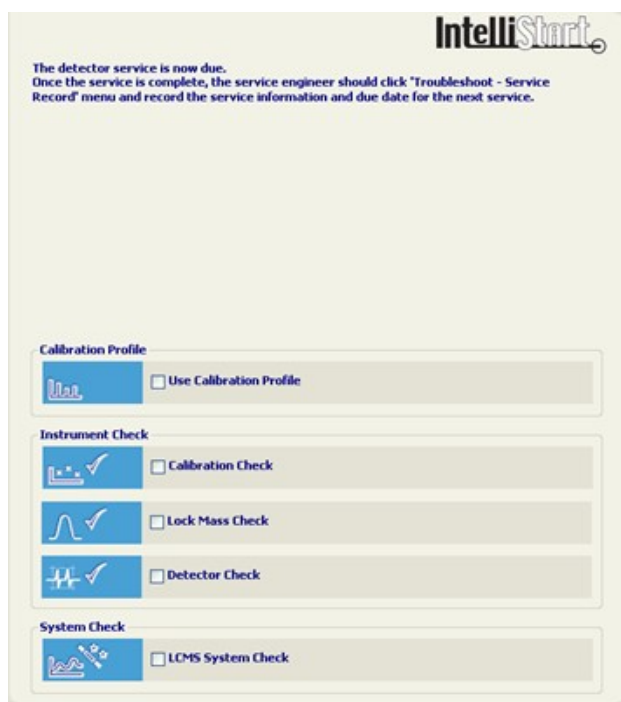1. In the MassLynx opening window, click MS Console 
2. Confirm that the Power and Status LEDs are showing green.
3. In left hand pane tree structure, select Instruments >SYNAPT G2-Si > Intellistart
4. In the Configure menu click Normal mode.



5. Select Detector Check, click Start , and run the wizard.

Tip: The detector requires checking only periodically. If the check fails, run Detector Setup (see *Related Topics*).

6. Select a calibration by selecting Use Calibration Profile, clicking .
7. In the wizard, select an active calibration profile from among those configured for the instrument.
8. Check the calibration by selecting Calibration Check, clicking , and following the wizard.

Recommendation: Check your calibration regularly and in particular, after performing these tasks:

- Cleaning
- Maintenance
- A change in TOF mass range

Tip: A change of mass range can also require a corresponding change of calibration.

If the calibration check fails, then calibrate the instrument to correct it.



9. Check the lock-mass intensity by selecting lock-mass Check, clicking , and running the wizard.

Tip: The lock-mass signal intensity and infusion flow rate required to achieve good stability depend on these conditions:

- Desolvation temperature
- Gas flow
- LC flow rate

- Solvent composition

Check the lock-mass setup specifying MS source conditions that mimic those used in the analysis. That is, with LC flow and, ideally, a flow rate and solvent composition at the midpoint of their gradient.

10. Check the LCMS system by selecting LCMS System Check, clicking , and running the wizard.
11. Select a system check type from the list.
12. Enter the reservoir bottle designation (that is, the position in the Auto sampler) of your sample and pre-run sample.
13. To run the check immediately, click Start .

Provided the checks are successful, the SYNAPT G2-Si is ready to run your analysis.

Related Topics

- [Calibration, Lockspray or LC system failure](#)
- [Configure a lock-mass for mass accuracy](#)
- [Check detector condition](#)

Schedule LCMS system Checks

- [MassLynx > MS Console > Instrument > Intellistart](#)

During a system check, MassLynx acquires data from an LC injection and compares the results with the value of acceptance parameters including retention time, peak area, peak height, peak width and the signal to noise ratio.

You can use IntelliStart to run regular tests to verify the performance of the LCMS system. You can activate an LCMS system check in two ways:

- Manually, from the console
- Automatically, according to a scheduled time and date

To enable operators to run LCMS system checks, you can save defined configurations for manual use later.

The OpenLynx Application Manager processes multiple injections and reports measurements for each specified parameter in all repeat analyses. If results fall within tolerances that you define, the LCMS system is ready to use.

The LCMS system check results are logged via the console, but they can appear in a printed report. The raw data and OpenLynx browser report are also stored for records.

Requirement: For the LCMS system check to work, OpenLynx must be installed with the MassLynx software (it is normally installed by default). If OpenLynx is not installed, rerun the installation DVD, click Modify, and select OpenLynx when prompted.

Before running a system check, you need a suitable tune file, method file, and inlet file. IntelliStart reads the tune, inlet, and method files from the Acqddb folder of the LCMS System Check project (SystemQC.pro). You must create an inlet file independently using the MassLynx Inlet Editor.

To set up a system check:

1. In the MassLynx Instruments window click MS Console .
2. In the left pane directory structure select Instruments > SYNAPT G2-Si> Intellistart.

3. Select LCMS System Check, and click Start .

Tip: This IntelliStart option is available in normal mode. If you cannot see the option, click Configure > Normal mode.

4. In the IntelliStart Setup window, LCMS System Check tab, click "Switch to advanced mode".
5. On the LCMS System Check tab, enter the number of pre-run and run injections, the reservoir bottle designation (position in the Autosampler of your sample), and the injection volume.

Tip: A pre-run sample cleans and stabilizes the system. You need not perform a pre-run prior to a run (pre-run sample injections are not included in the calculated results).

Recommendation: Make a minimum of three injections for tests where standard deviation calculations are required. If you select a single injection, acceptance parameters based on standard deviation are unavailable.

6. In the Method box, specify a tune file, inlet file, and MS file, or use the browse button to select a file.
7. Select the check boxes for the LC acceptance criteria you want to enable, and enter acceptance parameters for these the criteria:

LC acceptance parameters:

Parameter	Description
Retention Time	Enter the expected retention time in the Set Point box, and enter the tolerance in the Tolerance box. If the peak top falls outside the range for any samples, the test results in a failure. If you select three or more injections, the variation in retention time is calculated as the standard deviation (SD). Select the SD check box, and enter a tolerance above which the test results in a failure.
Peak Area	Enter the expected peak area in the Set Point box, and enter the tolerance in the Tolerance box. If the peak area falls outside the range for any sample, the test results in a failure. If you select three or more injections, the variation in peak area is calculated by relative standard deviation (%RSD) = (standard deviation/mean area)×100%. Select the %RSD check box, and enter a tolerance above which the test results in a failure.
Peak Height	Enter the expected peak height in the Set Point box, and enter the tolerance in the Tolerance box. If the peak height falls outside the range for any sample, the test results in a failure. The variation in peak height is calculated by relative standard deviation (%RSD) = (standard deviation/mean area)×100%. Enter a tolerance in %RSD above which the test results in a failure.
Peak Width	Enter the expected peak width in the Set Point box, and enter the tolerance in the Tolerance box. If the peak width falls outside the range for any sample, the test results in a failure.
Signal/Noise	In the Signal/Noise box enter the signal (peak height)-to-noise ratio below which the test results in a failure. In the Noise Range box, enter the range in retention time window over which the noise should be calculated. Select the method by which the noise is calculated: RMS or Peak-to-Peak

Tip: To determine suitable acceptance criteria, run an injection of your LCMS System sample, and review the chromatogram in MassLynx software.

8. Set your MS Acceptance Criteria with an expected Mass peak from your sample and specify an acceptable tolerance for your purposes in ppm or mDa.

9. Click Set Schedule, and choose Off, Daily, or Weekly. For Weekly, select the required days of the week.
10. Select Print Report, if you want a printed report of the results.

Tip: When you make changes, the Save As button, in the top, right-hand portion of the tab becomes active.

11. Click Save As, type a name for the system check type, and click Save.
12. Click Start, to run the system check immediately, or click OK, to close the dialog box, and run the LCMS system check according to your schedule.

Result: The status bar at the bottom of the IntelliStart window reports the progress of the system check.

Resolving system check failure

If the system check process fails, IntelliStart displays a warning light and the system becomes unavailable for use.

Click Resolve to clear the system check error and put the instrument back into a state ready for data acquisition.

The calculations do not include injections where no peak was found. For example, if only 8 of 10 injections detect peaks, the test results are based on 8 values rather than 10. However, the system check process is marked as failed because of the missing peaks. Investigate the cause of the system check failure, and run the system check again before proceeding with other analyses.

Related Topics

- [Prepare the instrument for use](#)
- [Troubleshoot the instrument](#)

Schedule LS MS system checks

- [→ MassLynx > MS Console > Instruments: SYNAPT G2-Si > Intellistart](#)

During a system check, MassLynx acquires data from an LC injection and compares the results with the value of acceptance parameters including retention time, peak area, peak height, peak width and the signal to noise ratio.

You can use IntelliStart to run regular tests to verify the performance of the LCMS system. You can activate an LCMS system check in two ways:

- Manually, from the console
- Automatically, according to a scheduled time and date

To enable operators to run LCMS system checks, you can save defined configurations for manual use later.


The OpenLynx Application Manager processes multiple injections and reports measurements for each specified parameter in all repeat analyses. If results fall within tolerances that you define, the LCMS system is ready to use.

The LCMS system check results are logged via the console, but they can appear in a printed report. The raw data and OpenLynx browser report are also stored for records.

Requirement: For the LCMS system check to work, OpenLynx must be installed with the MassLynx software (it is normally installed by default). If OpenLynx is not installed, rerun the installation DVD, click Modify, and select OpenLynx when prompted.

Before running a system check, you need a suitable tune file, method file, and inlet file. IntelliStart reads the tune, inlet, and method files from the Acqddb folder of the LCMS System Check project (SystemQC.pro). You must create an inlet file independently using the MassLynx Inlet Editor.

To set up a system check:

1. In the MassLynx Instruments window click MS Console .
2. In the left pane directory structure select Instruments > SYNAPT G2-Si > Intellistart.
3. Select LCMS System Check, and click Start .

Tip: This IntelliStart option is available in normal mode. If you cannot see the option, click Configure > Normal mode.

4. In the IntelliStart Setup window, LCMS System Check tab, click "Switch to advanced mode".
5. On the LCMS System Check tab, enter the number of pre-run and run injections, the reservoir bottle designation (position in the Autosampler of your sample), and the injection volume.

Tip: A pre-run sample cleans and stabilizes the system. You need not perform a pre-run prior to a run (pre-run sample injections are not included in the calculated results).

Recommendation: Make a minimum of three injections for tests where standard deviation calculations are required. If you select a single injection, acceptance parameters based on standard deviation are unavailable.

6. In the Method box, specify a tune file, inlet file, and MS file, or use the browse button to select a file.
7. Select the check boxes for the LC acceptance criteria you want to enable, and enter acceptance parameters for these the criteria:

LC acceptance parameters:

Parameter	Description
Retention Time	Enter the expected retention time in the Set Point box, and enter the tolerance in the Tolerance box. If the peak top falls outside the range for any samples, the test results in a failure. If you select three or more injections, the variation in retention time is calculated as the standard deviation (SD). Select the SD check box, and enter a tolerance above which the test results in a failure.
Peak Area	Enter the expected peak area in the Set Point box, and enter the tolerance in the Tolerance box. If the peak area falls outside the range for any sample, the test results in a failure. If you select three or more injections, the variation in peak area is calculated by relative standard deviation (%RSD) = (standard deviation/mean area)×100%. Select the %RSD check box, and enter a tolerance above which the test results in a failure.
Peak Height	Enter the expected peak height in the Set Point box, and enter the tolerance in the Tolerance box. If the peak height falls outside the range for any sample, the test results in a failure. The variation in peak height is calculated by relative standard deviation (%RSD) = (standard deviation/mean area)×100%. Enter a tolerance in %RSD above which the test results in a failure.
Peak Width	Enter the expected peak width in the Set Point box, and enter the tolerance in the Tolerance box. If the peak width falls outside the range for any sample, the test results in a failure.
Signal/Noise	In the Signal/Noise box enter the signal (peak height)-to-noise ratio below which the test results in a failure. In the Noise Range box, enter the range in retention time window over which the noise should be calculated. Select the method by which the noise is calculated: RMS or Peak-to-Peak

Tip: To determine suitable acceptance criteria, run an injection of your LCMS System sample, and review the chromatogram in MassLynx software.

8. Set your MS Acceptance Criteria with an expected Mass peak from your sample and specify an acceptable tolerance for your purposes in ppm or mDa.
9. Click Set Schedule, and choose Off, Daily, or Weekly. For Weekly, select the required days of the week.
10. Select Print Report, if you want a printed report of the results.

Tip: When you make changes, the Save As button, in the top, right-hand portion of the tab becomes active.

11. Click Save As, type a name for the system check type, and click Save.
12. Click Start, to run the system check immediately, or click OK, to close the dialog box, and run the LCMS system check according to your schedule.

Result: The status bar at the bottom of the IntelliStart window reports the progress of the system check.

Resolving system check failure

If the system check process fails, IntelliStart displays a warning light and the system becomes unavailable for use.

Click Resolve to clear the system check error and put the instrument back into a state ready for data acquisition.

The calculations do not include injections where no peak was found. For example, if only 8 of 10 injections detect peaks, the test results are based on 8 values rather than 10. However, the system check process is marked as failed because of the missing peaks. Investigate the cause of the system check failure, and run the system check again before proceeding with other analyses.

Related Topics

- [Prepare the instrument for use](#)
- [Troubleshoot the instrument](#)

Starting and stopping the instrument






- [MassLynx > MS Console > Instruments: > SYNAPT G2-Si Operate](#)


This topic describes the conditions required to start the SYNAPT G2-Si HDMS SYNAPT G2-Si. When you are not using the instrument, to conserve energy and reduce nitrogen consumption, stop the LC flow, and setting the instrument to [Source Standby](#).

To start the SYNAPT G2-Si:

The procedure for starting the SYNAPT G2-Si depends on the instrument's initial state.

Instrument state	Starting procedure
------------------	--------------------



Instrument state	Starting procedure
Powered-off	<p>Starting the Instrument</p> <ol style="list-style-type: none"> 1. Start the host (Work station) PC, wait for it to finish booting, and log in. 2. Power-on the instrument. <p>Turning instrument power on (with ACQUITY UPLC)</p> <p>! Notice: The ACQUITY UPLC module that contains the Ethernet switch must be turned on first to allow communications to be established.</p> <ol style="list-style-type: none"> 1. Turn on power to the ACQUITY module containing the Ethernet switch (for example the high temperature heater/cooler or the column manager). 2. Use the breakers at the rear of the instrument to switch on power to the instrument. 3. Check the instrument side panel illumination LED panel turns white. 4. Switch on power to the binary solvent manager. 5. Switch on power to any additional UPLC modules. 6. Allow up to three minutes for the EPC to boot successfully. <p>Turning instrument power on (with nanoACQUITY UPLC)</p> <ol style="list-style-type: none"> 1. Turn on power to the Ethernet switch. 2. Use the breakers at the rear of the instrument to switch on power to the instrument. 3. Check the instrument illumination LED panel turns white. 4. Switch on power to the binary solvent manager (if not already on). 5. Switch on power to any additional UPLC modules. 6. Allow up to three minutes for the EPC to boot successfully. <p>Starting the instrument (continued)</p> <ol style="list-style-type: none"> 3. Start the MassLynx software. 4. Click MS Console  5. In the Console window, click Operate  6. Depending on how long the instrument has been idle and in a vented state, it can require more than 24 hours to achieve sufficient vacuum.
Powered-on, but vented	<ol style="list-style-type: none"> 1. Click  2. Wait for the instrument to become ready. <p>Depending on how long the instrument has been resting in a vented state, it can require more than 24 hours to achieve sufficient vacuum.</p>
Instrument pumped down (evacuated) and in Standby mode	<ol style="list-style-type: none"> 1. Click  2. Wait 1 hour for the high voltages to stabilize and source temperatures to reach equilibrium.
Instrument in Source Standby mode	<ol style="list-style-type: none"> 1. Click  2. Wait 30 minutes for source temperatures to reach equilibrium.

Result: After you click , these events take place:

- The Operate icon flashes red.
- Pumping down (evacuation) commences.
- The API and collision gas start.
- Eventually, the Operate icon changes to green, and the instrument is ready for use.

Stop the instrument by selecting Standby status

When you are not using the instrument, to conserve energy and reduce nitrogen consumption, stop the LC flow, and put the instrument in Source Standby mode.

1. In the Instrument window select MS Tune  .
2. Click Source Standby  .

Result: These things happen when you click Source Standby:

- The high voltage in the source switches off.
- The sample flow is diverted to waste.
- The probe heater switches off.
- The ion block heater continues to operate.
- Source gases switch off.



Warning: To avoid burn injuries, exposure to biohazards, or toxic substances, see the *SYNAPT G2-Si Operator's Overview and Maintenance Guide* before carrying out any physical operations on the source.

Related Topics

- [Vent and shut down](#)
- [Prepare the for use](#)
- [Acquire data](#)

Vent and shutdown the instrument

- [MassLynx > MS Tune > Tune Window](#)

Venting and shutting down the SYNAPT G2-Si is an advanced/service user level task. Attempt this operation only if you understand its consequences. For further details, including procedures for restarting the instrument, see the *SYNAPT G2-Si Operator's Overview and Maintenance Guide*.

When it is not operating, allow the SYNAPT G2-Si to remain in Source Standby mode (see Related Topics). When it is necessary to completely switch off the instrument, vent it before disconnecting the power.

To vent the instrument:

1. If the instrument is in Operate mode, click Source Standby .
2. From the Tune window menu, click Vacuum > Vent.

When a pump or vent command is issued, status of the operation is shown.

To shut down the mass spectrometer:



Warning: To avoid electric shock, disconnect the power cable from the instrument's rear panel.

Disconnect the power cable from the electrical outlet.

Related Topics

- [Startup the SYNAPT G2-Si](#)

Calibrating the DESI zero position

- [MassLynx> Instrument > MS Tune >DESI](#)

All calculations for the automated DESI imaging are carried out based on a universal zero position.

Because the home position of the stage does not reflect the position of the sprayer, and therefore the analysis point on the surface, it is necessary to create a X-Y start point within MassLynx.

The start point is taken as the top left corner of a glass slide placed on the stage, into slot A.

Tip: To ease setting up the X-Y start point, a dot of red ink from the Sharpie® marker-pen provided in the installation kit, can be placed on the upper, left corner, of the glass slide before putting the slide into slot A.

1. On the tune page, click DESI source control.
2. Make sure the device status shows connected, if not, click connect device.
3. Click move to point.
4. Make sure that the sprayer is sampling the top left corner of the slide in slot A

Tip: If a dot of red ink from the Sharpie® marker-pen has been placed in the upper, left corner, an intense peak of m/z 433 in positive ion mode should be seen.

5. Make sure the sprayer is in its correct position. if not, alter the slot-start point X-Y and press move to point.
6. Make sure the device status shows connected, if not, click connect device.
7. Repeat step 6 as required.

Tip: The DESI spray position is now set and ready for automated image definition through HDI software.

Related Topics

- [Defining a DESI imaging experiment](#)
- [Selecting the DESI source](#)
- [Acquire with DESI](#)
- [Calibrate with DESI](#)

Defining a DESI imaging experiment

- [HDI imaging> HDI software](#)

You use HDI software to configure the DESI sample plate by positioning the tip of the DESI spray head above the plate. The tip of the DESI spray head must be accurately positioned above the plate and the start point parameters saved. And changes made to the angle of spray will also change the tip position.

You can export the HDI configuration to MassLynx software, for use in experiments.

HDI Imaging

Applying reference points

1. Open the HDI software.
2. On the Acquire page, click the Pattern tab.
3. Click create new plate
4. From the image file, select a saved plate.
5. Select the DESI sample source.
6. Select the appropriate slide from the holder.
7. Select the appropriate plate type.
8. Click next.
9. Select the four corners of the slide to establish reference points on the slide image.
10. Click next
11. Preview the slide.
12. Click finish.
 - **Recommendation:** Use the default pixel sizes.
14. Select the pixel sizes, if necessary.
 - **Recommendation:** Tick use single value for x and y pixels.
16. Select the rate and scan times, if necessary.
17. On the Acquire page, click Acquire.
18. Select the appropriate target instrument type.
19. Select the experiment as required.
20. In the active window, select an area of interest using the rectangle tool.
21. On the Acquire page, select the Process tab and make sure that the processing parameters are correct.
22. Save the metafile as a .pdm file.
23. Use the MassLynx button in the toolbar and export the software to MassLynx.

Related Topics

- [Calibrate with DESI](#)
- [Acquire with DESI](#)
- [Selecting the DESI source](#)
- [Calibrating the DESI zero position](#)

Selecting the DESI source

- [MassLynx](#) > [Instrument](#) > [MS Tune](#) > [DESI](#)

Before you can access the DESI source via MassLynx, you must first install the source on the instrument. See the *DESI Operator's Guide* (P/N 715004701) for details.

Selecting the DESI source

Requirement: The DESI source page controls the DESI control box. The USB cable from the DESI control box must be installed on the PC.

The axis parameters are set up in the topic, [Configure the DESI sample plate](#).

Requirement: The tip of the DESI spray head must be accurately positioned above the plate and the start point parameters saved. And changes made to the angle of spray will also change the tip position.

1. In the Tune window, click [Source](#) > [DESI](#).
2. In the DESI window, select the [DESI Source Control](#).
3. Select the [Desi motion parameters](#)
4. Verify and adjust, as necessary, the axis parameters and the general parameters.
5. Select [OK](#).

Related Topics

- [Defining a DESI imaging experiment](#)
- [Calibrate with DESI](#)
- [Acquire with DESI](#)
- [Calibrating the DESI zero position](#)

Operate the Maldi source

- [MassLynx](#) > [Instrument](#) > [MS Tune](#) > [Source](#) > [MALDI](#)

Before you can access the MALDI source through MassLynx, you must first physically install the source on the instrument. See the *MALDI SYNAPT G2-Si Operator's Overview and Maintenance Guide* for details.

Selecting the MALDI source

In the Tune window, click [Source](#) > [MALDI](#)

Result: The MALDI source page is displayed.

Requirement: The source activation cable must be connected before you can access the MALDI source page.

In addition to the tasks you can perform for a regular (non MALDI) SYNAPT G2-Si, you can also perform these additional functions from the MALDI source page:


- Load sample plates and select sample wells.
- Display a camera image of the sample.
- Fire the laser and optimize the spectrum.
- Save and recall settings for use with different samples.
- Set parameters on the MALDI source page.

Tip: Once installed, you can put the MALDI source into Standby and Operate modes, as with any other instrument source.


Using the camera viewer

Use the camera to view the sample plate in position. The camera image of the selected sample well shows cross hairs indicating the firing point. You can see the point at which the laser fired on the sample as a bare patch, indicating that the sample was consumed.

The MALDI camera viewer:

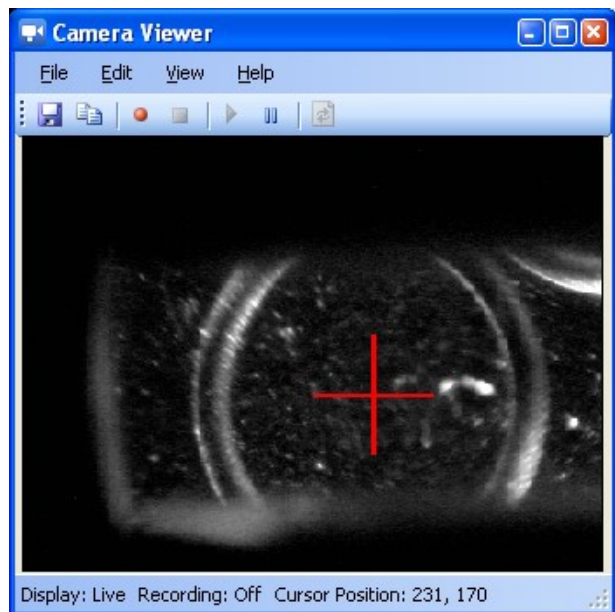

To operate the camera:

1. Click  to open the camera viewer.
2. Click the status bar to toggle between Live Image and Static Image modes.
3. To move the camera position, click the status bar to select Live Image mode, and then click and drag on the camera image.

Rule: In Static Image mode, you can move the cross hairs by clicking in the image. In Live Image mode, you cannot move the cross hairs.


You can also save a still image of the camera display by clicking the disc icon in the toolbar.

See also: MALDI Imaging online Help for information on aligning the camera image and aligning the laser and sample well.

Aligning the camera image with the laser position

When you fire the laser, the laser spot is not always perfectly aligned with the cross hairs. You can correct the misalignment from the camera viewer.

To align the camera image with the laser position:

1. In the camera viewer, click the status bar to select Live Image mode.
2. Move the camera position to view a prepared sample well.
3. In the Tune window, click .

Tip: Do not move the laser until enough material is removed from the sample plate so that you can see a laser spot.

4. Click the status bar to change the camera status to Static Image mode.
5. Click on the laser spot in the camera image to move the camera cross hairs to the correct position.
6. When the laser and cross hair positions are satisfactorily aligned with each other, click in the status bar to toggle back to Live Image mode.

Result: The camera software stores the cross hair position.



Aligning the camera image and the sample well

After you align the camera cross hairs to indicate the laser firing position, the center of the sample well does not always align with the camera cross hairs. You can align the sample well position relative to the camera image from the Tune window.

To realign the camera image with the sample well position:

1. From the Tune window, Click MALDI > Source Settings > Password
2. Type "access", and then click Apply.
3. In the Source Settings dialog box, click Sample Plate, and select Nudge Using Cross hairs.

Result: The plate moves to sample well A1.

4. Adjust the sample plate cross hairs until the center of the sample well aligns with the camera cross hairs.
5. Click  to store the current alignment.
6. Repeat steps 3 and 4 until the sample well center and the laser firing mark align satisfactorily.
7. Click  to accept the current alignment.

Result: The plate re indexes and moves to sample A1 using the latest alignment settings.

Related Topics

- [Calibrate using a MALDI source](#)
- [Create methods for MALDI](#)
- [Use MALDI pattern editor](#)
- [Perform a MALDI acquisition](#)

Controlling the MALDI sample plate

You can control the sample plate during an acquisition from the MALDI control window or via the MassLynx sample list. In the MALDI control window, you control the sample plate interactively, using a cross hair pointer or by following a predefined pattern. If you acquire data via the sample list, you can only use pattern control.

Creating patterns

A pattern comprises a series of nodes that define a series of 2D spatial positions between which the plate carrier moves. When the end of the pattern is reached, the movement stops. If you run an acquisition from the sample list, and the end of a pattern is reached, the acquisition stops. No limit applies to the number of nodes you can use in a pattern.

You can use three types of pattern file: PTN, PAT and PTO.

PTN files:

Pattern files, with the extension ".ptn" serve only non imaging experiments. You can load various predefined patterns, either spiral or straight-line, from the pattern editor. Additionally, you can define your own patterns. Pattern files are stored as simple comma-delimited text files that you can open and modify in any suitable application (such as a spreadsheet). Each pattern is scaled to fit the sample wells. For example, the same number of nodes are associated with both a 2.5-mm well and a 2-mm well.

PAT files:

Pattern files, with the extension ".pat" serve imaging experiments only. You can create these files with the MALDI Imaging Pattern Creator software. See the help system included with this software for more details.


PTO files:

Pattern files, with the extension ".pto" serve imaging experiments, although you can also use them in a non imaging context. With these files, you can run an imaging pattern that is offset by the position of the sample well. The software currently does not support creation of this pattern type. Nevertheless, several predefined .pto files are provided. If you add a .pto file into the sample list for a non imaging acquisition, it converts to the .ptn format before it runs.

Loading sample plates and selecting sample wells

Tip: You need not select Standby mode to load and unload plates.

To load plates:

1. Press the Open button.
2. Hold the sample plate by the sides, orienting it so that the corner indentation faces the bottom left-hand side, and place it into the recess in the plate carrier.
3. Close the plate carrier lid by lifting it manually until the magnetic clasps engage.
4. From the Tune window, click .

Result: The system evacuates the plate carrier air lock and transports the sample plate into the source.

5. Enter the sample well number in the sample box, to align the laser beam with a sample well.

Results:

- The selected sample well aligns with the laser beam.
- After the sample plate is loaded, the Unload button glows blue, to show that the unload function is operable.

Tip: Neither button glows when the instrument is acquiring data.

To unload a plate:


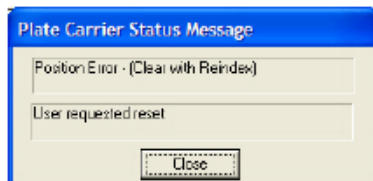
1. From the Tune window, click .
2. When the carrier stops moving, open the plate carrier lid, and then remove the plate, taking care not to touch the surface.

Plate carrier errors


The current status of the plate carrier appears in the Tune window's status bar. Error messages appear in the Plate Carrier Status Message dialog box. If you connect the source activation cable when a plate is already loaded, the following status message appears:



The message indicates that the sample plate control is restarting. If the plate carrier is in the unload position, this message does not appear.

Other errors give possible solutions in the dialog box and are indicated by the state of the Load/Unload buttons.

Load/Unload button states:

State	Meaning
Blue	Normal operation - Load/Unload
Red 	Recoverable error. The plate is out of position. Click the buttons to clear the error and reindex the plate.
Both disabled	Indicates one of the following conditions: <ul style="list-style-type: none"> • Acquiring data • Vented error • Plate carrier error

Sample plate formats

The MALDI SYNAPT G2-Si comes with several plate formats already defined. These are the standard 96-well MALDI sample plate, 384-well, BigSpot, and Imaging plate. See the MALDI Imaging online Help for information on using plate formats.

Viewing and setting plate formats

In the Tune window, click MALDI > Sample Plate Definition to open the Sample Plate Definition dialog box.

To select a predefined plate:

1. Select a plate from the Plate Type list.
2. Click OK.

To create a user defined plate:


1. Click New.
2. Enter relevant parameters for the plate type you are using.
3. In the Index Type drop-down list, select Numerical or Alpha Numeric addressing of the sample wells.
4. To define wells for lock mass correction, select the Use Lock Mass Wells check box, and specify the relevant parameter setting.
5. Click Save As, and save with an appropriate name.

To open a user defined plate definition:

1. Click Open.
2. From the Open dialog box select the required *.mtp file.
3. Check the settings are correct, and click OK.

Sample log sheet

A sample log sheet, MALDI_SAMPLE_LOG.pdf, is included with your MassLynx installation in the MassLynx folder.

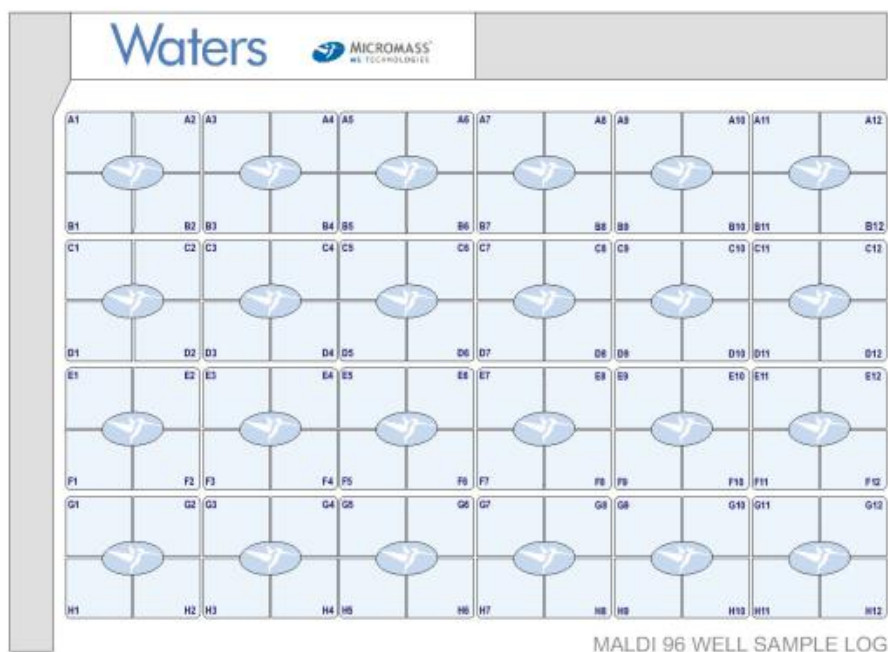Sample log sheet:

Using Adobe Reader, you can print log sheets, and record sample details.

Related Topics

- Operating the MALDI source
- Creating methods for MALDI
- Using the MALDI pattern editor
- Performing a MALDI acquisition

Selecting the REIMS source

- MassLynx > Instrument > MS Tune > REIMS

Before you can access the REIMS source via MassLynx, you must first install the source on the instrument. See the *REIMS Operator's Guide* (P/N 715004683) for details.

Heater select button.

The REIMS source uses an on-source heater. The heater is controlled from the REIMS tune window.

The heater will automatically default to the state it was in, before a standby selection, when the instrument is next put into operate.

When the REIMS source is in operation the heater, when switched on, is controlled automatically by software. The heater will automatically switch to OFF when the following conditions apply:

- The source door interlock is open.
- The backing pressure falls too low.
- The instrument is in "Standby".
- The instrument is vented.

Note: A caption is shown, stating that the heater is disabled because of insufficient backing pressure, indicating a possible blocked capillary.

Selecting the REIMS source

Note: The Heater Bias voltage is unique to REIMS. Both the Cone and the Heater Bias parameters are tunable

1. In the Tune window, click Source > REIMS.
2. As necessary, set the Cone and Heater Bias voltages using the sliders.
3. Select the Heater button to ON.

Related Topics


- [Calibration with REIMS](#)
- [Detector setup page and the REIMS source](#)

Add lockmass correction to methods

- [MassLynx > MS Method > LockSpray](#)

Where your source enclosure supports it, you can use a lock-mass. If you choose to acquire lock-spray data, you can apply corrections using the lock-mass as a reference. You can select only configured lock-mass compounds, those defined in the LockMass editor and with source conditions established.

To use lock-mass correction:

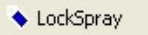1. Click  LockSpray and select one of these lock-mass correction options:
 - Acquire LockSpray - Apply correction
 - Acquire LockSpray - Do not apply correction
 - No LockSpray

Tip: If you do not want to use lock-mass correction, select No LockSpray, and click OK to exit back to the MS Method screen.

2. Select your LockSpray lock-mass compound and complete the LockSpray Acquisition Settings according to the table.

LockSpray Settings:

Parameter	Description	Typical short-run	Typical long-run
Scan Time	The dwell time for the LockSpray acquisition	0.3 s	1.0 s
Interval	The time between LockSpray acquisitions	20 s	60 s
Scans to average	Specifies the size of the rolling average used for the LockSpray mass correction	3	0.3 Da
Mass Window +/-	Adjusts the window size to locate the correct lock-mass peak in the LockSpray spectrum	1	0.3 Da

Related Topics

- [MS and MS/MS methods from functions](#)
- [Creating product ion discovery \(PID\) experiments](#)
- [MS^E method](#)
- [Fast DDA method](#)

Automatic calibration

- [Console System Tree > Instruments: SYNAPT G2-Si > Intellistart](#)

The SYNAPT G2-Si mass scale gives nominal mass measurement without a calibration. However, it must be calibrated before using it for exact-mass measurement. Most situations call for automatic calibration as described in this topic.

MassLynx software uses a polynomial equation to calibrate the system precisely over a wide mass range. Calibrating and checking calibrations always produces a report (stored in C:\MassLynx\IntelliStart\Results\Reports). You can optionally select Display report during the wizard to read the report when the calibration or calibration check finishes.

Requirement: For correct mass measurement performance, the TOF mass range during analysis should be within its calibration range. A warning is displayed during acquisition when this is not the case.

1. Set the Dynamic range to normal when calibrating the instrument. Extended and normal dynamic range acquisitions use the same calibrations.
2. You should first check any existing calibration to see if a calibration procedure is necessary.

These subjects are covered:


- [See "Checking your calibration"](#)
- [See "Calibration profiles"](#)
- [See "Performing a calibration"](#)

Checking your calibration

You should check your calibration regularly and in particular after performing these tasks:

- Cleaning
- Maintenance
- Changing the TOF mass range

To check the calibration:

1. In IntelliStart, click Configure > Normal Mode.
2. Select Calibration Check, then click Start , and follow the wizard.

Tip: If the calibration check fails your criteria, perform a calibration to correct it.

[Back to top](#)

Calibration profiles

Use the calibration profile editor to perform these calibration tasks:

- [Create new calibration profiles](#)
- [Edit existing calibration profiles](#)
- [Resetting existing calibration profiles](#)
- [Check the properties of existing calibration profiles](#)
- [Display calibration reports](#)

Creating a calibration profile

If the standard calibration profiles are not appropriate, you can create a new calibration profile for your calibrant, mass range, and ion mode.

To create a calibration profile:

1. In the Masslynx Instrument window click MS Console
2. In IntelliStart, click Configure > Configuration mode.
3. Click Configure > Calibration Profile Editor.
4. Click File > New.

Result: The New Calibration Profile dialog is displayed.

5. Enter a Calibration Profile Name.
6. Click Edit and specify a mass range.
7. Select the type of calibration.
8. Enable positive polarity, negative polarity or both (see [Enabling Positive Polarity](#) and [Enabling Negative Polarity](#) below).
9. Click OK to create the calibration profile ready for use during calibration.

Enable Positive Polarity

1. To enable Positive Polarity, click Edit.

Result: The Edit Calibration Reference - Positive dialog is displayed.

2. Select Enabled and then choose a reference compound from the drop-down menu.
3. Click OK to return to the New Calibration Profile dialog.

Enable Negative Polarity

1. To enable Negative Polarity, click Edit.

Result: The Edit Calibration Reference - Negative dialog appears.

2. Select Enabled and then choose a reference compound from the drop-down menu.
3. Click OK to return to the New Calibration Profile dialog.

[Back to top](#)

To edit an existing calibration profile:

1. In IntelliStart, click Configure > Configuration mode.
2. Click Configure > Calibration Profile Editor.
3. Select the profile that you want to edit.
4. Click Edit > Edit, or right-click on the profile and select Edit.

Result: The Edit Calibration Profile dialog appears.

5. Edit the mass range and calibration type using your required settings.

[Back to top](#)

Resetting a calibration profile

Rule: You cannot reset a profile that is currently loaded.

Calibration profiles can be reset to leave the instrument in an uncalibrated state.

To reset a calibration profile:

1. In IntelliStart, click Configure > Configuration mode.
2. Click Configure > Calibration Profile Editor.
3. Select the profile that you want to reset.
4. Click Edit > Reset, or right-click on the profile and select Reset.

Result: A warning message appears.

5. Click Yes to reset, then click File > Exit to leave the Calibration profile editor

[Back to top](#)

Checking the properties of a calibration profile**To check the properties of an existing calibration profile:**

1. In IntelliStart, click Configure > Configuration mode.
2. Click Configure > Calibration Profile Editor.
3. Select the profile that you want to check.
4. Click File > Properties, or right-click on the profile and select Properties.

Result: The Calibration Profile Properties dialog for your chosen profile appears.

[Back to top](#)

Displaying calibration reports

To view a summary report for an existing calibration profile:


1. In IntelliStart, click Configure > Configuration mode.
2. Click Configure > Calibration Profile Editor.
3. Select the profile that you want to check.
4. Click File > Properties.
5. Click the Display Report button (bottom left of the dialog box) and select Display Calibration Profile Report.

Result: A browser window opens displaying your chosen report.

[Back to top](#)

Performing a calibration

To perform a calibration:

1. In IntelliStart click Configure > Configuration Mode.
2. Select Create Calibration, then click , and follow the wizard.

The wizard displays the available calibration profiles. If a profile has been successfully calibrated, a green tick is shown next to it. Undefined calibrations have an orange triangle next to them. When both positive and negative are defined in a profile, they must use the same reference compound.

[Back to top](#)

Related Topics


- [Manual calibration](#)

Check and set-up the detector


→ MassLynx > MS Console > Instrument: SYNAPT G2-Si > Intellistart

As the detector ages it requires increasingly higher voltage to achieve the correct gain. When the required voltage reaches its maximum value, the detector needs replacing. Check your detector periodically to maintain the correct detector gain as the detector ages.

To check the detector:

1. Click Configure > Normal Mode
2. Select Detector Setup Check, click Start , and follow the wizard.

To set up the detector:

1. Click Configure > Configuration Mode
2. Select Detector Setup, click , and follow the wizard.

Related Topics

- Troubleshoot the instrument
- Prepare the instrument for use

Configure the APCI interface

- MassLynx > MS Console > SYNAPT G2-Si > MS Tune

This procedure describes how to configure the APCI interface. Note that it assumes the system settings are optimal.

! **Notice:** The Waters installation engineer leaves your instrument with optimal system settings. Do not adjust them. The system settings are password-protected to avoid unintentional changes. Changing them renders your calibrations invalid and reduces the performance of your system.

Note: If you have an ETD-enabled instrument, you must ensure that the glow discharge electrode is removed from the ion block and a blank plug fitted in its place. You must also ensure that the high voltage connector is attached to the ion block. For instructions, see *ETD Source for SYNAPT G2-Si Operator's Overview and Maintenance Guide*.

Configuring the APCI interface

Following the setup, the settings are saved in a tune (.ipr) file. The tune file saves values for all supported ionization modes at once, not just for the APCI mode.

Requirement: The APCI interface is installed on the SYNAPT G2-Si mass spectrometer (see the Waters *Atmospheric Pressure GC Operator's Guide Supplement*).

To set up the GC oven and reservoir bottles:

1. Set the transfer line temperature, controlled by thermal auxiliary 3, to 310.0 °C.
2. Set the nitrogen sheath gas, controlled by auxiliary column 3, to constant flow mode and with a flow rate of 350 mL/min.
3. Seat a reservoir bottle of heptacosane reference compound in one of the reservoir carriers. To prepare reference reservoirs, see the *Waters Atmospheric Pressure GC Operator's Guide Supplement*.


Requirement: Ensure the source enclosure is closed, the transfer line is connected, and the instrument is in Operate mode.

To set up the ionization:

1. From the Tune window, click File > Open, and open the default tune file.

Tip: At any time, you can click File > Revert, to revert to previously saved parameters.

2. If you have an ETD-enabled instrument, click Acq Mode > CID Mode.
3. From the Source menu, click APCI.
4. Select the required ion mode.

Example: click positive  ion mode.

5. Click MS mode  MS .

Tip: MS/MS configuration uses the same source parameter values as MS mode. However, besides the parameter values, you must specify the collision energy and set-mass.

6. On the APGC tab, specify these initial parameters settings.

Source values:

Parameter	Setting
Current mode	On
Corona current	1.3 μ A
Sampling Cone	30 V
Extraction Cone	3.0 V

Temperature settings suitable for sample infusion:

Parameter	Setting
Source	150 °C

Gas flow settings suitable for sample infusion:

Parameter	Setting
Cone Gas	20 L/hr
Auxiliary Gas (used to continually flush the source enclosure)	200 L/hr

7. Consult the Instrument tab to ensure the values are similar to the recommended values.

Instrument settings:

Parameter	Setting
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	MS/MS only

Requirement: Set the collision energy to Off, for MS mode. If setting up for MS/MS operation, you must then fragment a sample afterward, reset the collision energy to On, and specify the appropriate voltage.

8. Set the values in the non storage acquisition pane, to observe the spectrum.

Typical acquisition settings:

Parameter	Setting
Low mass	100 Da
High mass	1200 Da
Scan time	1.0 sec
Data type	Continuum or Centroid
MSMS Mass (only available in MS/MS mode)	Enter a precursor ion mass in Da

9. Zoom the peak display to a suitable range around the mass of interest.

Tip: On the spectrum, right-click, drag across the peak's full width, and release.

10. Optimize the corona current for signal intensity and noise.
11. Optimize the cone voltage for the required fragmentation.
12. Verify that the shape of the displayed peaks are correct.
13. Ensure that resolution, intensity, and stability are acceptable.
14. Click File > Save As, and save the settings using an appropriate name.

Starting with these settings, you can optimize the cone gas and auxiliary gas settings for your application, and save them in a separate tune file. Once the instrument is tuned for your application, you can acquire nominal mass data. To obtain exact mass data, calibrate the instrument, and configure and use a lock mass.

Related Topics

- [Automatic calibration](#)
- [Configure a lock mass](#)

Configure APPI for the optional source

- [MassLynx > MS Console](#)

This procedure explains how to configure the APPI for the optional source. It assumes optimally set system parameters.

! **Notice:** The Waters installation engineer leaves your instrument with optimal system settings. Do not adjust them. The system settings are password-protected to avoid unintentional changes. Changing them renders your calibrations invalid and reduces the performance of your system.

Requirement: If you have an ETD-enabled ion block fitted, ensure that the discharge electrode cable is connected to the ion block connector. For instructions, see *SYNAPT G2-Si Operator's Overview and Maintenance Guide* on the Documentation CD.

This procedure comprises the following two processes:

- [Setting up the beam initially by combining sample infusion](#)
- [Configuring specific applications](#)

Following the initial setup, save the settings as a default, sample-dependent tune file. You can then adjust the source parameters, as appropriate, for your application. To do so, you optimize temperatures and gas flows relative to the LC flow rate, sample type, and solvent type. You can then save these optimized source conditions as a separate tune (.ipr) file, ready for use.

The tune (.ipr) files save values for all supported ionization modes at once, not only the ionization mode selected.

For APPI, supplement the low flow-rate sample by specifying an LC flow that permits appropriate, high-flow, configuration settings. The sample must be typical for your application and taken from the vial mounted on the front of the instrument.

To set up the beam initially by combining sample infusion

Sample infusion is a low flow-rate infusion of a typical sample from the vial mounted on the front of the SYNAPT G2-Si. It allows initial source configuration before switching to the high-flow conditions generated by the LC system.

To set up the beam:

1. Install suitable sample, reference, and calibrant vials in positions A, B and C on the front of the SYNAPT G2-Si.

Tip: For installation instructions, see SYNAPT G2-Si Operator's Overview and Maintenance Guide on the Documentation CD.


Recommendation: Use reservoir bottles in the following order:

- A, for sample
- B, for reference
- C, for calibrant

Tip: Reference and calibrant solutions are described in their respective help topics. Base sample solutions on your eventual application. Filter or centrifuge them, to avoid blockages in the transfer line.

2. Open the Tune window.
3. Click File > Open, and open the default tune file.

Tip: At any time, you can click File > Revert to revert to previously saved parameters.

4. If the instrument is in Source Standby mode, click Operate , and confirm that the adjacent instrument status indicator shows green.
5. For an ETD-enabled instrument, click Acq Mode > CID Mode.
6. Click Source > APPI.

Tip: The Source tab is labeled with the source type and instrument polarity: for example, ES+ means electrospray source and positive ion mode. The Source tab controls source voltages, temperatures, and gas flows.

7. Click positive  ion mode in the toolbar. Negative  ion mode is not used.

Rationale: APPI is not generally useful in negative ion mode.

8. Click MS mode .

Requirement: MS/MS mode configuration uses the same source parameters as the MS mode configuration. In addition, however, you must specify trap and transfer collision energies and supply the set-mass.

9. On the Source tab, set the parameters to these initial values suitable for initial optimization.

Source values:

Parameter	Value
APPI Lamp	On
APPI Repeller	1.00 kV
Sampling Cone	40 V
Extraction Cone	4 V

Temperature values suitable for sample infusion:

Parameter	Value
Source	120 - 140 °C
Probe	600 °C

Gas flows suitable for sample infusion:

Parameter	Value
Cone Gas	0 L/hr
Desolvation Gas	800 L/hr

10. Consult the Instrument tab to ensure the values are similar to the recommended values:


Instrument values:

Parameter	Value
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	MS/MS only

Requirement: Set the collision energies to off, for MS mode, and then, if for MS/MS you need to fragment a sample afterward, turn them on, and set the appropriate voltages.

11. On the Fluidics tab in Sample Flow Control, select your sample vial (Reservoir A).

12. Purge the solvent delivery system for your sample.

- If you changed the vial designation in the previous step then the instrument performs an automatic purge.
- If the vial designation did not change then click Purge .

14. Enter an infusion flow rate of 10 $\mu\text{L}/\text{min}$.

15. Enter a Fill Volume, for example 250 μL .

16. Click Fill .

17. Select a Flow State of Combined.

Tip: APPI requires greater than 200 $\mu\text{L}/\text{min}$ during configuration.

18. Click  to start the flow.

19. Turn on the external (LC) flow.

Requirement In combined LC and sample-infusion mode, the sample is diluted by the ratio of the LC flow rate to infusion flow rate. If the mass of interest is not sufficiently intense for tuning, increase the intensity of the signal by increasing the sample flow or concentration.

20. Set the values in the non-storage acquisition pane to observe the spectrum.

Typical acquisition values:

Parameter	Value
Low mass	100 Da
High mass	1200 Da

Parameter	Value
Scan time	1.0 sec
Data type	Continuum or Centroid
MSMS Mass (only available in MS/MS mode)	Enter a precursor ion mass in Da

21. In the Tune window, zoom the peak displays to a suitable range around the mass of interest (on a spectrum, right-click, drag across the peak's full width, and release).

Tip: In the spectrum, right-click, drag across the peak's full width and release.

! **Notice:** Avoid contact between the APPI source and the APCI probe. The probe is hot and can melt the PEEK insulator.

22. Adjust the position of the APPI lamp to one of the two positions marked on the lamp assembly shaft.

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and more frequent source cleaning can result.

23. Use the probe adjusters to maximize the displayed peak intensity.

Tip: The Vernier adjuster on the mounting flange adjusts the probe tilt. The adjuster-ring on the probe-mount body adjusts the probe z-axis.

24. Adjust the Repeller voltage for maximum sensitivity.

Tip: Starting from a value of 1.0, increase the voltage. Set the voltage to the value giving the highest displayed peak intensity.

25. Optimize the cone voltage for maximum sensitivity:

26. Starting from a value of 15, increase the cone voltage by clicking and dragging the slider. Set the cone voltage to the value giving the highest displayed peak intensity on the sample peak of interest.

27. Confirm that the shape of the displayed peaks are correct.

28. Ensure that resolution, intensity, and stability are acceptable.

29. Click File > Save As and save the settings with an appropriate APPI sample-infusion related name.

Starting with APPI sample-infusion settings, you can optimize the source configuration for the flow-rates needed for your application(s), and you can save them as separate tune files.

Configuring for specific applications


1. Click File > Open and open your sample-infusion tune file.
2. On the Fluidics tab, change the flow state to Combined.
3. Set temperatures and gas flows for your application according to the following table:

Temperature and gas flow settings:

Flow rate (mL/min)	Source temperature (°C)	Desolvation temperature (°C)	Desolvation gas flow (L/h)
0.000 to 0.020	100	200	800

Flow rate (mL/min)	Source temperature (°C)	Desolvation temperature (°C)	Desolvation gas flow (L/h)
0.020 to 0.100	120	350	800
0.101 to 0.300	120	450	800
0.301 to 0.500	150	500	1000
> 0.500	150	600	1200

4. Allow the temperatures to stabilize.

5. Enter a sample Infusion Flow Rate of 10 µL/min and click  to start the flow.

6. Turn on your external (LC) flow.

Requirement: In combined LC and sample-infusion mode, the sample is diluted by the ratio of the LC flow rate to infusion flow rate. If the mass of interest is not sufficiently intense for tuning, the intensity of the signal can be increased as necessary by increasing the sample flow or concentration.

! **Notice:** Take care to avoid contact between the APPI source and the APCI probe. The APCI probe is hot and can melt the PEEK insulator.

7. Adjust the position of the APPI lamp to one of the two positions marked on the lamp assembly shaft.

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and more frequent source cleaning can result.

8. Use the probe adjusters to maximize the displayed peak intensity.

Tip: The Vernier adjuster on the mounting flange adjusts the probe tilt. The adjuster-ring on the probe-mount body adjusts the probe z-axis.

9. Adjust the Repeller voltage for maximum sensitivity:

10. Starting from a value of 1.0, increase the voltage. Set the voltage to the value giving the highest displayed peak intensity.

11. Optimize the cone voltage for maximum intensity.

12. If advised to do so by your applications specialist, adjust the cone gas flow for your sample requirements:

Starting from a value of 0, increase the cone gas flow in increments of 50; allow the flow rate to stabilize after each adjustment. To minimize the formation of solvent ion clusters, set the gas flow to the highest value that does not significantly reduce the peak intensity. For samples in involatile or complex matrix applications, using cone gas can also extend the cleaning interval for your source.

13. Click File > Save As, and rename and save the file.

Once tuned for your application, you can acquire nominal mass data. To obtain exact mass data, configure and use the lock-spray, and calibrate the instrument.

Related Topics

- [Automatic calibration](#)
- [Configure a lock-mass](#)

Configure the ETD mode

-
- [MassLynx](#) > [MS Console](#) > [Instrument: SYNAPT G2-Si](#) > [Intellistart](#)

This procedure describes how to configure the ETD mode. It assumes that the system has optimal system settings.



Notice: The Waters installation engineer leaves your instrument with optimal system settings. We recommend that you do not adjust them. The system settings are password-protected to avoid unintentional changes. Changing them renders your calibrations invalid and reduces the performance of your system.

Note: Before you use the ETD mode, you must ensure that the high voltage connector on the ion block is attached to the glow discharge electrode. For instructions, see SYNAPT G2-Si *Operator's Overview and Maintenance Guide* on the Documentation CD.

When you use the electron transfer dissociation (ETD) fragmentation mode, you must set the conditions for the glow discharge source, which generates the reagent ions, and specify the way in which refill scans are performed. Refill scans are used to introduce the reagent ions into the ETD collision cell, and are performed at selected intervals during tuning or acquisition to ensure a continuous supply of reagent ions.

All the settings defined on the Glow Discharge source tab and in the ETD Acquisitions dialog box are used only for refill scans.

After the initial setup is complete, the settings are saved as a default, sample-dependent tune file. The source parameters can then be adjusted as appropriate for your application. This is achieved by optimizing temperatures and gas flows relative to the LC flow rate, sample type and solvent type. These optimized source conditions can then be saved as a separate tune (.ipr) file, ready for use.


The tune (.ipr) files save values for all supported ionization modes at once, not just the ionization mode selected.

To configure the ETD mode:

1. Click Tune .

2. As a starting point, click File > Open and open the default tune file.

Tip: At any time you can click File > Revert, to revert to previously saved parameters.

3. If the instrument is in Source Standby mode, click Operate  and confirm that the adjacent instrument status indicator shows green.

4. Click Acq Mode > ETD Mode.

5. Click Setup > ETD Acquisition.

6. Enter your required values for refill time, refill interval, and set mass, then click Update.

7. Click Source > Glow Discharge.

Tips:

The source tab is labeled with the source type and instrument polarity, for example, ES+ means electrospray source and positive ion mode. The source tab controls source voltages, temperatures, and gas flows.

When you select the glow discharge source, the instrument is put into negative ion and MSMS mode. You can still select MS and positive ion modes, if required.

8. On the source tab, set the parameters to these initial values:

Source

Parameter	Value
Discharge current	15 μ A
Sampling Cone	0 V
Extraction Cone	4 V

Gas Controls

Parameter	Value
Cone Gas	25 L/hr
Purge Gas	100 L/hr
Make up Gas	10 ml/min

Refill parameters

Parameter	Value
Refill Time	0.1 sec
Refill Interval	1.0 sec
Set Mass	107.10 Da
	Note: You use a set mass of 107.10 with nitrobenzene as the ETD reagent.

Related Topics

- [Automatic calibration](#)
- [Adding lock-mass correction to methods](#)

Configure a SYNAPT G2-Si Mass Spectrometer

- [MassLynx > MS Console > Instruments:SYNAPT G2-Si > IntelliStart > Configure > Configuration mode](#)

Configure the SYNAPT G2-Si for your application. Doing so tunes the source and generates calibrations and LockSpray profiles for use during data acquisition. After completing the configuration, you can prepare the instrument for data acquisition.

To configure the SYNAPT G2-Si:

1. Optimize the source conditions for ionization and beam intensity
2. Define reference compounds for use as calibrants and lock-masses
3. Periodically, determine whether your detector is maintaining the correct gain as it ages
4. Define calibration profiles for each mass range used during acquisition
5. Create calibrations for the defined calibration profiles
6. Define LockSpray profiles for use in providing accurate mass acquisition
7. Set up lock-mass intensity for the defined LockSpray profiles
8. Optionally, adjust the quadrupole RF settings, to optimize the quadrupole mass transmission window

9. Prepare for acquisition.

Related Topics

- [Optimize LockSpray ESI](#)
- [Optimize nanoflow ESI](#)
- [Define reference compounds](#)
- [Check and set up the detector](#)
- [Define calibration profiles](#)
- [Automatic mass calibration](#)
- [Configure a lock-mass for mass accuracy](#)
- [Optimise Quadrapole](#)
- [Prepare the instrument for use](#)

Configure a lockmass

- [MassLynx > MS Console > Instrument: SYNAPT G2-Si > Intellistart](#)

For the LockSpray and nanoLockSpray sources, you can configure lock-mass compounds to use during data acquisition.

Note: In the LockSpray source, you cannot adjust the lock-spray probe position. The lock-spray and sample sprayers share a common gas supply. The nebulizer gas flow is fixed at approximately 75 L/hr by the aperture in the ESI probes.

For the nanoLockSpray source, positional adjustment for the lock-spray probe is achieved using the Vernier adjuster.

The lock-mass sprayer runs from the instrument's integrated solvent delivery system at a flow rate that gives a stable spray, typically in the range 10 to 50 $\mu\text{L}/\text{min}$. To run sample lists overnight or over a weekend, you can use an external lock-spray reservoir. The built-in reservoirs have a capacity of 30 mL giving the following acquisition times as a function of the lock-spray flow rate in use:

- 10 $\mu\text{L}/\text{min}$ gives 50 hours
- 20 $\mu\text{L}/\text{min}$ gives 25 hours
- 30 $\mu\text{L}/\text{min}$ gives 16 hours
- 40 $\mu\text{L}/\text{min}$ gives 12 hours
- 50 $\mu\text{L}/\text{min}$ gives 10 hours

Choose a concentration of lock-mass solution that gives a suitable ion intensity.

Setting up long-term mass accuracy during your sample runs

To set up long-term mass accuracy during your sample runs:

1. [Configure MassLynx for your compound](#)
2. [Check the lock-mass signal intensity](#)
3. [Perform LockSpray source setup.](#)

If a lock-mass entry does not already exist for your chosen reference compound, you must create one.

Prerequisite: To configure your lock-mass, you must know the relative intensities of the mass peaks that you want to use for your reference compound. Before you start, run the lock-mass compound, and view the spectrum in the Tune window. Make sure the peaks you choose are spread across the spectrum and that the highest mass you choose is at least double the lowest mass.

To configure MassLynx for your compound:

1. From IntelliStart, click Configure > Configuration Mode, and then Configure > Lock Mass Editor.
2. In the editor window, click File > New.
3. In the New LockMass dialog, enter a unique LockMass Name (usually based on the name of your lock-mass compound).
4. Choose to create a lock-mass based on a reference file and select a compound from the drop-down list.

Alternative: Specify your own compound.

Tip: If you choose to specify your own compound, to add a line to the Mass table, right-click, click Add and enter a lock-mass value. To remove items, right-click on the item and select Remove.

5. Select MS mode or MSMS mode.

Requirement: The lock-mass mode must match the mode you intend using for your acquisition.

6. Against either Positive Polarity or Negative Polarity, select collision energy, On or Off, and set a suitable voltage.

For your reference compound, the lock-mass polarity you set up is optional. You can supply either positive only, negative only, or both positive and negative lock-mass settings. Repeat these instructions for each polarity.

7. Select two or more mass values to use for the multipoint lock-mass calculations by placing a check in the Include column.


Rule: If you are defining lock-mass from an existing reference file, the intensity range of your selected peaks must not be more than 4:1. For example, select peaks between 100% and 25% intensity or between 40% and 10% intensity.

Requirement: The highest mass must be at least double the lowest mass.


Tip: If you are defining your own masses, then no intensity information is required.

8. Click OK to save the lock-mass file.

Checking the lock-mass signal intensity

1. From IntelliStart, click Configure > Normal Mode.
2. Select Lock Mass Check, click Start , and follow the wizard.

Performing LockSpray source setup

1. From IntelliStart, click Configure > Configuration mode.
2. Select LockSpray Source Setup, click , and follow the wizard.

Related Topics

- [Automatic calibration](#)

- [Adding lock-mass correction to methods](#)

Configure a lockspray with an APCI probe

- [MassLynx > MS Console](#)

This procedure describes how to configure the LockSpray source and APCI probe. It assumes that the system has optimal system settings.

! **Notice:** The Waters installation engineer leaves your instrument with optimal system settings. Do not adjust them. The system settings are password-protected to avoid unintentional changes. Changing them renders your calibrations invalid and reduces the performance of your system.

Note: If you have an ETD-enabled ion block fitted, you must ensure that the discharge electrode cable is connected to the ion block connector and that the discharge needle is connected. For instructions, see *SYNAPT G2-Si Operator's Overview and Maintenance Guide* on the Documentation CD.

This procedure comprises the following two processes:

- [Setting up the beam initially by combining sample infusion](#)
- [Configuring specific applications](#)

Sample infusion is a low flow-rate infusion of a typical sample from the vial mounted on the front of the SYNAPT G2-Si. It allows initial source configuration before switching to the high-flow conditions generated by your LC system.

Following the initial setup, save the settings as a default, sample-dependent tune file. You can adjust the source parameters as appropriate for your application. To do so, you optimize temperatures and gas flows relative to the LC flow rate, sample type and solvent type. You can save these optimized source conditions as a separate tune (.ipr) file, ready for use.

The tune (.ipr) files save values for all supported ionization modes at once, not just the ionization mode selected.

For APCI, supplement the low flow-rate sample by specifying an LC flow that permits appropriate, high-flow, configuration settings. The sample must be typical for your application and is taken from the vial mounted on the front of the instrument.


To set up the beam initially by combining sample infusion:

1. Install suitable sample, reference and calibrant vials in positions A, B and C on the front of your SYNAPT G2-Si. For installation instructions, see *SYNAPT G2-Si Operator's Overview and Maintenance Guide* on the Documentation CD.
2. **Recommendation:** Use vials in this order:
 - A, for sample
 - B, for reference
 - C, for calibrant




Reference and calibrant solutions are described in their respective help topics. Ensure that sample solutions are based on your eventual application and that they have been filtered or spun in a centrifuge to avoid blockages in the transfer line.

3. Open the Tune window.
4. Click File > Open and open the default tune file.

Tip: At any time you can click File > Revert, to revert to previously saved parameters.

- If the instrument is in Source Standby mode, click Operate , and confirm that the adjacent instrument status indicator shows green.
- If you have an ETD-enabled instrument, click Acq Mode > CID Mode.
- Click Source > APCI.

Tip: The source tab controls source voltages, temperatures, and gas flows. It is labeled with the source type and instrument polarity, for example, ES+ means electrospray source and positive ion mode.

- Click positive  or negative  ion mode in the toolbar.
- Click MS mode .

Tip: MS/MS configuration uses the same source parameters as MS but in addition, trap and transfer collision energies, and set-mass must be supplied.

- On the source tab, set the parameters to the following initial values.

Source values:

Parameter	Value
Current mode	On
Corona current	3.0 μ A
Sampling Cone	40 V
Extraction Cone	4 V

Temperature values suitable for sample infusion:

Parameter	Value
Source	120 - 140 °C
Probe	600 °C

Gas flows suitable for sample infusion:


Parameter	Value
Cone Gas	0 L/hr
Desolvation Gas	800 L/hr

- Consult the Instrument tab to ensure the values are similar to the recommended values:

Instrument values:

Parameter	Value
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	MS/MS only

Tip: Set the collision energies to off, for MS mode, and then, if for MS/MS you need to fragment a sample afterward, turn them on, and set the appropriate voltages.

12. On the Fluidics tab in Sample Flow Control, select your sample vial (Reservoir A).
13. Purge the solvent delivery system for your sample.
 - If you changed the vial designation in the previous step then the instrument performs an automatic purge.
 - If the vial designation did not change then click Purge .

14. Enter an infusion flow rate of 10 $\mu\text{L}/\text{min}$.

15. Enter a Fill Volume, for example 250 μL .

16. Click Fill .

17. Select a Flow State of Combined.

Tip: APCI requires greater than 200 $\mu\text{L}/\text{min}$ during configuration.

18. Click  to start the flow.

19. Turn on the external (LC) flow.

Tip: In combined LC and sample-infusion mode, the sample is diluted by the ratio of the LC flow rate to infusion flow rate. If the mass of interest is not sufficiently intense for tuning, increase the intensity of the signal by increasing the sample flow or concentration.

20. Set the values in the non-storage acquisition pane to observe the spectrum.

Typical acquisition values:

Parameter	Value
Low mass	100 Da
High mass	1200 Da
Scan time	1.0 sec
Data type	Continuum or Centroid
MSMS Mass (only available in MS/MS mode)	Enter a precursor ion mass in Da

21. In the Tune window, zoom the peak displays to a suitable range around the mass of interest (on a spectrum, right-click, drag across the peak's full width, and release).

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and more frequent source cleaning can result.

22. Use the Vernier probe adjuster on the probe's mounting flange to maximize the peak intensity.

23. In Current mode, adjust the corona current to the minimum value required for acceptable intensity.

Rationale: Using the minimum value possible for the desired effect prolongs the life of your corona pin.

24. Optimize the cone voltage for maximum sensitivity:


25. Starting from a value of 15, increase the cone voltage by clicking and dragging the slider. Set the cone voltage to the value giving the highest displayed peak intensity on the sample peak of interest.
26. Confirm that the shape of the displayed peaks are correct.
27. Ensure that resolution, intensity, and stability are acceptable.
28. Click File > Save As and save the settings with an appropriate APCI sample-infusion related name.
29. Starting with APCI sample-infusion settings, you can optimize the source configuration for the flow-rates needed for your application(s) and save them as separate tune files.

To configure for specific applications:

1. Click File > Open and open your sample-infusion tune file.
2. On the Fluidics tab, in Sample Flow Control, change the flow state to Combined.
3. Set temperatures and gas flows for your application according to the following table.

Temperature and gas flow settings:

Flow rate (mL/min)	Source temperature (°C)	Desolvation temperature (°C)	Desolvation gas flow (L/h)
0.000 to 0.020	100	200	800
0.020 to 0.100	120	350	800
0.101 to 0.300	120	450	800
0.301 to 0.500	150	500	1000
> 0.500	150	600	1200

4. Allow the temperatures to stabilize.
5. Enter a sample Infusion Flow Rate of 10 µL/min and click  to start the flow.
6. Turn on your external (LC) flow.

Tip: In combined LC and sample-infusion mode, the sample is diluted by the ratio of the LC flow rate to infusion flow rate. If the mass of interest is not sufficiently intense for tuning, the intensity of the signal can be increased as necessary by increasing the sample flow or concentration.

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and more frequent source cleaning can result.

7. Use the Vernier probe adjuster on the probe's mounting flange to maximize the peak intensity.
8. In Current mode, adjust the corona current to the minimum value required for acceptable intensity.

Rationale: Using the minimum value possible for the desired effect prolongs the life of the corona pin.
9. Optimize the cone voltage for maximum intensity.
10. If advised to do so by your applications specialist, adjust the cone gas flow for your sample requirements.
11. Starting from a value of 0, increase the cone gas flow in increments of 50; allow the flow rate to stabilize after each adjustment. To minimize the formation of solvent ion clusters, set the gas flow to the highest value that

does not significantly reduce the peak intensity. For involatile or complex matrix applications, using cone gas can also extend the cleaning interval for your source.

12. Click File > Save As and save the file with a new name.

Once tuned for your application, you can acquire nominal mass data. To obtain exact mass data, configure and use the lock-spray and calibrate the instrument.

Related Topics

- [Automatic calibration](#)
- [Configure a lock-mass](#)

Configure the lockspray for ESCi operation

- [MassLynx > MS Console](#)

This procedure explains how to configure the LockSpray source for ESCi operations. It assumes optimally set system parameters.

! **Notice:** The Waters installation engineer leaves your instrument with optimal system settings. Do not adjust them. The system settings are password-protected to avoid unintentional changes. Changing them renders your calibrations invalid and reduces the performance of your system.

Requirement: If you have an ETD-enabled ion block fitted, ensure that the discharge electrode cable is connected to the ion block connector. For instructions, see ETD Source for SYNAPT G2-Si *Operator's Overview and Maintenance Guide* on the Documentation CD.

Rule: It is possible to select ESCi in ETD mode, but this option is not applicable for ETD experiments.

This procedure comprises the following two processes:

- Setting up the beam initially using sample infusion
- Configuring specific applications

Following the initial setup, save the settings as a default, sample-dependent tune file. You can then adjust the source parameters as appropriate for your application. To do so, you optimize temperatures and gas flows relative to the LC flow rate, sample type, and solvent type. You can save the optimized source conditions as a separate tune (.ipr) file, ready for use.

The tune (.ipr) files save values for all supported ionization modes at once, not only the ionization mode selected.

Setting up the beam initially using sample infusion

Sample infusion is a low flow-rate infusion of a typical sample from the vial mounted on the front of the SYNAPT G2-Si. It allows initial source configuring before switching to the high-flow conditions generated by the LC system.

To set up the beam initially using sample infusion:

1. Install suitable sample, reference, and calibrant vials in positions A, B, and C on the front of the SYNAPT G2-Si.

Tip: For installation instructions, see *SYNAPT G2-Si Operator's Overview and Maintenance Guide* on the Documentation CD.

Recommendation: Use reservoir bottles in the following order:


- A, for sample

- B, for reference
- C, for calibrant




Tip: Reference and calibrant solutions are described in their respective help topics. Base sample solutions on your eventual application. Filter or centrifuge them, to avoid blockages in the transfer line.

2. Open the Tune window.
3. As a starting point, click File > Open, and open the default tune file.

Tip: At any time you can click File > Revert, to revert to previously saved parameters.

4. If the instrument is in Source Standby mode, click Operate , and confirm that the adjacent instrument status indicator shows green.
5. For an ETD-enabled instrument, click Acq Mode > CID Mode.
6. Click Source > ESCi.

Tip: The Source tab is labeled with the source type and instrument polarity; for example, ES+ means electrospray source and positive ion mode. The Source tab controls source voltages, temperatures, and gas flows.

7. Click positive  or negative  ion mode in the toolbar.
8. Click MS mode .

Requirement: MS/MS mode configuration uses the same source parameters as the MS mode configuration. Nevertheless, you must specify trap and transfer collision energies and supply the set-mass.

9. On the Source tab, set the parameters to these initial values, which are suitable for optimization by infusion.

Source values:

Parameter	Value
Capillary	3.00 kV
Corona current	2.0 μ A
Sampling Cone	40 V
Extraction Cone	4 V

Temperature values suitable for sample infusion:

Parameter	Value
Source	100 $^{\circ}$ C
Desolvation	200 $^{\circ}$ C

Gas flows suitable for sample infusion:




Parameter	Value
Cone Gas	0 L/hr
Desolvation Gas	800 L/hr

10. Consult the Instrument tab to ensure the values are approximate to the recommended values.

Instrument values:

Parameter	Value
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	MS/MS only

Requirement: Set the collision energies to off, for MS mode, and then, if for MS/MS you must fragment a sample afterward, turn them on, and set the appropriate voltages.

11. On the Fluidics tab, select the sample vial (Reservoir A).
12. Purge the solvent delivery system for your sample:
 - If you changed the vial designation in the previous step, the instrument performs an automatic purge.
 - If the vial designation did not change, click Purge .
13. Enter an infusion flow rate of 10 $\mu\text{L}/\text{min}$.
14. Enter a Fill Volume, for example 250 μL .
15. Click Fill .
16. Select a flow state of infusion.
17. Click  to start the flow.
18. Set the values in the non-storage acquisition pane to observe the spectra.

Typical acquisition values:

Parameter	Value
Low Mass	100 Da
High Mass	1200 Da
Scan Time	1.0 sec
Data Type	Continuum or Centroid
MSMS Mass (only available in MS/MS mode)	Enter a precursor ion mass in Da

19. In the Tune window, zoom the peak displays to a suitable range around the mass of interest (on a spectrum, right-click, drag across the peak's full width, and release).
20. Use the Vernier probe adjuster on the probe's mounting flange to maximize the peak intensity.

To ascertain its position, set the capillary voltage to zero. If the signal intensity on the ESI spectrum fails to drop to less than 10% of that observed with optimum capillary voltage, move the probe away from the sample cone slightly. Return the capillary voltage to its previous setting, and observe the signal intensity. Repeat as necessary.

Requirement: This probe position is appropriate for a 10 $\mu\text{L}/\text{min}$ flow rate. You must adjust it later, however, for higher flow rates (see later).

21. For ESI, adjust the Capillary voltage for maximum sensitivity.

Tip: Most compounds are optimized at between 1.0 and 3.0 kV.

22. For APCI operation, in Current mode, adjust the corona current to the minimum value required for acceptable intensity.

Rationale: Using the minimum value possible for the desired effect prolongs corona pin life.

23. Optimize the capillary protrusion using the knurled knob atop of the ESI probe.

24. Optimize the cone voltage for maximum sensitivity:

Starting from a value of 15, increase the cone voltage by clicking and dragging the slider. Set the cone voltage to the value giving the highest displayed peak intensity on the sample peak of interest.

25. Confirm that the shape of the displayed peaks are correct.

26. Ensure that resolution, intensity, and stability are acceptable.

27. Click File > Save As and save the settings with an appropriate sample-infusion related name.

Configuring for specific applications

Starting with sample-infusion settings, you can optimize the source configuration for the higher flow-rates needed for your application(s), and save them as separate tune files.


To configure for specific applications:

1. Click File > Open and open your sample-infusion tune file.
2. On the Fluidics tab, change the flow state to Combined.
3. Set temperatures and gas flows for your application according to the following table.

Settings required for specific flow rates:

Flow rate (mL/min)	Source temperature (°C)	Desolvation temperature (°C)	Desolvation gas flow (L/h)
0.000 to 0.020	100	200	800
0.020 to 0.100	120	350	800
0.101 to 0.300	120	450	800
0.301 to 0.500	150	500	1000
> 0.500	150	600	1200

4. Allow the temperatures to stabilize.

5. Enter a sample infusion flow rate of 10 μ L/min, and click  to start the flow.

6. Start the external (LC) flow.

Requirement In combined LC and sample-infusion mode, the sample is diluted by the ratio of the LC flow rate to infusion flow rate. If the mass of interest is not sufficiently intense for tuning, increase intensity of the signal as necessary, by increasing the sample flow or concentration.

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and more frequent source cleaning can result

7. Use the Vernier probe adjuster on the probe's mounting flange to maximize the peak intensity.

To ascertain its position, set the capillary voltage to zero. If the signal intensity fails to drop to less than 10% of that observed with optimum capillary voltage, move the probe away from the sample cone slightly. Return the capillary voltage to its previous setting and observe the signal intensity. Repeat as necessary.

8. For ESI, adjust the Capillary voltage for maximum sensitivity.

9. For APCI, in Current mode, adjust the corona current to the minimum value required for acceptable intensity.

Rationale: Using the minimum value possible for the desired effect prolongs the life of the corona pin.

10. Optimize the capillary protrusion using the knurled knob on the top of the ESI probe.

11. Optimize the cone voltage for maximum intensity.

12. If advised to do so by your applications specialist, adjust the cone gas flow for your sample requirements.

Starting from a value of 0, increase the cone gas flow in increments of 50; allow the flow rate to stabilize after each adjustment. To minimize the formation of solvent ion clusters, set the gas flow to the highest value that does not significantly reduce the peak intensity. For involatile or complex matrix applications, using cone gas can also extend the cleaning interval for your source.

13. Click File > Save As, and rename and save the file.

Once tuned for your application you can acquire nominal mass data. To obtain exact mass data, configure and use the lock-spray, and calibrate the instrument.

Related Topics

- [Automatic calibration](#)
- [Configure a lock-mass](#)

Define calibration profiles

- [MassLynx > MS Console > Instruments:SYNAPT G2-Si > Intellistart > Configure > Configure mode](#)

Calibration profiles store a mass range, polarity, mode (sensitivity or resolution), and reference compound. Once defined, calibration profiles are used to calibrate the instrument for the conditions present at the time of the calibration.

The pusher frequency is a significant instrument condition that changes according to the acquisition mass range. The pusher is the metal plate that accelerates the ions into the flight tube of the TOF region of the instrument. It is connected to a switched, high-voltage, power supply and generates an impulse electrical field across the ion beam that fires the ions into the flight tube.

For lower mass ranges, the pusher accelerates groups of ions thousands of times per second. For higher mass ranges, to compensate for longer flight times, the pusher frequency is automatically set to a lower value as it exceeds each of the mass boundaries listed below:

- 1200 Da
- 2000 Da
- 5000 Da

- 8000 Da
- 14000 Da
- 32000 to 100000 Da

Tip: You must create calibration profiles that are compatible with your desired acquisition mass range. Compatible mass ranges use the same pusher frequency. Where the calibration is not compatible, subsequent acquisitions are not calibrated correctly.

For example, when the range of the active calibration is between 50 and 1300 Da, the calibration passes above the 1200-Da mass boundary and uses the lower pusher frequency. An acquisition from 500 to 1200 Da using this calibration does not pass above the 1200 Da boundary. Thus it uses a higher pusher frequency, so the results are not correctly calibrated.

Once created, calibrations based on the calibration profile are stored on the instrument.

To add a calibration profile:

1. Click Configure > Calibration Profile Editor.
2. In the Calibration Profile Editor, click File > New.
3. Enter a name for the calibration profile name and, optionally, a description.
4. Enter a value for the start mass of the mass range.
5. For the end mass, select a standard mass, or enter a custom value.

Tip: The standard mass values take account of the pusher frequency boundaries.

6. Select Automatic or Assisted.

Tip: An assisted calibration includes the opportunity to select and deselect matching peaks in the Calibration dialog box. Automatic does not.

7. Enable the polarities, and select reference compounds from the lists.
8. Click OK, to create the calibration profile.

Related Topics

- [Prepare the instrument for use](#)
- [Define reference compounds](#)
- [Automatic mass calibration](#)

Define reference compound use

- [MassLynx > MS Console > Instruments: SYNAPT G2-Si > Intellistart > Configure > Configuration mode](#)

Use the reference compound editor to define compounds for use as calibrants or lock-masses.

Tips:

- Set up calibrants that are suitable for the mass ranges over which acquisitions run.
- Set up LockSpray compounds that produce peaks you can use as references at the high and low mass ends of the range.

A list of reference compounds appears in a table along with their ion mode and a list of their mass peaks. For these ionized compounds, each mass peak represents the ion peak position on the mass axis of the acquired spectra equivalent to the ionic mass-to-charge ratio.

Select a row to display further details, including the mass spectrum. To make changes, right-click, and click an option.

To add a reference compound:

1. Click Configure > Reference File Editor.
2. In the editor, click File > New.
3. Enter a unique name, like the compound's name.
4. Select the polarity required.
5. Select MS mode or MSMS mode.

Tip: If you select MSMS and specify a collision energy, the masses you specify relate to collision-induced fragments of the compound.

Requirement: For lock-mass compounds, the mode selected must match the mode used for acquisition.

6. Select ToF or Mobility mode.
7. Provide further details in the description.
Tip: Such details can include the intended use, compound supplier details, and order codes.
8. Enter the cone and capillary voltages necessary for the ionization.
9. For MSMS mode, enter a set mass, enable the collision energy, and set a collision voltage.
10. Add a row for each expected monoisotopic mass peak.

Tip: The monoisotopic mass is the sum of the masses of the atoms in a molecule calculated using the principle (most abundant) isotope for each element, not the average mass of each element. For typical organic compounds, this calculation results in selection of the lightest isotopes. Nevertheless, for some heavier atoms, such as iron and argon, the principle isotope is not the lightest isotope.

Requirement: The SYNAPT G2-Si provides mass accuracies for small molecules for which the mass of the charge carrier is significant. For the monoisotopic mass peaks you provide, be sure to consider the mass of the charge carriers (including electrons).

11. Click Save.

Related Topics

- [Prepare the instrument for use](#)
- [Configure a lock-mass](#)
- [Define LockSpray profiles](#)

Manual calibration

- [Console System Tree > Instruments: SYNAPT G2-Si > Intellistart](#)

Most situations call for automatic calibration, See " Calibration (Automatic)" . This manual method is provided for users who prefer to use it, or prefer to use as an assisted method when fault finding or experimenting.

Note: Following a manual calibration, MassLynx software does not inform you when the acquisition mass range falls outside the calibration mass range. For APGC operation, you must perform a manual calibration.



The SYNAPT G2-Si mass scale gives nominal mass measurement without a calibration. However, it must be calibrated before using for exact-mass measurement. MassLynx software uses a polynomial equation to calibrate the system precisely over a wide mass range. Calibrating and checking calibrations always produces a report (stored in C:\MassLynx\IntelliStart\Results\Reports). You can optionally select Display report during the wizard to read the report when the calibration or calibration check finishes.

Tip: You should first check any existing calibration to see if a calibration procedure is necessary.

To calibrate manually:

1. Specify the conditions necessary for acquiring data from the reference compound.


Tip: For the APGC, place a reservoir bottle of heptacosane in one of the reservoir bottle carriers. You can also lower the temperature of the transfer line to 150 °C, to reduce interference.

2. From IntelliStart, click Tune , and select Acquire > Start.
3. Acquire 2 minutes of TOF MS continuum data, with a normal dynamic range for your compound, across the desired mass range.
4. Process the data ready for calibrating:
 - From the MassLynx main window, click Chromatogram.
 - Click File > Open, and open the recently acquired chromatogram.
 - Click Display > Remove, and remove unwanted chromatograms from the window.
 - Right-click at the zero end of the trace and drag across the chromatogram to combine 2 minutes of spectra. Then right click again to display a marker and the spectrum window
 - In the spectrum window, click Process > Automatic Peak Detection, and click OK to create a centroid spectrum.
 - Click Display > Remove, and remove the continuum spectrum.
 - From the File menu, click Save Spectrum.
5. Return to IntelliStart and click Configure > Configuration Mode.
6. Click Create Calibration and then Start .
7. On the Create Calibration dialog box, click Next.
8. Click Calibration Profile Editor.
9. On the Calibration Profile Editor, click File > New.
10. On the New Calibration Profile dialog box, provide a name for the calibration.
11. For type of calibration, select Manual.
12. Below the required polarity, select the mode in which the data was acquired and then click Edit.
13. Enable the calibration polarity, and select a reference compound from the drop-down list.
14. Browse to the raw data file.

- In the Select Raw File dialog box, select the file, and click History.
- Select the saved, centroid spectrum.

15. Click OK.

Result: The Calibration dialog appears.

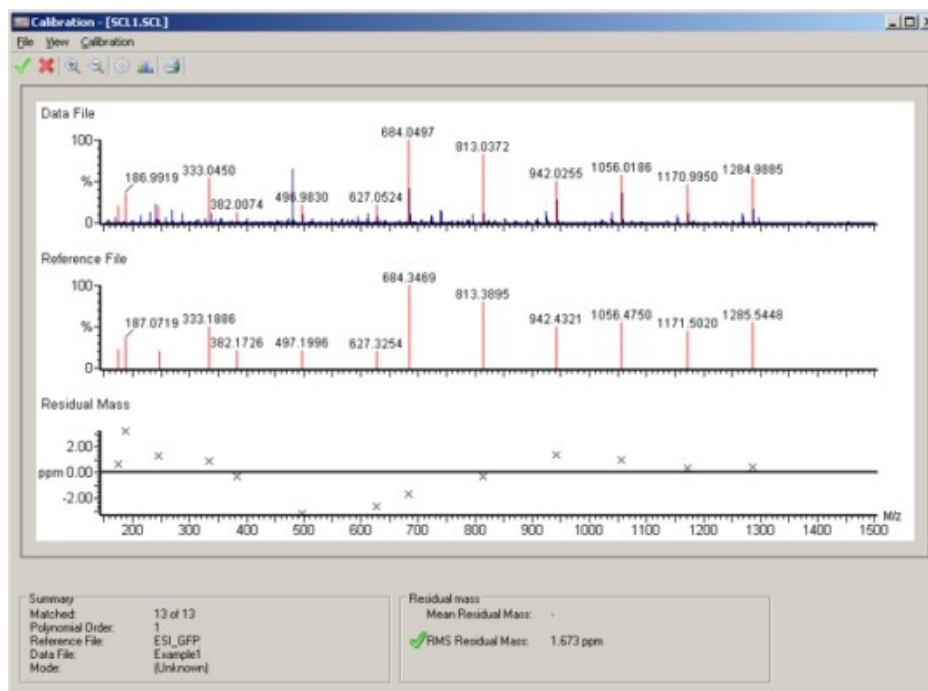

Tips:

- Use the mouse wheel to zoom the data.
- For an alternative view of the residuals, click View > Display Residuals as PPM or View > Display Residuals as AMU.

16. Click Calibration > Quality of fit, and set a residual mass.

Tip: Set a threshold, for example 10 ppm, then match peaks until the Residual Mass readout below the data display changes to a green tick symbol.

17. Check for outliers and remove them from the calibration curve by right-clicking the outlying peak in the reference file spectrum and right-clicking the associated peak on the data file spectrum.
18. Try alternative polynomial orders in Calibration > Parameters for a better fit of calibration curve.

Requirement: After changing the polynomial order remove your outliers again.

19. Accept or reject the calibration.

- If satisfied with the results, click File > Accept Calibration.
- If dissatisfied with the results, click File > Reject Calibration.

Result: Once accepted here, the calibration becomes the active calibration on the instrument.

For positive polarity calibrations, where you selected Both, the process continues for Negative.

Related Topics

- [Automatic calibration](#)

Optimize lockspray for_ESI

- [MassLynx > MS Console > Instruments: SYNAPT G2-Si > MS Tune](#)

Perform the following tasks, in sequence, to optimize the electrospray ionization and ESI source for your application:

The Waters installation engineer leaves your instrument with optimal system settings. Waters recommends you leave them unchanged. The system settings are password-protected, to avoid unintentional changes. Changing the settings renders your calibrations invalid and hinders the performance of your system.

- [Set up the beam initially using sample infusion](#)
- [Fine-tune the beam for your specific application](#)

After completing the initial setup, which involves tuning the beam by infusion, save the settings as a default, sample-infusion tune file. You can then adjust the source parameters, as appropriate for your application, by optimizing temperatures and gas flows relative to your LC flow rate, sample type, and solvent type. You can save these optimized source conditions as a separate tune (.ipr) file, ready for use.

The tune (.ipr) files save values for all supported ionization modes at once, not just the ionization mode selected.

Sample infusion is a low-flow-rate infusion of a typical sample from the reservoir bottle mounted on the front of the SYNAPT G2-Si. It allows initial source configuration before switching to the high-flow conditions generated by your LC system.

To set up the beam initially, using sample infusion:

1. Install suitable sample, reference (lock-mass), and calibrant reservoir bottles in positions A, B, and C on the front of your SYNAPT G2-Si. See the *Operator's Overview and Maintenance Guide* for installation instructions.


Recommendation: Use reservoir bottles in this order:

- A, for sample
- B, for reference
- C, for calibrant.


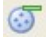

Reference and calibrant solutions are described in their respective help topics. Base sample solutions on your eventual application. Filter or centrifuge them, to avoid blocking the transfer line.

1. Open the Tune window , click File > Open, and open the default tune file.

Tip: At any time you can click File > Revert, to revert to previously saved parameters.

2. If the instrument is in source standby mode, click Operate  and confirm that the instrument's status indicator shows green.
3. From the Source menu, click ESI.

Tip: The Source tab is labeled with the source type and instrument polarity. For example, ES+ means electrospray source and positive ion mode. The source tab controls source voltages, temperatures, and gas flows.

- Click positive  or negative  ion mode in the toolbar.
- Click MS mode .

Tip: MS/MS configuration uses the same source parameter settings as those used for the MS configuration. However, you must in addition provide collision energy and a set-mass.

- On the Source tab, set the parameters to these initial values shown in the table:

Source values:

Parameter	Value
Capillary	3.00 kV
Sample cone	40 V
Extractor cone	4 V

Temperature values suitable for sample infusion:

Parameter	Value
Source temperature	100 °C
Desolvation temperature	200 °C

Gas flows suitable for sample infusion:

Parameter	Value
Cone gas flow	0 L/hr
Desolvation gas flow	800 L/hr

Refill:

Note: The Refill parameter is only available if you have chosen ETD as the acquisition mode.

- Consult the Instrument tab to ensure the values reflected there approximate the recommended instrument values:


Instrument tab values:

Parameter	Value
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	MS/MS only

Tip: For MS mode, set the collision energy to off. If, for MS/MS mode, you must fragment a sample, set the collision energy to on, and specify an appropriate voltage.

- On the Fluidics tab, select the sample reservoir bottle (Reservoir A) in Sample Flow Control.
- Select Infusion as the flow path.
- Enter an infusion flow rate of 10 µL/min.

5. Select a fill volume: for example, 250 μL .

6. Click Purge , to wash the pump, and then refill the pump with sample solution.

Tip: Each wash cycle represents four refills of the pump with wash solution.

7. Click  to start the flow.

8. Set the appropriate values to observe the spectrum:

Typical acquisition values:

Parameter	Value
Low Mass	100 Da
High Mass	1200 Da
Scan Time	1.0 sec
Data Type	Continuum or Centroid
MSMS Mass (only available in MS/MS mode)	Enter a precursor ion mass in Da

9. In the plot control pane, click the drop-down, click Mode Settings, set the mode as TOF mode, and the scanning mode as MS or MS/MS.

Rationale: This mode is used in combination with the view settings to determine how the data is presented in the trace display. Quad mode allows displays of Intensity (chromatographic view), Spectrum view, Resolution view, Signal-to-Noise view. TOF mode allows a Spectrum view, Spectrum + Multi, Spectrum + Intensity, Spectrum + Resolution, and Spectrum + Mass Position.

10. Zoom the peak display to a suitable range around the mass of interest.

Tip: On the spectrum, click and drag across the desired range, and then release.

11. Use the vernier probe adjuster on the probe's mounting flange to maximize the peak intensity.

Tips:

- Some compounds can give a high signal, the result of ionization other than electrospray, when the probe is too close to the sample cone or pointing directly at the aperture. Such a signal is undesirable.
- To test for excessive non-electrospray ionization, set the capillary to zero. If the signal intensity fails to drop to less than 10% of that observed with optimum capillary voltage, move the probe away from the sample cone slightly. Return the capillary voltage to its previous setting, and observe the signal intensity. Repeat as necessary.
- This probe position is appropriate for 10 $\mu\text{L}/\text{min}$. You must adjust it later for higher flow rates.

12. Adjust the capillary voltage for maximum sensitivity.

Tip: Most compounds are optimized at between 1.0 and 3.5 kV.

13. Optimize the capillary protrusion for maximum intensity using the knurled knob atop the ESI probe.

14. Optimize the sample cone voltage, for maximum sensitivity:

Tip: Starting from a value of 15, increase the cone voltage by clicking and dragging the slider. Set the cone voltage to the value giving the highest displayed peak intensity on the sample peak of interest.

15. Ensure that resolution, intensity, and stability are acceptable (see "Define calibration profiles").
16. Click File > Save As and save the Tune file (.ipr) with a new name.

Starting with sample-infusion settings, you can optimize the source configuration for the higher flow-rates needed for your specific application.


[Back to top](#)

Fine-tune the beam for your specific application:

1. Click File > Open, and open the sample-infusion tune file.
2. Click the Fluidics tab.
3. Change the sample flow path to Combined.
4. Set temperatures and gas flows for your application.

Recommended temperatures by flow rate:

Flow rate (mL/min)	Source temperature (°C)	Desolvation temperature (°C)	Desolvation gas flow (L/h)
0.000 to 0.020	100	200	800
0.020 to 0.100	120	350	800
0.101 to 0.300	120	450	800
0.301 to 0.500	150	500	1000
> 0.500	150	600	1200

5. Allow the temperatures to stabilize.
6. Enter a sample infusion flow rate of 10 µL/min, and start the flow .
7. Start the external (LC) flow.

Tip: In combined (LC and sample-infusion mode), the sample is diluted by the ratio of the LC flow rate to the infusion flow rate. If the mass of interest is not sufficiently intense for tuning, increase the intensity of the signal, as necessary, by increasing the sample flow or concentration.

! **Notice:** If you position the probe too close to the sample cone, nonlinear data, sample suppression, and the need for frequent source cleaning can result.

8. Maximize the peak intensity by turning the vernier probe adjuster on the probe's mounting flange.

Tip: To ascertain the probe's position, set the capillary voltage to zero. If the signal intensity fails to drop to less than 10% of that observed with optimum capillary voltage, move the probe away from the sample cone slightly. Return the capillary voltage to its previous setting, and observe the signal intensity. Repeat as necessary.

9. Adjust the capillary voltage for maximum sensitivity.
10. Optimize the capillary protrusion using the knurled knob atop the ESI probe.
11. Optimize the sample cone voltage for maximum intensity.
12. If advised to do so by your applications specialist, adjust the cone gas flow for your sample requirements.

Tip: Starting from a value of 0, increase the cone gas flow in increments of 50; allow the flow rate to stabilize after each adjustment. To minimize the formation of solvent ion clusters, set the gas flow to the highest value that does not significantly reduce the peak intensity. For samples in involatile or complex matrix applications, using cone gas can also extend the source's cleaning interval.

13. Click File > Save As, and save the file with a new name.

[Back to top](#)

Having tuned the instrument for your application, you can acquire nominal mass data. To obtain exact mass data, calibrate the instrument, and then configure and use the LockSpray.

Related Topics

- [Prepare the instrument for use](#)
- [Automatic calibration](#)
- [Configure a lock-mass](#)

Optimise nanospray for ESI

- [MassLynx > MS Console > Instruments: SYNAPT G2-Si > MS Tune](#)

Perform the tasks presented in this topic to configure the nanoLockSpray source for ESI operation with the Universal Sprayer. See also the *SYNAPT G2-Si Operator's Overview and Maintenance Guide*, for information about selecting, preparing, and positioning the sample sprayer.

The Waters installation engineer leaves your instrument with optimal system settings. Waters recommends you leave them unchanged. The system settings are password-protected, to avoid unintentional changes. Changing the settings renders your calibrations invalid and hinders the performance of your system.


We recommend that you optimize the source in steps, as follows:

- [Enable the video camera](#)
- [Set up the beam initially, using sample infusion](#)
- [Fine-tune the beam for your specific application](#)

After completing the initial setup, which covers tuning the beam, save the settings as a default NanoFlow tune file. You can then adjust the source parameters, as appropriate for your application, by optimizing temperatures and gas flows relative to your flow rate, sample type, and solvent type. You can save these optimized source conditions as a separate tune (.ipr) file, ready for use.

The tune (.ipr) files save values for all supported ionization modes at once, not just the ionization mode selected.


The sample sprayer receives gas from a mass-flow controller mounted beneath the stage on the front of the nanoLockSpray source to which it is connected by a 1/16-inch PTFE tube. Positional adjustment for the sprayer is achieved using the x-, y-, z-axis position adjusters. A camera monitors the position of the probe and spray.

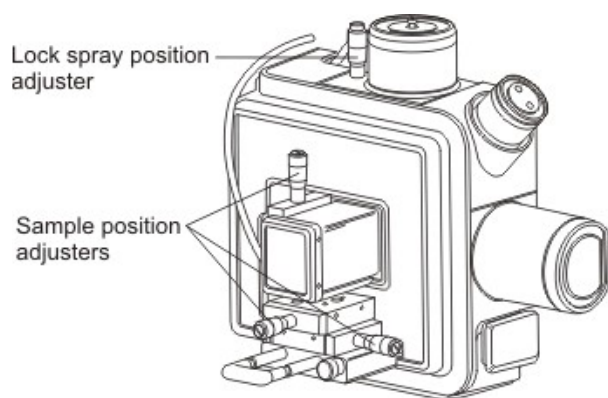

To enable the video camera:

Requirement: The nanoLockSpray source must be fitted to the instrument for this option to be available.

Tip: When a new source is fitted, an automatic source recognition dialog is displayed. Confirm the source type, and click OK.

1. Click Camera, .
2. Manipulate the camera mount while viewing the video in the Camera Control dialog box.
Tip: To aim, move the camera mount. To focus, rotate the mount.
3. Position the camera so that the tip of the sample sprayer and the cone aperture are in view.


[Back to top](#)

To set up the beam initially, using sample infusion:

1. Temporarily connect 75 μm , fused-silica tubing between the grounded union and the universal sprayer (see the *Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide*).

Rationale: For tuning, reconfigure the nanoflow solvent delivery to introduce your sample spray using the reference fluidics plumbing. Tune the source using the reference fluidics pump. This pump performs better at the low flow rates required for nanoflow.

Front panel of SYNAPT G2-Si showing the grounded union (A):

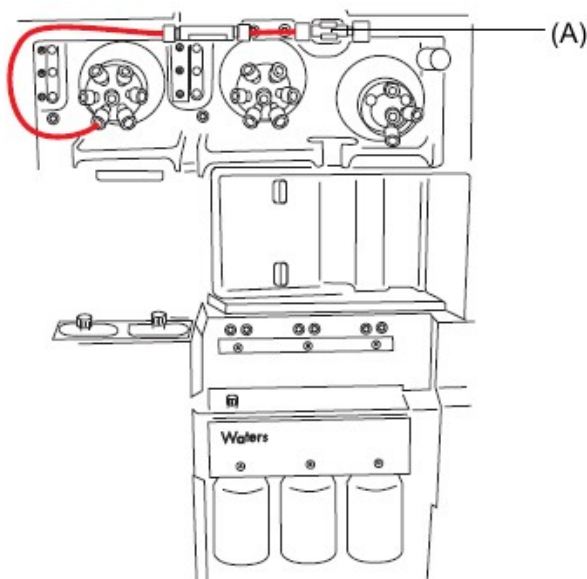


2. Install suitable sample, reference, and calibrant reservoir bottles in positions A, B, and C on the SYNAPT G2-Si. See also: *The Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide*, for installation instructions.

Recommendation: Use reservoir bottles in this order: A, for sample; B, for reference; C, for calibrant. Reference and calibrant solutions are described in their respective Help topics. Base sample solutions on your eventual application. Filtered or centrifuge them to avoid blocking the transfer line.

3. Open the Tune window, click File > Open, and open the default tune file.

Tip: At any time you can click File > Revert, to revert to previously saved parameters.

4. If the instrument is in source standby mode, click Operate , and confirm that the instrument's status indicator shows green.

5. From the Source menu, click NanoFlow.

Tip: The Source tab is labeled with the source type and instrument polarity. For example, NanoFlow+ means NanoLockSpray source and positive ion mode. The source tab controls source voltages, temperatures, and gas flows.

6. Click positive  or negative  ion mode in the toolbar.

7. Click MS mode .

Tip: MS/MS configuration uses the same source parameter settings as those used for the MS configuration. However, you must in addition provide collision energy and a set-mass.

8. On the source tab, set the parameters to these initial values, which are suitable for tuning by infusion:

Temperature values:

Parameter	Value
Source Temperature	80 °C

Gas flows:

Parameter	Value
-----------	-------

Parameter	Value
Cone Gas Flow	40 L/hr
NanoFlow Gas	0.3 Bar
Purge Gas Flow	100 L/hr

Refill:

Parameter	Value
Discharge current	0 micro amp

Note: The refill parameter is available only if you have chosen ETD as the acquisition mode.

- Consult the Instrument tab to ensure the values reflected there approximate the recommended instrument values:

Instrument tab values:

Parameter	Value
Trap Collision Energy	Off for MS mode
Transfer Collision Energy	Off for MS mode
Gas Controls	Automatic
Resolving Quadrupole	4.7 (MS/MS only)

- For MSMS mode, in the settings pane, click the drop-down and click Instrument Settings and set Collision Energy to 6 V.



Tip: For MS mode, set the collision energy to off. If, for MS/MS mode, you must fragment a sample, set the collision energy to on, and specify an appropriate voltage.

- On the Fluidics tab, in LockSpray fluidics control, select the sample reservoir bottle (Reservoir A).

Rationale: During this procedure, the reference fluidics line is plumbed to the sample sprayer and not the reference sprayer.

- If you did not change the reservoir bottle designation, click Purge .

Tip: If you changed the reservoir bottle designation in the previous step, the instrument performs an automatic purge.

- Enter a flow rate of 0.5 μ L/min.
- Enter a fill volume, for example 50 ml
- Click fill .
- Click  to start the flow.
- Set the appropriate values to observe the spectrum:

Typical acquisition values:

Parameter	Suggested Value
-----------	-----------------

Parameter	Suggested Value
Low Mass	100 Da
High Mass	1200 Da
Scan Time	1.0 sec
Data Type	Continuum or Centroid
MSMS Set Mass (only in MS/MS mode)	Enter a precursor ion mass in Da

18. In the plot control pane, click the drop-down, click Mode Settings, set the mode as TOF mode, and the scanning mode as MS or MS/MS.

Rationale: This mode is used in combination with the view settings to determine how the data is presented in the trace display. Quad mode allows displays of Intensity (chromatographic view), Spectrum view, Resolution view, Signal-to-Noise view. TOF mode allows a Spectrum view, Spectrum + Multi, Spectrum + Intensity, Spectrum + Resolution, and Spectrum + Mass Position.

19. Zoom the peak display to a suitable range around the mass of interest.

Tip: On the spectrum, click and drag across the desired range, and then release.

20. Use the x-, y-, z-axis position adjusters and video display to maximize the peak intensity.

21. Adjust the capillary voltage for a stable spray with maximum sensitivity.

Tip: Most compounds are optimized at between 1.0 and 3.0 kV.

22. Confirm that the shape of the displayed peaks is correct.

23. Adjust the purge gas to reduce or remove peaks caused by the ingress of external atmospheric contaminants.

Rationale: The purge gas provides positive pressure in the source to prevent the ingress of contaminants. Start with a value of 100 and tune according to your application. Setting the purge gas too high, can adversely affect spray stability.

24. Ensure that resolution, intensity, and stability are acceptable.

25. Click File > Save As and save the Tune file (.ipr) with a new name.

26. Reconnect the grounded union to the reference sprayer (see the *Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide*).

Starting with these settings, you can optimize the source configuration for your application, and store them as separate tune files (.ipr).

[Back to top](#)

Fine-tune the beam for your specific application:

1. Click File > Open, and open your default nanoflow tune file.
2. Start the external (LC) flow.
3. Use the x-, y-, z-vernier probe adjuster and video display to maximize the peak intensity.
4. Adjust the capillary voltage for a stable spray with maximum sensitivity.
5. Optimize the cone voltage for maximum intensity.

- Click File > Save As, and save the file with a new name.

Having tuned the instrument for your application, you can acquire nominal mass data. To obtain exact mass data, calibrate the instrument, and then configure and use the LockSpray.

Related Topics

- Prepare the instrument for use
- Automatic calibration
- Configure a lock-mass

Optimize Quadrupole Profile

- MassLynx > MS Console > Instruments: SYNAPT G2-Si > MS Tune > Setup > Quad Profile

The default quadrupole automatic profile is suitable for general mass spectroscopy. The ability to adjust the quadrupole radio frequency (RF) voltage configuration is available for experts, for fine tuning.

The transmission window of the quadrupole is a function of the quadrupole set-mass. A fixed quadrupole set-mass results in a limited mass range for MS acquisitions. To overcome this limitation, the quadrupole set-mass scans during acquisitions. Enabled by a default setting of Auto Profile, this scan function is automatic.


Where appropriate for your application, it is also possible to manually define a scanning profile or set a fixed quadrupole set-mass. The effect of changing the parameter settings can be observed in real time on the peak display.

Select from Fixed mass, Auto profile, or Manual profile as follows:

Fixed mass

The graph below shows a typical quadrupole transmission window for a fixed mass. For a mass-to-set-mass ratio of 0.8, transmission cuts off sharply and, at the high-mass end, as the mass approaches seven times the quadrupole set-mass, the cutoff is more gradual. Thus for a typical fixed set-mass of 500 Da, no ions are transmitted below 400 Da ($0.8 \times$ set-mass), the transmission remains high, at 2500 Da ($5 \times$ set-mass) and diminishes above 3500 Da ($7 \times$ set-mass).

Typical fixed mass transmission profile:

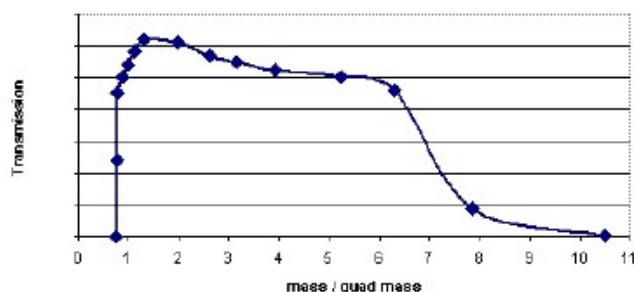


Auto profile

In Auto Profile mode the RF voltage on the quadrupole is scanned, giving an optimal transmission profile over the mass range of interest. The RF scanning range is a function of the acquisition start and end mass values.

Manual profile

When you select the Manual Profile mode, you can ramp the quadrupole set-mass, enabling the acquisition of mass ranges wider than the quadrupole transmission window. Specify appropriate values for mass, dwell time, and ramp time to suit the acquisition mass range.

Manual quad RF profile:

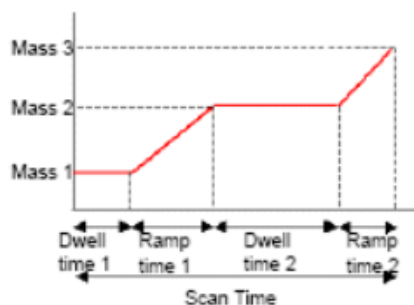

Typical manual quad RF profile parameter values:

Parameter	Value
Mass 1	100
Mass 2	300
Mass 3	500
Dwell Time (% Scan Time) for Mass 1	20
Dwell Time (% Scan Time) for Mass 2	20
Ramp Time (% Scan Time) for Mass 1	20
Ramp Time (% Scan Time) for Mass 2	40

Requirement: Dwell and Ramp times must add to 100%.

Status information from the front of the instrument:

Status LED strip	Instrument State
Off	Standby
Green	Operate
Amber	Source Standby
Flashing Green	Operate with pump override on
Flashing Amber	Standby with pump override on
Flashing Red	Not at vacuum
Red	RF Trip

Related Topics

- MS and MS/MS methods
- MSe method

Tune the mass spectrometer manually for ionKey/MS operation

- MassLynx > MS Console > Instruments: SYNAPT G2-Si > MS Tune



The ionKey/MS system integrates micro-scale UPLC separation into the source of the mass spectrometer by combining the ACQUITY UPLC M-Class and the SYNAPT G2-Si.

To enable an ionKey/MS system with an existing instrument, see the *SYNAPT G2-Si Operator's Overview and Maintenance Guide*

See also: The ACQUITY UPLC M-Class System Guide and the ionKey/MS System Guide.



Note: The Tune window is for advanced operation. Adjust its settings only if you are familiar with its features. The easiest and fastest way to tune the source is to allow IntelliStart to do it automatically.


To tune manually for ionKey/MS operation:

- In the Tune window, click Operate , and confirm that the adjacent instrument status indicator shows green.
- Click Setup, and then select Low Flow Mode.
- To select the source ionization mode, click Ion Mode, and select IonKey+ from the drop-down list.
- Click MS mode .
- On the ionKey+ Source tab, set the parameters to the recommended values, shown below.

Recommended ionKey+ Source tab parameter values:

Parameter	Description
Capillary	3.0 to 4.0 kV (typically lower for the Infusion Tile and Higher for the Analytical Tile).
Sampling Cone	40 V (tunable value, sample-dependent)
Source Offset	30 V (or as installed value for Std ESI)
Source Temperature	120 °C
Cone Gas	0 L/h (up to 150 L/h depending on background)
Nano Flow Gas	0.2 (tunable)
Trap Cooling Gas	0 (not configured for ionKey)

- Click API Gas  to start the nitrogen flow.
- Click Collision Gas  to start the trap gas flow.
- On the Vacuum tab, ensure the vacuum for the Trap gauge is approximately 8.0e-3 mbar.
- Insert the Ikey infusion tile and engage electronics/gas connections by rotating the handle 180 degrees in an anti-clockwise rotation on the right of the source enclosure.
- On the Fluidics tab, from the Reservoir list, select the bottle containing your tuning solution, and specify an infusion flow rate of 3 µL/min.

11. In the peak editor, on the right-hand side of the Tune window, select Function 1, and select a type of scan from the drop-down list.
12. Specify the mass for your compound, with an appropriate span and gain depending on the compound of interest (typically a span of 5, and a gain of 25).
13. Observe the peak in the Tune window.
14. To enable the camera, click 
Result: The Camera Control dialog box opens.
15. Using the camera controls, focus on the analyte sprayer.
16. Adjust the sprayer's X, Y, and Z controls to optimize the spray position for maximum signal intensity.
17. Adjust the capillary voltage for maximum signal intensity.
18. Adjust the cone gas flow to optimise the signal to background ratio.
Tip: Allow the pressure to stabilize after each adjustment.
19. Adjust the NanoFlow gas pressure, in increments of 0.1 bar, allowing the pressure to stabilize after each adjustment.
Tip: Making a small incremental change, such as from 0 bar to 0.1 bar will not be effective. It is better to make a large change, and then put it back to a 0.1 increment, and then allow the change to take effect.
20. Repeat steps 16 to 19 (inclusive), for maximum intensity and signal stability.

Related Topics

Related Topics

- [Prepare the instrument for use](#)
- [Automatic calibration](#)
- [Configure a lock-mass](#)

Tuning for Signal Intensity

- [MassLynx > MS Console > Instruments:SYNAPT G2-Si > MS Tune](#)

The SYNAPT G2-Si is optimally tuned by your installation engineer. In addition to this, you should optimize the source conditions for your application. For the appropriate instructions, click one of the links below.

Related Topics

- [Optimize lockspray for ESI](#)
- [Optimize for nanoflow ESI](#)

The Tune window

- [MassLynx > MS Tune](#)

The Tune window is the interface for instrument fine control. Use it for these purposes:

- Tuning the instrument
- Starting and stopping the flow of gasses
- Monitoring vacuum pressures
- Monitoring acquisitions

The Tune window includes parameter tabs, a status bar, and peak displays.

Peak display

The peak display shows the real-time signal on the detector, using the MS as the mass analyzer, without acquiring or making a data file. Resolution is shown for continuum data, and intensity of ions per push (IPP) is shown for centroid data.

Instrument control pages

You access the controls necessary to set up the system via the Tune window from the toolbar or instrument control tabs.

The source, Instrument and Fluidics tabs are always visible:

- Source tab – Controls source voltages, temperatures and gas flows.
- Instrument tab – Controls the quadrupole, collision energies (Trap and Transfer), and gas controls.
- Fluidics tab – Controls the fluidics system for sample, reference, and calibrant introduction.

Tip: The Source tab is labelled with the source type and instrument polarity, for example, ES+ means electrospray source and positive ion mode.

In addition, when you are using the ETD mode, the following tab is visible:

- TriWave tab - Controls

Tip: Click File > Revert, to revert to previously saved parameters.

Related Topics

- [Tuning for signal intensity](#)
- [Acquire data](#)

Acquiring data using the DESI source

- [MassLynx >File](#)

To acquire data

1. On the MassLynx home screen, select File - Import Worksheet.
2. Select the required experiment.
3. Save the file name in the file name tab.
4. In MS Tune window, select DESI source control.
5. Make sure that the DESI device is connected.
6. Specify the appropriate DESI configuration start point as defined in the topic, "Selecting the DESI source".
7. Click Start run.
8. Save the sample list.
9. Make sure that acquire sample and auto process samples are ticked.
10. Click OK

Result: The acquisition parameters are set up in the HDI imaging page.

Related Topics

- [Selecting the DESI source](#)
- [Defining a DESI imaging experiment](#)
- [Calibrate with DESI](#)
- [Calibrating the DESI zero position](#)

Calibrating with a DESI source

- [MassLynx >File](#)

Calibrating the instrument for DESI operation is similar to calibrating it in other, non-DESI modes. DESI calibrations must encompass the entire mass range you want to acquire. Acquiring data beyond the calibrated range diminishes the degree of mass accuracy returned by the instrument.

For mass data of consistent accuracy, calibrate the DESI SYNAPT G2-Si regularly against reference compounds appropriate to your project.

Calibrations created from MS acquisitions are valid for MS, MS/MS, and survey experiments.

Certain restrictions apply when calibrating a DESI source:

- Manual calibration is supported (you must identify peaks manually, and use previously acquired centroid data when you create or edit a calibration profile).
- Automatic/assisted calibration is supported and ensures correct attenuation of the beam before acquisition.
- The software does not support calibration checking.

Related Topics

- [Defining a DESI imaging experiment](#)
- [Selecting the DESI source](#)
- [Acquire with DESI](#)
- [Calibrating the DESI zero position](#)

Calibrate with a MALDI source

Calibrating for MALDI operation is similar to calibrating the instrument in non-MALDI mode. MALDI calibrations must encompass the entire mass range you want to acquire. Acquiring data beyond the calibrated range limits the mass accuracy returned by the instrument.

Tuning in Mobility TOF mode, to optimize ion mobility separation for ES sources, also applies to the MALDI source.

For mass data of consistent accuracy, calibrate the MALDI SYNAPT G2-Si regularly against reference compounds appropriate to your project. Calibrate daily, before long, automated data acquisitions or for each sample plate.

Calibrations created from MS acquisitions are valid for MS, MS/MS, and survey experiments.

Certain restrictions apply when calibrating a MALDI source:

- Only manual calibration is supported (you must identify peaks manually and use previously acquired centroid data when you create or edit a calibration profile).

- The software does not support calibration checking.

Related Topics

- Automatic calibration
- Operate the MALDI source


Use the MALDI pattern editor

- MassLynx> Instrument > MS Tune > Source >Maldi > Pattern Editor

A pattern file is a series of 2D-coordinates (nodes) between which the laser moves. You can also turn the laser on or off at each node. Using the pattern editor, you define a series of movements the laser makes within the area of the sample well, and save them to a file.



Use the MALDI pattern editor to create new non-imaging pattern files, or modify the example files provided with the MassLynx software. Non imaging pattern files are denoted by the ".ptn" extension.

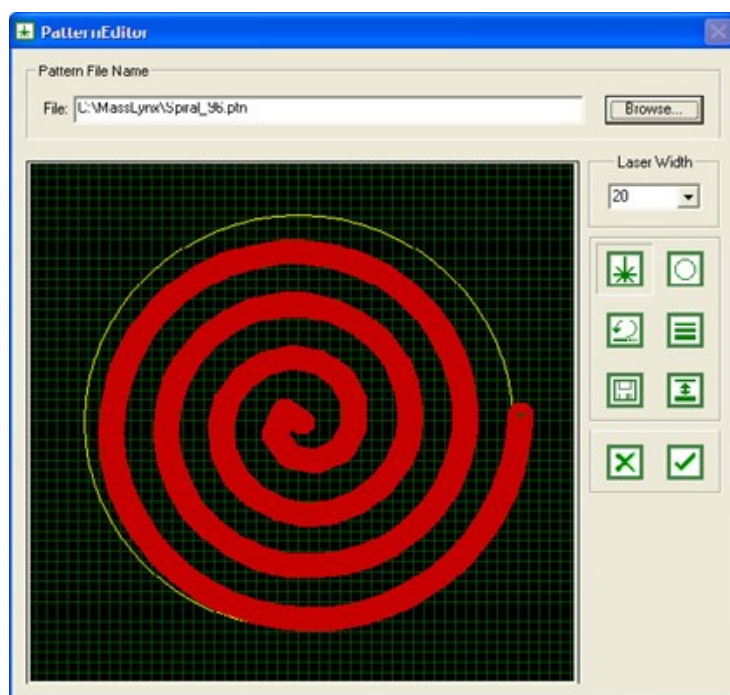The MALDI pattern editor:



To create a pattern:


In the MALDI Pattern Editor,



1. Select a laser width from the drop-down list.
2. Click Laser Width .
3. Click the grid (inside the yellow circle) to add the start node.
4. Click Fire Laser  to toggle the laser on.
5. Click the grid at the next desired node position.

6. Continue clicking the grid to add more nodes.

7. Click Save .

Result: A prompt appears asking you to enter a name for the pattern file.

You can also use Undo  to remove the last node added or New , to clear all pattern data from the grid.

Tip: The laser-width setting refers only to the width of the drawing lines, not the physical width of the laser beam. You can vary the laser width between nodes by selecting a new value from the drop down list and clicking Laser Width . Using a finer width , you can create more complex patterns within the sample well area. You can also apply the currently set laser width to all lines in the drawing by clicking Modify All Laser Widths .

Related Topics

- [Create methods for MALDI](#)
- [Operate the MALDI source](#)

Calibration using the REIMS source

- [MassLynx >File](#)

Calibrating with REIMS source

Calibrating the instrument for REIMS operation is similar to calibrating it in ESI mode. REIMS calibrations must encompass the entire mass range you want to acquire. Acquiring data beyond the calibrated range diminishes the degree of mass accuracy returned by the instrument.

For mass data of consistent accuracy, calibrate the REIMS SYNAPT G2-Si regularly against reference compounds appropriate to your project.

Calibrations created from MS acquisitions are valid for both MS and MS/MS.

Certain restrictions apply when calibrating the instrument using a REIMS source:

- Manual calibration is supported (you must identify peaks manually, and use previously acquired centroid data when you create or edit a calibration profile).
- Automatic/assisted calibration is supported and ensures correct attenuation of the beam before acquisition.
- The software does not support calibration checking.

Related Topics

- [Detector setup page and the REIMS source](#)
- [Selecting the REIMS source](#)

Detector setup page and the REIMS source

- [MassLynx >File](#)

The Detector Setup Page, when used with REIMS source, represents the selections in the Fluidics section as unavailable. The Current Tune Page Settings are therefore selected by default.

Related Topics

- Calibration with REIMS
- Selecting the REIMS source

Trouble shooting advice

Troubleshooting information provides advice for resolving issues that you can experience while operating the SYNAPT G2-Si. Waters recommends that you perform regular maintenance to prevent instrument problems. For maintenance procedures, see the *SYNAPT G2-Si Operator's Overview and Maintenance Guide*.

The instrument provides self-check facilities that identify some of these problems automatically. Others, however, such as loss of sensitivity, are caused by dirt or erosion from samples. You can also experience other instrument problems, like heater failure or vacuum leaks. For cleaning procedures, see the *Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide*.

Approach

For sample or lock-mass problems, determine source conditions via the software first. If all the settings appear correct, proceed to the physical instrument. Following are common physical instrument problems:

- The sample or reference compound concentration is too weak.
- The probe or capillary is leaking or blocked.
- Air is in the tubing.
- Gas pressure is too low.
- The probe is misaligned.

For further help and support, contact your local [Waters representative](#).

Related Topics

- Calibration, Lockspray or LC system failure
- Mass accuracy problems
- Troubleshooting poor-quality or absent beams
- Decreases in sensitivity
- Video problems

Calibration, Lockspray or LC system failure

This troubleshooting topic describes conditions that cause IntelliStart checks to return failure results.

Symptoms:

See "Calibration check returns a failing result:"

See "LockSpray check returns a failing result:"

See "LCMS system check failed:"

See "Detector check returns a fail:"

Calibration check returns a failing result:

Problem	Solution
Error message: Mean or RMS residual mass exceeds bounds	<p>Calibration acceptance criteria are configured in IntelliStart during calibration creation. When the acceptance criteria are exceeded, the check returns a failed result.</p> <p>Recalibrate the instrument (see "<i>Calibrate the SYNAPT G2-Si</i>").</p> <p>Review recent data for accuracy, using the Mass Difference Calculator, by clicking Windows Start > All Programs > MassLynx > Spectrum Tools > MassDiff.</p> <p>Tip: The Mass Difference Calculator compares peak information obtained experimentally with standard reference peaks. The calculator calculates the mass difference between each reference peak and the nearest entry in the peak list, reporting the difference as the absolute error (mDa) or fractional error (ppm).</p>

LockSpray check returns a failing result:

Problem	Solution
The intensity of the lock-mass beam is outside of the necessary bounds.	Run Intellistart lock-mass setup, see " <i>Configure a lock-mass</i> " in <i>Related Topics</i> .
The lock-mass reservoir bottle is empty	Refill the reservoir bottle

LCMS system check failed:

During a system check, MassLynx acquires data from an LC injection and compares the results with the value of acceptance parameters including retention time, peak area, peak height, peak width and the signal to noise ratio. Open the report to determine which criteria failed.

Note: The calculations do not include injections where no peak was found. For example, if only 8 of 10 injections detect peaks, the system check test results are based on 8 values rather than 10. Nevertheless, the system check process indicates a failure because of the missing peaks.

Problem	Solution
The report indicates that all checks returned failing results	<p>Where there is no ion beam, the wrong compound is in use, or the chromatographic peak falls outside the expected retention window, all criteria fail.</p> <ol style="list-style-type: none"> 1. Open the chromatogram via the MassLynx software, click Chromatogram > Realtime Update, to ensure that the chromatographic peak arrives at the expected time. (see Poor or absent beam in <i>Related Topics</i>) 2. Click Resolve to clear the system check error, and return the instrument to a state ready for data acquisition. 3. Before proceeding with other analyses, run the system check again.
Retention time incorrect	<p>Verify that the peak arrives inside the LCMS System Check retention time window.</p> <p>Adjust the chromatographic conditions, or adjust the LCMS System Check parameters.</p>
Chromatographic peak found but mass peaks are too weak or exhibit poor mass measurement	Determine the MS source conditions, and ensure that the calibration is valid.

Problem	Solution
Reproducibility - peaks are found consistently but their intensity varies	Determine the MS source conditions. Determine the chromatographic method, particularly the injector setup, see <i>ACQUITY online Help</i> .

Detector check returns a fail:

Problem	Solution
Detector check returned a failing result	Run the Detector Setup function.

Related Topics

- Calibrate the SYNAPT G2-Si
- Configure a lock-mass
- Check and set-up the detector
- Troubleshooting poor-quality or absent beams

Troubleshooting poor-quality or absent beams

A variable beam or no beam is generally caused by either a lack of sample entering the probe or ionization problems.

Symptoms:

See "Sample flow inconsistent:"

See "No ESI ionization:"

See "Blocked sample cone:"

Sample flow inconsistent:

Problem	Solution
Air in tubing	Disconnect the PEEK tubing from the probe, and observe liquid flow.
Solvent delivery is blocked	Inspect the probe to determine whether it is blocked. If it is blocked, replace the probe capillary. See <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .

No ESI ionization:

Problem	Solution
No capillary voltage	Contact your Waters service engineer.

Blocked sample cone:

Problem	Solution
The sample cone is blocked	Determine the backing pressure in the Tune window. Normal backing pressure is approximately 2 mbar. As the cone becomes blocked, the backing pressure decreases. If the pressure falls below 1.5 mbar, the cone is substantially blocked. Note: A high backing pressure indicates a leak or a pump problem. Clean the cone, see <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .

Problem	Solution
---------	----------

Decreases in sensitivity

A reduction in the ion count below the expected count suggests reduced sensitivity, which can mean that expected peaks go undetected or compounds go unidentified.

Symptom: Reduced sensitivity:

Problem	Solution
Source components became contaminated by sample, over time	Clean the source components. See the <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .
Probe tip became contaminated by sample, over time	Clean or replace the ESI probe tip. See the <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .
Probe tip is eroded by samples	Research solvent compatibility and replace the probe tip. See the <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .
The heater fails to heat the ion block	Replace the ion block heater cartridge. See the <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .
Solvent delivery is blocked	Replace the ESI probe capillary. See the <i>Waters Operator's Overview and Maintenance Guide</i> . SYNAPT G2-Si

Inaccurate peaks

Problems

TOF instruments rely on lock-mass correction to provide accurate mass data.

Peaks do not appear at their accurate mass

Problem	Solution
Without lock mass correction, peaks can fail to appear at their accurate masses.	If lock mass correction is not available, for instance on the APPI/APCI source, monitor the instrument's calibration regularly, and recalibrate when necessary.
Laboratory temperature varies widely through the sample run.	Use lock mass, and maintain the machine's physical environment as close to a set temperature as possible.

Mass accuracy problems

When mass peaks fail to appear at their accurate mass, several tests can help resolve the issue.

Symptoms:

See "Mass scale exhibits a permanent offset:"

See "Erratic mass readings across spectra:"

See "Poor mass accuracy:"

See "Mass values not indicated:"

Mass scale exhibits a permanent offset:

Problem	Solution
Lock-mass is not functioning	<p>If lock-mass data are not being acquired, then mass peaks are not corrected to their accurate mass over time.</p> <ol style="list-style-type: none"> 1. Confirm that the lock-mass peaks lie within the selected calibration mass-range. <p>Rationale: Where the calibration mass-range is not suitable for the lock-masses, the lock-mass peaks are not recognized.</p> <ol style="list-style-type: none"> 2. Confirm that the lock-mass reservoir bottle contains sufficient solution. 3. Run a lock-mass setup check.
Calibration is not valid for this acquisition	<p>Changes to instrument and source configuration can affect the calibration.</p> <p>Run a calibration check, to ensure the calibration is valid.</p>

Erratic mass readings across spectra:

Problem	Solution
Lock-mass intensity unstable or too low	<p>As the intensity of the lock-mass peaks pass below the usable threshold, the lock-mass correction switches off. When the intensity is useable again, the correction is applied. Where the lock-mass intensity is close to the useable boundary, the correction can switch on and off for subsequent spectra, giving erratic mass readings.</p> <p>From the MassLynx main window, click Chromatogram > Display > TIC, and select the LockSpray function from the available channels. Click the chromatogram, and tab through the spectra using the arrow button. Where the lock-masses were used, the peaks are colored blue. If the lock-masses were not used on some spectra, they are missing or colored red.</p> <p>Increase the concentration or flow rate of the reference compound.</p>
Mass value flips between two similar values	<p>Where interference from nearby peaks exists, the calculation of center-mass for the Gaussian peak can change which mass is identified between spectra.</p> <p>Where interference is due to lock-mass peaks, try a more suitable lock-mass compound.</p>

Poor mass accuracy:

Problem	Solution
Poor mass accuracy at high or low mass ends of spectra	<p>The acquisition can be outside the calibration mass range.</p> <p>For manual calibrations, the mass range is not stored, so it is possible to acquire data that are outside the calibration range.</p> <p>Create a calibration for the required mass range.</p>

Problem	Solution
<p>The acquisition mass range implies a different pusher frequency from the calibration mass range</p>	<p>Within the TOF region of the instrument, the pusher frequency decreases, in steps, to compensate for the longer flight times of larger molecules. The frequency decreases for the whole acquisition as the mass range crosses each of the following higher mass boundaries:</p> <ul style="list-style-type: none"> • 1200 Da • 2000 Da • 5000 Da • 8000 Da • 14000 Da • 32000 to 100000 Da <p>This pusher behavior also applies to calibrations. Calibrations are valid only for a specific pusher frequency, as decided by the highest mass boundary included in their mass range.</p> <p>For example, an acquisition from 50 to 1000 Da uses the highest pusher frequency. If the active calibration has a range of 50 to 1300 Da, it is not valid for the acquisition because this calibration crosses the first mass boundary.</p>

Mass values not indicated:

Problem	Solution
<p>Peaks of interest are too small or appear split</p>	<p>Although the spectrum can show the mass peaks of interest, sufficient data are necessary to calculate a statistically significant Gaussian curve. The mass value provided is then the center-mass of the calculated Gaussian distribution.</p> <p>Increase the concentration or flow rate of the analyte compound.</p>


Related Topics

- [Prepare the instrument](#)
- [Automatic calibration](#)
- [Configure a lock-mass](#)

Video problems

The video camera on the nanoLockSpray source is connected to the workstation PC through the mass spectrometer. A BNC connection on the instrument's back panel connects to a USB converter box. The USB cable then connects the converter box to the workstation PC.

Camera viewer

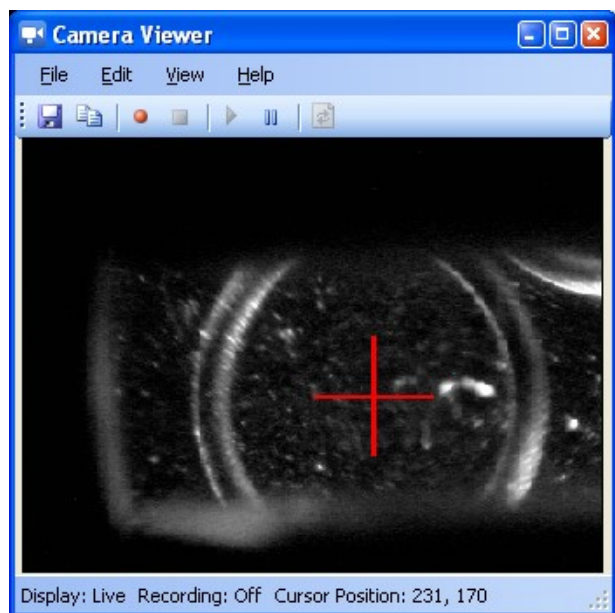

Symptoms:

- No video (nanoLockSpray source)
- Focus and position problems

No video:

Problem	Solution
Video not connected	Confirm the video connections are correct and secure. See the <i>Waters SYNAPT G2-Si Operator's Overview and Maintenance Guide</i> .
Windows video driver not functioning	Try the following measures, in turn: Close and reopen the video window. Move the USB lead to a different USB port on the workstation PC. Close the MassLynx software and restart the workstation PC.

Focus and position problems:

Problem	Solution
The camera is pointing away from the sprayer tip	Manipulate the camera mount to view the sample sprayer. Tip: To aim, move the camera mount. To focus, rotate the mount.
The image is out of focus	Turn the camera mount to adjust the focus.

Error messages

The instrument control computers contain many error messages that may display from time to time. In general, the message are self explanatory and provide advice on what action to take or reason for failure (usually incorrect settings or missing information).

In many cases problems are caused by communication breakdown between the instrument and control software running on the embedded PC. Rebooting the PC could provide a solution. See the SYNAPT G2-Si *Operator's Overview and Maintenance Guide* for reboot instructions.

This table contains a list of errors that need supplementary advice.

Message	Action or description
Source pressure monitoring is indicating a possible exhaust problem.	See the SYNAPT G2-Si <i>Operator's Overview and Maintenance Guide</i> Appendix B - Connecting the Nitrogen exhaust line
The source pressure test failed. Possible source leak.	See the SYNAPT G2-Si <i>Operator's Overview and Maintenance Guide</i> Chapter 5 - Cleaning the source components
The source pressure test failed to complete due to high exhaust pressure.	See the SYNAPT G2-Si <i>Operator's Overview and Maintenance Guide</i> Appendix B - Check rotary pump
The source pressure test failed to complete due to low API gas pressure.	See the SYNAPT G2-Si <i>Operator's Overview and Maintenance Guide</i> Appendix B
A leak has been detected in the MS fluidics system.	See the SYNAPT G2-Si <i>Operator's Overview and Maintenance Guide</i> Chapter 5

Prepare to acquire data

- [MassLynx > MS Console > Instruments: SYNAPT G2-Si > Intellistart](#)

Once you have set up the instrument, you then need to prepare for data acquisition. Before you perform the preparation tasks detailed below, ensure you set up the calibration and set up the lockspray (see [Related Topics](#))

Perform the following tasks to prepare the SYNAPT G2-Si for data acquisition.

Requirements: The instrument is configured for your application. The Waters engineer who installed your equipment will have configured and calibrated it for your applications but if you are changing the type of use, you may need to get these settings altered. If this is the case, contact Waters.

To prepare for data acquisition:

Before starting the SYNAPT G2-Si, you must complete the following tasks:

- Ensure that the instrument is pumped down (evacuated) and the ionization source for your application fitted. See the *Waters SYNAPT G2-Si Mass Spectrometry System Operator's Overview and Maintenance Guide* on the SYNAPT G2-Si Documentation CD.
- Visually inspect the ion probe position.
- Ensure that the wash bottle contains a sufficient volume of a suitable wash solution.

Tip: The composition of the wash solution is sample-dependent. However, for most applications, Waters recommends an 80:20 methanol/water mixture.

- Ensure that your sample, lock-mass reference, and calibrant solutions are ready and mounted on the instrument.

Recommendation: Use reservoir bottles in this order: A for sample, B for reference, C for calibrant. Compose the sample solution according to the application, filtering it or centrifuging to avoid blocking the transfer line.

- Ensure the waste line is correctly plumbed.
- Ensure the LC system is configured and ready to use.

Related Topics

- Start and stop the instrument
- Vent and shut down the instrument
- Prepare the instrument for use
- Automatic calibration
- Configure a lockmass


Acquire data

- [MassLynx > MS Tune > Tune window](#)

You normally acquire data using methods selected in the MassLynx sample list. However, to manually start acquisition of data for a single sample, use the Tune window.

To acquire data from the Tune window:



1. Click .
2. Click Acquire > Start.
3. In the Start Acquisition dialog box, choose the function type, lock-mass, data format, and dynamic range.
4. Enter a run time and mass range.
5. Click Start.

Result: The mass spectrometer begins to acquire data and save them to disk.

Related Topics

- [Acquisition Methods](#)

Acquisition methods

- [MassLynx > MS Method](#)

Your SYNAPT G2-Si makes a variety of functions and experiments available for you to customize.

Related Topics

- [Fast-DDA method](#)
- [Product Ion Discovery \(PID\) method](#)
- [MS^E method](#)
- [Control fluidics and external devices](#)



MS and MS/MS methods

- [MassLynx > MS Method > MS|MSMS](#)

The MS and MS/MS functions enable you to create a custom experimental method. You define the duration of MS and MS/MS functions during set up. The MS/MS acquisition is targeted at known or expected molecular species of interest by specification of m/z and retention time values.

The manual MS/MS method set up does not provide automated operations such as the MS to MS/MS switching provided by Survey, Fast DDA and Mobility DDA methods.

To create methods with MS and MS/MS functions:

1. From the main MassLynx window, click MS Method.
2. Combine multiple MS and MS/MS functions in a single method.
 - For MS functions – Click MS  MS to acquire spectra from precursor ions.
 - For MS/MS functions – Click MS/MS  MSMS to acquire spectra of product ions formed from selected precursors in the collision cell.
3. Specify your MS or MS/MS function on the following tabs:
 - **Acquisition** – Configure the initial acquisition parameters.
 - **TOF MS** – Configure the MS parameters (MS tab).
 - **TOF MS/MS** – Configure the MS/MS parameters (MS/MS tab).
 - **Trap CE Control** – Configure the trap collision energy (CE) parameters.
 - **Transfer CE Control** – Configure the transfer collision energy (CE) parameters.

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Fragmentation: Select the fragmentation mode. Polarity: Select the polarity. Analyser Mode: Select the type of optics required. Dynamic range: Select either Normal or Extended. Extended is not available for ETD fragmentation. Sensitivity: Normal or Enhanced Duty Cycle (EDC) for a selected mass.

[Back to top](#)

TOF MS tab

On the TOF MS tab, you configure the following parameters:

- m/z range over which to acquire data
- Conditions at which the instrument scans
- Instrument conditions

Acquiring TOF MS over the m/z range

Parameter	Description
Low Mass	The start m/z value
High Mass	The end m/z value

Scanning Conditions

Parameter	Description
Scan Time	Specifies the duration over which TOF ion-detections accumulate in each individual spectrum presented to MassLynx.
Data Format	Specifies the type of data collected (centroid or continuum).

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in the tune file	Select this override and enter the preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear ramp between these two values.

[Back to top](#)

TOF MS/MS tab

Using the TOF MS/MS tab, you can configure the following parameters:

- m/z range over which to acquire data (described below)
- Conditions at which the instrument scans
- Instrument conditions

Acquiring MS/MS over the m/z range

Tip: When analyzing ions with multiple charges, ensure that the mass range is sufficient to include product ions with a lower charge but higher m/z value than the precursor.

Parameter	Description
Low Mass	The start m/z value over which data are acquired.
High Mass	The end m/z value at which data acquisition stops.
Set Mass	Specifying the precursor mass selected by the quadrupole for fragmentation.

Scanning Conditions

Parameter	Description
Scan Time	Specifies the duration of each scan in seconds.
Data Format	Specifies the type of data collected (centroid or continuum). See "Creating methods with MS method editor".

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in tune file	Select this override, and enter the preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear voltage ramp between these two values.

[Back to top](#)

Trap CE and Transfer CE Control tabs

The default collision energies are set in the Tune window.

Note: The Trap CE Control settings are not available when you choose the ETD fragmentation mode.

Alternatively, you can control the trap and transfer collision energies for this method using the following options:

Fixed Collision Energy Value

Using this option you can specify a single value that is used throughout the function.

Collision Energy Profile

Using this option you can select up to five collision energies. The collision energy cycles around the selected values, using one collision energy per accumulated spectrum, to produce optimal fragmentation. The following energy values can be specified:

- Collision Energy 1 - For survey scans
- Collision Energy 2 - Low mass for MS/MS scans
- Collision Energy 3 - High mass for MS/MS scans
- Collision Energy 4 - Start for low mass ramp for MS/MS scans
- Collision Energy 5 - End for low mass ramp for MS/MS scans

Collision Energy Ramp

For the duration of the scan time, while ion detections are accumulating into each individual spectrum, the collision energy follows a linear voltage ramp between these two values.

[Back to top](#)

Related Topics

- [Adding lock-mass correction methods](#)
- [Creating Mobility MS and Mobility MS/MS methods](#)

Mobility MS and MS-MS methods



- [MassLynx > MS Method > Functions > Mobility MS | Mobility MS/MS](#)

The mobility MS and MS/MS functions enable you to employ ion mobility separations to increase the selectivity of the experiment for compounds of interest by separating them on the basis of their mobility as well as providing a route to cleaner fragment ion spectra that are easier to interpret (supported in the MSE viewer tool of MassLynx). The Mobility MS and Mobility MS/MS functions enable creation of a custom experimental method where MS/MS

acquisitions are targeted using specification of m/z and retention time values of known or expected molecular species of interest.

The manual MS/MS method set up does not provide automated operations such as the MS to MS/MS switching provided by Survey, Fast DDA and Mobility DDA methods.

To create methods with Mobility MS and Mobility MS/MS functions:

- From the main MassLynx window, select MS Method.
- Combine multiple Mobility MS and Mobility MS/MS functions in a single method.
 - For HDMS MS functions – Click MS  MS to acquire spectra from precursor ions.
 - For HDMS MS/MS functions – Click MS/MS  MSMS to acquire spectra of product ions formed from selected precursors in the collision cell.
- Specify your MS or MS/MS function on the following tabs:
 - Acquisition** – Configure the initial acquisition parameters.
 - TOF MS** – Configure the MS parameters (MS tab).
 - TOF MS/MS** – Configure the MS/MS parameters (MS/MS tab).
 - Mobility** - Configure data collection data acquisition parameters.
 - Trap CE Control** – Configure the trap collision energy (CE) parameters.
 - Transfer CE Control** – Configure the transfer collision energy (CE) parameters.

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Fragmentation: Select the fragmentation mode. Polarity: Select the polarity. Analyser Mode: Select the type of optics required. Sensitivity: Normal, Enhanced Duty Cycle (EDC) or High Duty Cycle (HDC)

[Back to top](#)

TOF MS tab

On the TOF MS tab, you configure these parameters:

- m/z range over which to acquire data
- Conditions at which the instrument scans
- Instrument conditions

Acquiring TOF MS over the m/z range

Parameter	Description
Low Mass	The start m/z value
High Mass	The end m/z value

Scanning Conditions

Parameter	Description
Scan Time	Specifies the duration over which TOF ion-detections accumulate in each individual spectrum presented to MassLynx.
Data Format	Specifies the type of data collected (centroid or continuum).

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in tune the file	Select this override and enter the preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear voltage ramp between these two values.

[Back to top](#)

TOF MS/MS tab

Using the TOF MS/MS tab, you can configure these parameters:

- m/z range over which to acquire data (described below)
- Conditions at which the instrument scans
- Instrument conditions

Acquiring MS/MS over the m/z range

Parameter	Description
Low Mass	The start m/z value over which data are acquired.
High Mass	The end m/z value at which data acquisition stops.
Set Mass	Specifying the precursor mass selected by the quadrupole for fragmentation. <div style="display: flex; align-items: center;"> <div style="font-size: 2em; margin-right: 10px;">!</div> <div> <p>Notice : When analyzing ions with multiple charges, ensure that the mass range is sufficient to include product ions with a lower charge but higher m/z value than the precursor.</p> </div> </div>

Scanning Conditions

Parameter	Description
Scan Time	Specifies the duration of each scan in seconds.
Data Format	Specifies the type of data collected (centroid or continuum). See "Creating methods with MS method editor".

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in tune file	Select this override, and enter the preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear voltage ramp between these two values.

[Back to top](#)

Mobility tab

Using the Mobility tab, you can configure data collection data acquisition parameters.

Data Collection

Select the Add Drift Time Function to save collapsed retention time data to a separate function in the Raw file.

Data Acquisition

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave settings will be used.
Apply Rule File for Charge State / Drift Time Stripping	Select this option to filter the acquired data using a rule file (.rul).

[Back to top](#)

Trap CE and Transfer CE Control tabs

Note: The Trap CE Control settings are not available when you choose the ETD fragmentation mode.

The collision energy is set on the Tune window. However, you can choose to control collision energy directly for this function using the following options:

Fixed Collision Energy Value

Using this option you can specify a single value that is used throughout that function. This value will override the value specified in your instrument parameter file.

Collision Energy Profile

Using this option you can select up to five collision energies. The collision energy cycles around the selected values, using one collision energy per accumulated spectrum, to produce optimal fragmentation. The following energy values can be specified:

- Collision Energy 1 - For survey scans
- Collision Energy 2 - Low mass for MS/MS scans
- Collision Energy 3 - High mass for MS/MS scans
- Collision Energy 4 - Start for low mass ramp for MS/MS scans
- Collision Energy 5 - End for low mass ramp for MS/MS scans

Collision Energy Ramp

For the duration of the 'scan time' while ion detections are accumulating into each individual spectrum, the collision energy follows a linear voltage ramp between these two values.

Collision Energy From Look Up Tables

Use a file of values to specify the collision energy during the scan.

[Back to top](#)

Related Topics

- [Adding lock mass correction to methods](#)

MS^E_method

- [MassLynx > MS Method > MS^E Continuum | MS^E Centroid](#)



The MS^E acquisition mode is an extremely useful generic, untargeted data acquisition method where high quality quantitative and qualitative information is obtained from every detectable component in a chromatographic run. In other words, the MS^E method enables comprehensive exact-mass precursor and product ion information to be acquired from all ion species with a single injection in a data-independent manner. As such, utilization of this method requires no prior "knowledge" of the sample to drive the fragment ion analysis (fragmentation is conducted in the StepWave or TRANSFER T-wave regions) . It provides simultaneous precursor, product ion, and neutral-loss detection, and when combined with the relevant MassLynx application manager, quick and easy interpretation and reporting of results can be obtained.

Note: This experiment type is conducted in continuum or centroid mode and supports only CID fragmentation (not ETD).

MS^E data can be acquired in two formats, continuum and centroid:

- The continuum data format preserves peak shape and is used in protein identification methods. For continuum data, the detector acquires several thousand complete TOF spectra per second (depending on the mass range). These are mass-corrected and summed in the embedded PC for the defined "Scan Time". One complete spectrum is saved to disk per scan.
- The centroid data format is used for small-molecule identification. The raw spectra from the detector are summed for the duration of the scan time, and then the data are converted to "stick" form, preserving mass and intensity information, but losing peak shape. The file size is smaller than for continuum data.

To create an MS^E method:

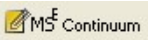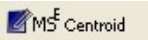1. Click  , or  .
2. Select the source mode, polarity, and specify the acquisition times.

Tip: Acquisition times are measured, in minutes, from the start of the method.

3. Select an analyser mode of either Resolution or Sensitivity.
4. Specify your scan parameters on the TOF MS tab according to the descriptions in the following tables.

TOF MS parameters:

Parameter	Description
Mass Range	Enter the start and end values, in Da.

Parameter	Description
Scan Time	Specify the duration over which ion-detections accumulate in each individual spectrum recorded.
Data Format	Continuum or Centroid as chosen by the experiment type.

- Enter the collision energies.

Tip: An MS^e acquisition involves recording two MS functions as the instrument cycles between the low and high collision energies.

Low and high collision energy parameters:

Parameter	Description
Low Energy	Set a low energy (typically 4 V), during which no fragmentation occurs to the precursor ions.
High Energy Ramp	Set the elevated energy ramp, in MS mode, to linearly ramp the energy (typically between 15 V and 40 V). Linear ramping induces fragmentation of any species present in the gas cell at that time.

- If required specify an alternate cone voltage.
- Click OK.
- In the Experiment Setup window, click LockSpray.

Requirement: MS^E acquisitions require exact mass, so you must enable LockSpray ionization.

Tips:

- Where your source enclosure supports it, you can use a lock-mass. The LockSpray sources support lock-mass, through use of an electrospray probe and baffle. The settings presented in the tables below are necessary for acquiring lock-mass data from any source that supports lock-mass acquisition.
- The LockSpray source runs at a flow rate that gives a stable spray, typically in the range 5 to 20 $\mu\text{L}/\text{min}$. If you must run sample lists over a weekend, use an external reservoir. These are the calculated times for a 30-mL reservoir at various lock-mass flow rates:

Lock-mass flow rates and running times (30 mL reservoir):

Flow rate ($\mu\text{L}/\text{min}$)	Running time (hours)
5	100
10	50
20	25

- Choose to acquire LockSpray with or without correction applied.

Tip: Some data applications need correction data to be separate.

- Select a reference compound.
- Specify the lock-mass acquisition settings.

Requirement: Also complete these values for APGC sources.

LockSpray Settings:

Parameter	Description	Typical short-run	Typical long-run
Scan Time	The dwell time for the lock-mass acquisition	0.3 s	1.0 s
Interval	The time between lock-mass acquisitions	20 s	60 s
Mass Window +/-	Adjusts the window size to locate the correct lock-mass peak in the LockSpray spectrum	0.3 Da	0.3 Da

12. Click OK.

Related Topics

- Prepare the instrument for use

Mobility MS^E method

- [→ MassLynx > MS Method > Mobility MS^E > Continuum](#)

The High Definition MS^E, or HDMS^E, method provides all of the advantages of the generic MS^E methodology described in the MS^E method section above however it is particularly useful for use when the selectivity of the analysis needs to be maximized for e.g with the use of very complex samples where many molecular species can co-exist at the same point on the m/z scale. Ion mobility is used to separate ionic species prior to MS or MS^E analysis (where fragmentation is conducted in the TRANSFER T-wave only).

Note: This experiment type is not supported for ETD fragmentation mode.

Creating a Mobility MS^E method:

1. Click Mobility MS^E Continuum .
2. Specify your function on the following tabs:
 - Acquisition
 - ToF MS
 - Cone Voltage
 - Mobility

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter settings for the relevant function.

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Polarity: Select the polarity. Analyser Mode: Select the type of optics required.

Tof MS tab

On the TOF MS tab, configure the following parameters:

- Da range
- Scanning conditions
- Da range
- Acquire Mobility MS^E over the range specified by the start and end values.
- Scanning conditions

Parameter	Description
Scan Time	Specifies the duration over which TOF ion detections accumulate into each individual spectra and are presented to MassLynx software.
Data Format	Specifies the type of data collected (centroid or continuum).
Collision Energy tab	When acquiring Mobility MS ^E data, two MS functions are used in an alternating fashion. During the course of an experiment, the instrument cycles between the low and elevated collision energies.

Parameter	Description
Function 1 -	Low Energy Trap Collision Energy: Can switch on or off, and specify voltage. Transfer Collision Energy: Can switch on or off and specify voltage.
Function 2 -	High Energy Ramp Trap Collision Energy: Can switch on or off and specify a voltage range. Ramp Transfer Collision Energy: Can switch on or off and specify a voltage range. Use collision energy from lookup tables: Select files created in the lookup table editor (see the link in Related Topics).

Cone Voltage tab

If required, override the cone voltage value specified in the tune file with an alternative voltage.

Mobility tab

Use the mobility tab to specify rule files for charge state and drift time stripping for both low and high collision energy functions. The stripping lets you filter the data acquired to the RAW data file for charge state and drift time.

Notes:

The stripping is not applied to data acquired to the Mobility data file. The full data are always acquired to the Mobility data file so that filter operations can be performed using the DriftScope application.

Filtering is not possible as a post process step for the RAW data file because the mobility information is not stored.

Parameter	Description
-----------	-------------

Parameter	Description
Function 1 -	Low Energy Add a rule file for charge state / drift time stripping at low energy. Tip: Rules are created in DriftScope and saved in files with a .rul extension.
Function 2 -	High Energy Add a rule file for charge state / drift time stripping at low energy.

Creating a LookUp table

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Fragmentation: Select the fragmentation mode. Polarity: Select the polarity. Analyser Mode: Select the type of optics required. Sensitivity: Normal, Enhanced Duty Cycle (EDC) or High Duty Cycle (HDC)

[Back to top](#)

Tof MS tab

On the TOF MS tab, you configure these parameters:

- m/z range over which to acquire data
- Conditions at which the instrument scans
- Instrument conditions

Acquiring TOF MS over the m/z range

Parameter	Description
Low Mass	The start m/z value
High Mass	The end m/z value

Scanning Conditions

Parameter	Description
Scan Time	Specify the duration over which TOF ion-detections accumulate in each individual spectrum presented to MassLynx.
Data Format	Specify the type of data collected (centroid or continuum).

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in tune the file	Select this override and enter the preferred cone voltage.

Ramp the Cone Voltage during the scan	Select this option and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear voltage ramp between these two values.
---------------------------------------	--

[Back to top](#)

TOF MS/MS tab

Using the TOF MS/MS tab, you can configure these parameters:

- m/z range over which to acquire data (described below)
- Conditions at which the instrument scans
- Instrument conditions

Acquiring MS/MS over the m/z range

Parameter	Description
Low Mass	The start m/z value over which data are acquired.
High Mass	The end m/z value at which data acquisition stops.
Set Mass	Specifies the precursor mass selected by the quadrupole for fragmentation. <div style="display: flex; align-items: center;"> <div style="font-size: 2em; margin-right: 10px;">!</div> <div> <p>Notice: When analyzing ions with multiple charges, ensure that the mass range is sufficient to include product ions with a lower charge but higher m/z value than the precursor.</p> </div> </div>

Scanning Conditions

Parameter	Description
Scan Time	Optionally specify the duration of each scan in seconds.
Data Format	Specify the type of data collected (centroid or continuum). See "Creating methods with MS method editor".

Instrument Conditions

Parameter	Description
Override Cone Voltage Specified in tune file	Select this override, and enter the preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. While ion-counts are accumulating into each individual spectrum acquisition, the cone voltage follows a linear voltage ramp between these two values.

[Back to top](#)

Mobility

Using the Mobility tab, you can configure data collection data acquisition parameters.

Data Collection

Select the Add Drift Time Function to save collapsed retention time data to a separate function in the Raw file.

Data Acquisition

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave settings will be used.
Apply Rule File for Charge State / Drift Time Stripping	Select this option to filter the acquired data using a rule file (.rul)

[Back to top](#)

Trap CE and Transfer CE Control tabs

Rule: The Trap CE Control settings are not available when you choose the ETD fragmentation mode.

The collision energy is set on the Tune window. However, you can choose to control collision energy directly for this function using the following options:

Fixed Collision Energy Value

Using this option you can specify a single value that is used throughout that function. This value will override the value specified in your instrument parameter file.

Collision Energy Profile

Using this option you can select up to five collision energies. The collision energy cycles around the selected values, using one collision energy per accumulated spectrum, to produce optimal fragmentation. The following energy values can be specified:

- Collision Energy 1 - For survey scans
- Collision Energy 2 - Low mass for MS/MS scans
- Collision Energy 3 - High mass for MS/MS scans
- Collision Energy 4 - Start for low mass ramp for MS/MS scans
- Collision Energy 5 - End for low mass ramp for MS/MS scans

Collision Energy Ramp

For the duration of the 'scan time' while ion detections are accumulating into each individual spectrum, the collision energy follows a linear voltage ramp between these two values.

Collision Energy From Look Up Tables

Use a file of values to specify the collision energy during the scan.

[Back to top](#)

HD_DDA_method

▪ [MassLynx](#) > [MS Method](#) > [HD-DDA](#)

SYNAPT G2-Si also features a new HD-DDA method which supersedes the previous Mobility Survey method.

HD-DDA features a smarter algorithm (like FastDDA) that provides more efficient timing of MS/MS acquisition, and more sensitive ion detection through the integration of the (ion-mobility enabled) Wideband Enhancement mode of operation. The combination of these leads to higher quality data from low abundance species and the greatest number switches possible per unit time.

Alternatively, ion mobility can be used to remove interferences by precursor drift time matching.

It can be used to perform one of the following tasks:


- Enhance quality of MS/MS information from weak samples or
- To improve identification of compounds from mixtures of simple/medium complexity.

You assign the criteria for selecting ions for an MS/MS acquisition in this order:

1. Threshold - MS Survey tab
2. Peak Selection - Peak Detection tab
3. Exclude list - Exclude tab
4. Include list - Include tab

During the survey, once the peaks are selected, the MS/MS cycles through them until the Stop MS/MS criteria are met (MS/MS tab).

To create a survey experiment:

1. Select  HD-DDA.
2. Specify parameter values on these tabs:
 - Acquisition
 - MS Survey
 - MS/MS
 - Peak Detection
 - Exclude
 - Include
 - Collision Energy
 - Mobility

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Acquisition settings:

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Polarity: Select the polarity. Analyser Mode: Select the type of optics required. Improved Sensitivity: When Improved Sensitivity is selected, MSMS data is acquired in Sensitivity mode. This only applies in Resolution mode.

[Back to top](#)

MS Survey tab

Using the MS Survey tab, you can configure the following parameters:

- MS survey Da range – The mass range (m/z) over which to acquire data.
- Stop the MS survey – The conditions at which the survey scan stops and the acquisition switches to MS/MS mode.
- MS survey scanning conditions – The conditions for which to acquire data.

MS survey Da range:

Parameter	Description
Start	The lower m/z limit
End	The upper m/z limit
Max No. of Ions	Select 1 - 8.

MS Survey Switching Threshold:

Allows the threshold of intensity to be set for switching to MS/MS acquisition.

MS survey scanning conditions:

Parameter	Description
Scan Time	Specifies the duration over which TOF ion-detections accumulate in each spectrum presented to MassLynx.
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

MS/MS tab

In the MS/MS tab, you can configure these parameters:

Parameter	Description
Stopping MS/MS	The conditions at which an MS/MS acquisition stops and switches back to an MS survey acquisition
Scan rates and instrument conditions for MS/MS	The conditions under which MassLynx acquires data

Tip: The parameter values you specify to stop an acquisition depend upon the type of experiment, survey or PID.

Stop MS/MS:

MS/MS stops set by accumulated threshold or by time out.

Scan rate for MS/MS:

Parameter	Description
Scan Time	The dwell time of each MS/MS scan.

Parameter	Description
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

Peak Detection tab

To prevent the selection of more than one isotopic peak from a single chromatogram, you specify a peak detection window. The window's size determines which peak of an isotopic cluster is marked as a potential candidate for an MS/MS experiment. Using a small peak-detection window, 0.5 Da in the case of a doubly charged species, you can select all the peaks in a cluster. If you increase the peak window, only the most abundant peak within the window is selected for MS/MS examination. For peptides, the typical values for a peak detection window are 2 to 3 Da.

! **Notice:** These peak detection parameters must be used only by experienced operators with the necessary training and experience. Others should adopt the suggested settings below.

Suggested peak detection parameters:

Parameter	Description
Apply no criteria other than intensity (as specified on the MS Survey tab) for peak selection	The default setting. (Peak selection occurs by intensity only.)
Charge State Peak selection	When selected, the Charge State Peak detection parameters find peaks of interest. See below: Charge state peak selection parameters
Isotope Pattern selection	When selected, isotope pattern detection algorithms find peaks of interest in the MS scan. See below: Isotopic Pattern
Deisotope Peak selection	When selected, a deisotope peak filter is added to the isotopic pattern algorithm. See below: Deisotope peak selection

Charge state peak selection parameters

The charge state peak selection parameters are available only when you select Charge State Peak selection. They are typically used for multiply charged proteolytic peptide mixtures.

Parameter	Description
Select Charge States of interest	Selects the required charge states. (Selected charge states are indicated by depressed buttons)
Use advanced options for...	When selected, enables the advanced parameters. When cleared, the default values apply. Number of Components: Specifies the maximum number of components you allow identified from a single MS survey spectrum. Typically, a value such as 60 is used. Tolerance Window +/-: Specifies the tolerance window for charge-state determination. Typically, 0.2 Da is used. Peak Extraction Window: Recommendation: Specify 3 Da for peptide work.

Isotopic Pattern

Specify an isotope pattern in the Isotopic Pattern box, to only select components with the correct isotopic pattern for an MS/MS experiment.

Tip: Selecting isotopic pattern parameters is particularly useful when the targeted precursors exhibit a specific distribution, such as that caused by chlorinated or brominated molecules.

Isotopic pattern parameters:

Parameter	Description
1st Mass Difference	Specifies the first mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the first mass difference.
2nd Mass Difference	Specifies the second mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the second mass difference

Detection parameters

The detection parameters provide additional thresholds, ratios, and tolerances for detecting the presence of the specified isotope pattern and subsequently switching to perform an MS/MS experiment.

Tolerances:

Parameter	Description
Intensity above	Considers peaks where intensity is above the specified value.
Mass within +/-	Considers peaks where mass is within the specified range.
Only consider peaks within a ratio of	When selected, enables the Intensity Ratio Tolerance parameter.
Intensity Ratio Tolerance	Considers peaks within the specified intensity ratio tolerance range.

Deisotope Peak Selection

Selecting the Deisotope Peak Selection check box extracts an area of a spectrum and subjects it to the charge state recognition function. The extraction window must be wide enough to extract all peaks in a cluster, otherwise a false answer results. You set the minimum and maximum number of charges considered by selecting the appropriate charge state. For proteolytic peptides, 2X-, 3X- and 4X- charged components are normally considered. The charge state algorithm calculates theoretical isotope distributions for the identified component mass for all charge states considered.

The Peak Extraction window must encompass all the peaks in a cluster. When the extraction window is too small, the deisotope function returns an incorrect mass.

Deisotope peak selection parameters:

Parameter	Description
Use advanced options for	When selected, enables the advanced parameters. When cleared, the default values are used.
Tolerance Window +/-	Specifies the tolerance window for finding the best fit for the isotopes to determine their charge state. Usually set at half the extraction window, 3.5 Da.
Peak Extraction Window	For Deisotope Peak Selection, a typical value is 7 Da (double the overall peak cluster size).

[Back to top](#)

Exclude tab

On the Exclude tab you specify predefined masses that you want to exclude from MS/MS analysis. These are the two types of peak exclusion:

- Dynamic.
- Fixed.

Dynamic peak exclusion:

Parameter	Description
Enable real time exclusion of masses from MS/MS	When selected, enables the real-time exclusion option.
Acquire once then always exclude for the rest of the acquisition	When selected, automatically excludes masses from MS/MS scans after their first acquisition.(Used for infusion type experiments.)
Acquire and then exclude for (seconds)	When selected, allows excluded masses for a specific amount of time. Using the Exclude Mass dialog box, you can specify the amount of time for which a mass is excluded.(Used for LCMS type experiments.)
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Exclude list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Exclude list.

Fixed peak exclusion:

Parameter	Description
Exclude from range	When selected, you can specify mass ranges so that any masses detected within the ranges are excluded. Separate individual masses with commas. Use an underscore to denote a range of masses: for example, 510, 520, 550_600, 700. Tip: Excluding all masses in this way works for only a small number of components. For larger numbers use "Exclude from file".
Exclude from file	When selected, specifies the file name of an excluded masses file. Tip: This parameter provides a useful way of excluding large numbers of components. For example, you can use it to exclude components or impurities that are not permanently present. The list is also a file editor. Use the Save and Save As buttons to create files of excluded components.

To edit the Exclude list:

1. In the Fixed Peak Exclusion frame, click Add.

Tip: To edit an existing entry in the Exclude list, double-click the entry.

2. In the Exclude Mass dialog box, specify settings, as required.

The parameters in the Exclude Mass dialog box are described below.

Parameter	Description
Add/Modify	Exclude Mass (m/z) specifies the mass to exclude from the acquisition.


Parameter	Description
Exclude Time	When selected, specifies the excluded time, in minutes. If a mass is within the Exclude tab Retention Time Window value, it is excluded from the acquisition.

[Back to top](#)

Include tab

In the Include tab, you can predefine m/z values for including in the MS/MS scan. Then you can select ions from the available list, or select them preferentially.

Include File List frame, with file loaded:

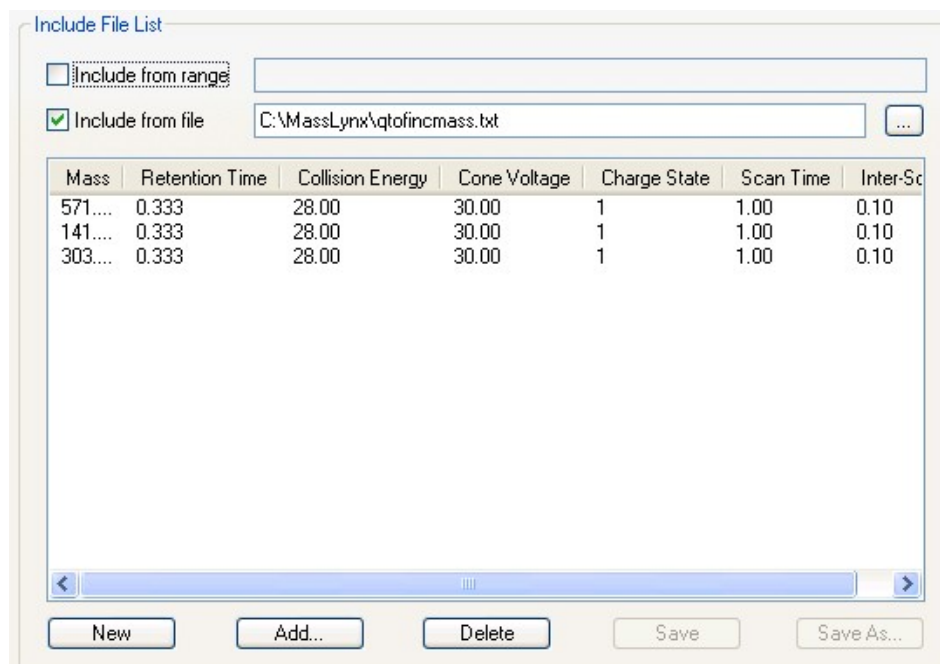

Include Mass parameters:

Parameter	Description
Enable MS/MS peak selection using an include list	Select to activate the Include list.
Only select peaks that are in the include list	Include peaks in the list and exclude all others.
Preferentially select peaks that are on the include list	When selected, masses on the Include list have priority over any other mass detected for MS/MS
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Include list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Include list.

Include File list parameters:

Parameter	Description
-----------	-------------

Parameter	Description
Include from range	<p>When selected, you can specify mass ranges so that any masses detected within the ranges are included. Separate individual masses with commas. Use an underscore to denote a range of masses, for example, 510, 520, 550_600, 700.</p> <p>Tip: Including all masses in this way works for only a small number of components. For larger numbers use Include from file.</p>
Include from file	<p>When selected, specifies the name of a file containing included masses file.</p> <p>Tip: This parameter provides a useful way of including large numbers of components.</p> <p>The list is also a file editor. Use the Save and Save As buttons to create files containing included mass. When you add an entry, you can also enter additional configuration values.</p>

To edit the Include File List:

1. Click Add.

Tip: To edit an existing entry in the Include list, double-click the entry.

2. In the Include Masses dialog box, select the check box for the required option, and enter a value in the adjacent text box.

Rule: The text box becomes available only when you select an option.

Include Masses dialog box:

Parameter	Description
Include Mass	Specifies the mass added to the mass list.
Fragmentation mode	Choose the fragmentation mode used for the mass.
Retention Time	Specifies the associated retention time, in minutes. If the specified mass elutes at the specified retention time and is within the Include tab's retention time (+/-) value, the mass is considered to be on the Include File list. Masses eluting from a column are included as masses of interest.
Collision Energy	Specifies the collision energy used to fragment the detected mass for the MS/MS scan when you choose the CID fragmentation mode.
Cone Voltage	Specifies the cone voltage value applied during MS/MS acquisitions.
Charge State	Include the mass by the charge state.
Scan Time	The scan time, in seconds, used during MS/MS scans for the detected mass.

[Back to top](#)

Collision Energy

Using the Collision Energy tab, you can:

- Use Trap CE or Transfer CE
- Use Collision Energy Profile
- Use Charge State Recognition
- Use Collision Energy Ramp

Use Charge State Recognition

Allows for the software to calculate the required collision energy for each MS/MS mass, depending on its charge state.

Use Collision Energy Profile

This mode sets the collision energy according to the charge state z and m/z of the precursor as determined by the software. Collision energy files for charge states from 1 to 6 may be created and selected. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so each m/z across the mass range will have different energies applied.

Example: When a mass of interest is detected in an MS experiment, its charge state is calculated. Using the mass and its charge state, you can obtain a collision energy value from the entered charge state table. The value is then used in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files in the CE Control tab, click Modify, and select the appropriate tab. When you have finished, click Save and OK.

Use Collision Energy Ramp

The collision energy (CE) ramps between start and end voltages. The voltages are specified for low-mass (LM) and high-mass (HM), with intermediate MS/MS scans using linearly interpolated values. The low and high masses are set on the MS Survey tab.

[Back to top](#)

Transfer CE Control tab

On the Transfer CE Control tab, you can select these options:

- Use Fixed Collision Energy
- Use Collision Energy Ramp
- Use Charge State Recognition

Use Collision Energy Ramp

For MS/MS scans, enter the start and end energy values.

Use Charge State Recognition

The Use Charge State Recognition option sets the collision energy according to the charge state z and m/z of the precursor, as determined by the software. Collision energy files for charge states 1 through 6 can be created and selected. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so different energies are associated with each m/z across the mass range.

Example: When a mass of interest is detected in an MS experiment, its charge state is calculated. Using the mass and its charge state, you can obtain a collision energy value from the charge state specified in the table. The value is then used in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files in the CE Control tab, click Modify, and select the appropriate tab. When you finish, click Save and OK.

Mobility

Data Collection

Instrument Conditions:

Parameter	Description
Override IMS Wave velocity specified in tune file	Select this option to input a velocity in m/s
Override Transfer Wave velocity specified in tune file	Select this option to input a velocity in m/s

MS Survey:

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave setting is used.
Apply Rule File for Charge State / Drift Time Stripping	Select this option to filter the acquired data using a rule file (.rul).
Store Stripped Data	Save to raw original IMS data without Rule File stripping

MS/MS:

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave setting is used.
Use Precursor Drift Times Only	This option automatically applies a rule to MSMS that filters out drift times when the precursor ion is not present
Store Stripped Data	Save to raw original IMS data without Dynamic Rule stripping.
Sensitivity	<p>Normal:</p> <ul style="list-style-type: none"> Select Normal for normal sensitivity. <p>Target Enhancement:</p> <ul style="list-style-type: none"> Select Target Enhancement (Enhanced Duty Cycle) to boost the sensitivity for the selected masses. For this option you have only 5000.0 or 100000.0 end mass selection options. <p>Wideband Enhancement:</p> <ul style="list-style-type: none"> Select to use settings in the HDC lookup file on the MSMS data.

[Back to top](#)

Related Topics

- [Add lockmass correction to methods](#)

Fast_DDA_method

- [MassLynx > MS Method > Fast DDA](#)

SYNAPT G2-Si features improvements to the FastDDA method.

FastDDA now features a smarter algorithm that provides more efficient timing of MS/MS acquisition, leading to more efficient LC/MS/MS and a greater number of switches per unit time.

FastDDA uses CIDMS/MS data and the acquisition is based on a variety of user-defined switching criteria. FastDDA uses up to 30 precursors, selected from the MS 'survey' scan, for MS/MS. FastDDA has the same intelligent decision making nature of 'Survey' DDA (for ETD and ETD/CID) and uses an embedded PC, not the host PC. MS to MS/MS switching can be done faster and this helps capture MS/MS data more successfully from lower intensity components or those components with very narrow chromatographic peaks.


Note: FastDDA supports CID fragmentation but does not support HDMS mode of operation or ETD fragmentation.

The experiment assesses the ion selection criteria for an MS/MS acquisition in the following order:

1. Threshold – MS Survey tab
2. Peak Selection – Peak Detection tab
3. Exclude list – Exclude tab
4. Include list – Include tab

Once the peaks of interest are selected, the MS/MS cycles through them until the stop criteria are met (MS/MS tab).

To create a survey experiment:

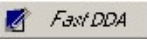1. From the main MassLynx window, click MS Method.
2. Click .
3. In the Function dialog box, specify parameter settings on the following tabs:
 - Acquisition
 - MS Survey
 - MS/MS
 - Peak Detection
 - Exclude
 - Include
 - Collision Energy
4. Click OK.
5. In the Experiment Setup window, click the LockSpray button.

Tips:

- Where your source enclosure supports it, you can use a lock-mass. The LockSpray sources support lock-mass through use of an electrospray probe and baffle. The settings described below are necessary for acquiring lock-mass data from any source that supports lock-mass acquisition.
- The LockSpray source runs at a flow rate that gives a stable spray, typically in the range 5 to 20 $\mu\text{L}/\text{min}$. For extra capacity, for instance when you must run sample lists overnight or over a weekend, you can use an external reservoir.

For your reference, here are the calculated times for a 30 mL reservoir at different lock-mass flow rates:

Flow rate ($\mu\text{L}/\text{min}$)	Running time (hours)
5	100

Flow rate ($\mu\text{L}/\text{min}$)	Running time (hours)
20	25

6. Choose to acquire data with or without correction applied.

Tip: Some data applications require correction data to be separate.

7. Select a reference compound.

8. Specify the lock-mass acquisition settings.

Requirement: Also complete these values for APGC sources.

LockSpray Settings:

Parameter	Description	Typical short-run	Typical long-run
Scan Time	The dwell time for the lock-mass acquisition	0.3 s	1.0 s
Interval	The time between lock-mass acquisitions	20 s	60 s
Scans to average	Specifies the size of the rolling average used for the lock-mass correction	3	1
Mass Window +/-	Adjusts the window size to locate the correct lock-mass peak in the LockSpray spectrum	0.3 Da	0.3 Da

9. Click OK.

10. Double click on the experiment bar.

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Acquisition settings:

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Select the polarity, and analyser mode. Improved Sensitivity: When Improved Sensitivity is selected, MSMS data is acquired in Sensitivity mode. This only applies in Resolution mode

[Back to top](#)

MS Survey tab

When selected, acquire regularly spaced surveys at specified intervals.

Using the MS Survey tab, you can configure the following parameters:

Parameter	Description
-----------	-------------

Parameter	Description
MS Survey Da Range	The mass range (m/z) over which to acquire data
MS Survey Switching Threshold	Set the intensity of individual ion
MS survey scanning conditions	The conditions for which to acquire data

MS survey Da range:

Parameter	Description
Start	The lowest mass in the spectrum
End	The highest mass in the spectrum

MS survey switching threshold:

Parameter	Description
Threshold	Sets the intensity/time when switch to MS/MS acquisition will occur

MS survey scanning conditions:

Parameter	Description
Scan Time	Specifies the duration over which TOF ion detections accumulate in each individual spectrum presented to MassLynx.
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

MS/MS tab

In the MS/MS tab, you can configure these parameters:

Parameter	Description
MS/MS Da range	The acquisition mass range where the range must differ from the MS acquisition mass range
Select (1-30)	Maximum number of precursor ions from each MS survey scan.
Scan rate for MS/MS	The dwell time for each MS/MS acquisition and the data format
Stop MS/MS	The conditions at which the MS/MS acquisition stops and the acquisition switches back to MS survey
Accumulated TIC threshold	MS/MS of a given precursor will stop when the TIC threshold is reached in the MS/MS data and unless an exceptional maximum time has changed

MS/MS Da range:

Parameter	Description
Low mass	The lowest mass in the MS/MS spectrum.
High mass	The highest mass in the MS/MS spectrum.
Max No. of Ions.	Select 1 - 8

Select (1-30):

Select from 1 to 30 ions as the maximum number of precursor ions for MS/MS from each MS survey scan.

Scan rate for MS/MS:

Parameter	Description
Scan Time	The dwell time of each MS/MS scan.
Data Format	The data format used to save the scan, centroid or continuum.

Stop MS/MS:

MS/MS spectra are acquired until either of the following two stop criteria are met:

Parameter	Description
Accumulated TIC Threshold	<p>Stop the MS/MS scan when the TIC or BPI falls below or rises above the specified intensity (normalized to one second).</p> <p>To estimate the required intensity normalized to one second (intensity/sec), read the intensity from the chromatogram, and divide it by the scan time.</p> <p>Tip: Where your selected scan time is not a whole second, hover over the units to display an active tool tip showing the corresponding intensity per scan.</p>
Time	<p>Stop the MS/MS scan when the specified number of seconds has elapsed.</p> <p>Tip: Calculate the optimum time by considering the chromatographic peak width and the typical number of ions selected by each survey scan.</p>

[Back to top](#)

Peak Detection tab

Set a peak detection window to prevent the selection of more than one isotopic peak from a single chromatogram. The size of this window determines which peak of an isotopic cluster is marked as a potential candidate for an MS/MS experiment. In a small peak detection window—0.5 Da, in the case of a doubly charged species—you can select all the peaks in a cluster. If you increase the peak window, only the most abundant peak within the window is selected for MS/MS scanning. Typical values for a peak detection window for peptides are 2 to 3 Da.



Notice: Waters advises that only experienced operators with the necessary training and experience use these peak detection parameters. Others must adopt the suggested settings, below.

Suggested peak detection parameters:

Parameter	Description
Apply no criteria other than intensity (as specified on the MS Survey tab) for peak selection	The default setting. Peak selection is by intensity only.
Charge State Peak selection	<p>When selected, the Charge State Peak selection parameters are used to find peaks of interest.</p> <p>See below: Charge state peak selection parameters</p>
Isotope Pattern selection	<p>When selected, isotope pattern detection algorithms are used to find peaks of interest in the MS scan.</p> <p>See below: Isotopic Pattern</p>

Parameter	Description
Deisotope Peak selection	When selected, a deisotope peak filter is added to the isotopic pattern algorithm. See below: Deisotope peak selection

Charge state peak selection parameters

The Charge State Peak selection parameters are typically used for multiply charged proteolytic peptide mixtures.

Parameter	Description
Select Charge States of interest	Selects the required charge states. (selected charge states are indicated by depressed buttons)
Number of Components	Specifies the maximum number of components you allow identified from a single MS survey spectrum. Typically a value such as 60 is used.
Tolerance Window +/-	Specifies the tolerance window for charge state determination. Typical value of 0.2 Da is used.
Peak Extraction Window	For peptide work a value of 3 Da is recommended.

Isotopic Pattern

Specify an isotope pattern in the Isotopic Pattern box, to only select components with the correct isotopic pattern for an MS/MS experiment.

Tip: Selecting isotopic pattern parameters is particularly useful when the targeted precursors exhibit a specific distribution, such as that caused by chlorinated or brominated molecules.

Isotopic pattern criteria:

Parameter	Description
1st Mass Difference	Specifies the first mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the first mass difference.
2nd Mass Difference	Specifies the second mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the second mass difference

Detection parameters

The detection parameters allow the use of tolerances to detect the presence of the specified isotope pattern.

Tolerances:

Parameter	Description
BPI above	Considers peaks where intensity is above the specified value.
Mass within +/-	Considers peaks where mass is within the specified range.
Only consider peaks within a ratio of	When selected, enables the Intensity Ratio Tolerance parameter.
Intensity Ratio Tolerance	Considers peaks within the specified intensity ratio tolerance range.

Deisotope Peak Selection

When you select Deisotope Peak Selection, an area of a spectrum is extracted and passed to the Charge State Recognition algorithm. The extraction window must be wide enough to extract all the peaks within the cluster, otherwise a false answer results. You set the minimum and maximum number of charges considered by selecting the appropriate charge state. For proteolytic peptides the charged components, 2+, 3+, and 4+, are normally considered. The charge state algorithm first calculates theoretical isotope distributions for the identified component mass for all charge states that must be considered.

In the case of Deisotope Peak Selection, the Peak Extraction window must reveal all the peaks in the cluster. If the extraction window is too small, the deisotope algorithm returns an incorrect mass.

Deisotope peak selection parameters:

Parameter	Description
Use advanced options for	When selected, enables the advanced parameters. When cleared, the default values are used.
Tolerance Window +/-	Specifies the tolerance window for finding the best fit for the isotopes to determine their charge state. Usually set at half the extraction window, 3.5 Da.
Peak Extraction Window	For Deisotope Peak Selection, a typical value is 7 Da (double the overall peak cluster size).

[Back to top](#)

Exclude tab

On the Exclude tab you specify predefined masses that you want to exclude from MS/MS analysis. These are the two types of peak exclusion:

- Dynamic.
- Fixed.

Dynamic peak exclusion:

Parameter	Description
Enable real time exclusion of masses from MS/MS	When selected, enables the real-time exclusion option.
Acquire once then always exclude for the rest of the acquisition	When selected, automatically excludes masses from MS/MS scans after their first acquisition.(Used for infusion type experiments.)
Acquire and then exclude for (seconds)	When selected, allows excluded masses for a specific amount of time. Using the Exclude Mass dialog box, you can specify the amount of time for which a mass is excluded.(Used for LCMS type experiments.)
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Exclude list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Exclude list.

Fixed peak exclusion:

Tip: The list accepts up to 100 entries. Further entries are not used.

Parameter	Description
-----------	-------------

Parameter	Description
Exclude from range	When selected, you can specify mass ranges so that any masses detected within the ranges are excluded. Separate individual masses with commas. Use an underscore to denote a range of masses? for example, 510, 520, 550_600, 700. Tip: Excluding all masses in this way works for only a small number of components. For larger numbers use "Exclude from file".
Exclude from file	When selected, the list of masses becomes a file editor. Use the Save and Save As buttons to create files or click browse and load a preexisting file. When you add an entry, you can also optionally enter further configuration values.

To edit the Exclude list:

1. In the Fixed Peak Exclusion frame select the Exclude from file check box.
2. Select Add to open the Exclude Mass dialog box.

Requirement: To edit an existing entry in the Exclude list, double-click the entry.

3. Enter the values as required.

Exclude mass dialog box parameters:

Parameter	Description
Add/Modify	Exclude Mass (m/z) specifies the mass to exclude from the acquisition.
Exclude Time	When selected, specifies the excluded time, in minutes. If a mass is within the Exclude tab Retention Time Window value, it is excluded from the acquisition.

[Back to top](#)

Include tab

In the Include tab, you can predefine m/z values for including in the MS/MS scan. Then you can select ions from the available list, or select them preferentially.

Include Mass:

Parameter	Description
Enable MS/MS peak selection using an include list	Select to activate the Include list.
Only select peaks that are in the include list	Include peaks in the list and exclude all others.
Preferentially select peaks that are on the include list	When selected, masses on the Include list have priority over any other mass detected for MS/MS
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Include list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Include list.

Include File List:

Tip: The list accepts up to 100 entries. Further entries are not analyzed.

Parameter	Description
-----------	-------------

Parameter	Description
Include from range	When selected, you can specify mass ranges so that any masses detected within the ranges are included. Separate individual masses with commas. Use an underscore to denote a range of masses, for example, 510, 520, 550_600, 700. Tip: Including all masses in this way works for only a small number of components. For larger numbers use the "Include from file" parameters.
Include from file	When selected, the list of masses becomes a file editor. Use the Save and Save As buttons to create files, or click browse, and load a preexisting file. When you add an entry, you can also optionally enter further configuration values.

To edit the Include File List:

1. Check the box for Include from file.
2. Click Add to open the Include Masses dialog box

Tip: To edit an existing entry in the Include list, double-click it.
3. Select the check box for the required option, and enter the value in the adjacent text box.

Rule: The text box becomes available only when you select an option.

Include Masses parameters:

Parameter	Description
Include Mass	Specifies the mass added to the mass list.
Retention Time	Specifies the associated retention time, in minutes. If the specified mass elutes at the specified retention time and is within the Include tab's retention time (+/-) value, the mass is considered to be on the Include list. Masses eluting from a column are included as masses of interest.
Collision Energy	Specifies the collision energy used to fragment the detected mass for the MS/MS scan.
Cone Voltage	Specifies the cone voltage value applied during MS/MS acquisitions.
Charge State	Include the mass by the charge state.
Scan Time	The scan time, in seconds, used during MS/MS scans for the detected mass.

[Back to top](#)

Collision Energy tab

Using the Collision Energy tab, you can:

- Use iTRAQ[®] Mode
- Use Trap CE
- Use Transfer CE
- Use Collision Energy Profile
- Use MSMS Collision Energy

- Use Collision Energy Map

Use iTRAQ® Mode

Checking the iTRAQ Mode button allows the Elevated collision energy to be adjusted.

Use Trap CE or Use Transfer CE

Checking Use Trap CE brings up Use Change State Recognition button, Use Collision Energy Ramp button, Use Collision Energy Profile button and the Collision Energy Profile window.

Checking Use Transfer CE brings up Use a Fixed Collision Energy Value (in Volts) and Use Collision Energy Ramp buttons.

Use Collision Energy Profile

The collision energy profile provides a range of different collision energies, depending on the m/z value of the selected precursor. In addition, you can use between 1 and 5 collision energies for each precursor m/z to attempt to produce optimal fragmentation.

For ranges not specified in the profile, this option uses the collision energy taken from the Tune window unless you override it and specify an Override Collision Energy.

Use Charge State Recognition

This mode sets the collision energy according to the charge state z and m/z of the precursor as determined by the software. You can create and select collision energy files for charge states from 1 to 6. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so each m/z across the mass range will have different energies applied.

Example: When the software detects a mass of interest in an MS experiment, it calculates the mass's charge state. Using the mass and its charge state, you can obtain a collision energy value from the charge state table, and then use the charge state in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files, in the CE Control tab, click Modify, and select the appropriate tab. When you have finished, click Save and OK.

Use Collision Energy Ramp

The collision energy (CE) ramps between start and end voltages. The voltages are specified for low-mass (LM) and high-mass (HM), with intermediate MS/MS scans using linearly interpolated values. The low and high masses are set on the MS Survey tab.

[Back to top](#)

Related Topics

- [Prepare the instrument for use](#)

DDA Survey method (Data directed analysis for ETD)

→ [MassLynx](#) > [MS Method](#) > [MS Survey](#)

This method is particularly useful for generation of high quality sequence information from modified (phospho- or glyco-) peptide species to enable identification of site of attachment as in addition to the core peptide sequence. 'Survey' is a DDA method where MS/MS data is generated from precursors of interest using Electron Transfer Dissociation (ETD) alone or in combination with CID, and the acquisition based on a variety of user-defined switching criteria. The method allows the user to conduct a 'Survey' (MS) scan to identify (up to 8) precursors for MS/MS analysis based on the ions detected in that survey scan. The 'real time' decision making for switching from

MS to MS/MS is made on the host PC therefore the switching time differs from that of FastDDA. The available MS/MS analysis time (between MS survey scans) is user defined and can be distributed between the target precursors.


Recommendation: This topic describes the standard MS Survey experiment. However, we recommend that you use FastDDA where possible. FastDDA is a data directed survey experiment that is controlled on the instrument rather than in the MassLynx software. This makes the analysis and switching more rapid and avoids missing narrow chromatographic peaks.

The criteria for selecting ions for an MS/MS acquisition are assessed in the following order:

1. Threshold – [MS Survey tab](#)
2. Peak Selection – [Peak Detection tab](#)
3. Exclude list – [Exclude tab](#)
4. Include list – [Include tab](#)

Once the peaks of interest are selected, the MS/MS cycles through them until the Stop MS/MS criteria are met (MS/MS tab).

To create a survey experiment:

1. From the main MassLynx window, click MS Method.
2. Click Survey  Survey, to add a survey experiment.
3. Specify your parameters on the following tabs:
 - [Acquisition](#)
 - [MS Survey](#)
 - [MS/MS](#)
 - [Peak Detection](#)
 - [Exclude](#)
 - [Include](#)
 - [CE Control](#)
 - [Cone Voltage](#)
4. Click OK
5. In the Experiment Setup window, click the LockSpray button.

Tips:

- Where your source enclosure supports it, you can use a lock-mass. The LockSpray sources support lock-mass through use of an electrospray probe and baffle. The APGC source supports lock-mass through the use of reference vials in the source. The settings described below are necessary for acquiring lock-mass data from any source that supports lock-mass acquisition.
- The LockSpray runs at a flow rate that gives a stable spray, typically in the range 5 to 20 $\mu\text{L}/\text{min}$. For extra capacity, for instance when you need to run sample lists overnight or over a weekend, you can use an external reservoir. For your reference, here are the calculated times for a 30 mL reservoir at different lock-mass flow rates:

Flow rate ($\mu\text{L}/\text{min}$)	Running time (hours)
5	100
10	50
20	25

6. Choose to acquire data with or without correction applied.

Tip: Some data applications need correction data to be separate.

7. Select a reference compound.

8. Specify the lock-mass acquisition settings.

Requirement: Also complete these values for APGC sources.

LockSpray Settings:

Parameter	Description	Typical short-run	Typical long-run
Scan Time	The dwell time for the lock-mass acquisition	0.3 s	1.0 s
Interval	The time between lock-mass acquisitions	20 s	60 s
Scans to average	Specifies the size of the rolling average used for the lock-mass correction	3	1
Mass Window +/-	Adjusts the window size to locate the correct lock-mass peak in the LockSpray spectrum	0.3 Da	0.3 Da

9. Click OK.

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Select the polarity, and dynamic range.

[Back to top](#)

MS Survey tab

Using the MS Survey tab, you can configure the following parameters:

- MS survey Da range – the mass range (m/z) over which to acquire data.
- Stopping the MS survey – the conditions at which the survey scan stops and the acquisition switches to MS/MS mode.
- MS survey scanning conditions – the conditions for which to acquire data.

MS survey Da range

Parameter	Description
Start	The start m/z value.
End	The end m/z value.
Discard MS Survey Scans that do not contain ions selected for MS/MS analysis	Where no MS/MS scans are performed discard the MS scan data. This helps to reduce file size.

Stopping the MS survey:

Parameter	Description
TIC rises above	When selected, the MS scan stops and the MS/MS scan starts when the total ion current (TIC) of the scan rises above the specified threshold value.
Intensity of individual ion rising above	Threshold above which the software recognizes an ion and selection passes to the Peak Detection tab.
Threshold	The value at which the switch to MS/MS acquisition is made.

MS survey scanning conditions:

Parameter	Description
Scan Time	Specifies the duration over which TOF ion-detections accumulate in each individual spectrum presented to MassLynx.
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

MS/MS tab

In the MS/MS tab, you can configure these parameters:

- MS/MS m/z range – the m/z range over which to acquire data.
- Stopping MS/MS – the conditions at which the MS/MS acquisition stops and the acquisition switches back to MS survey.
- Scan rates and instrument conditions for MS/MS – the conditions under which you will acquire the data.

Tip: The parameter values you specify to stop an acquisition depend upon the type of experiment: Survey or PID.

MS/MS m/z range:

Parameter	Description
Acquire MS/MS over the range	
Start	The m/z at which the MS/MS scan starts. The lowest mass in the MS/MS spectrum.
End	The m/z at which the MS/MS scan stops. The highest mass in the MS/MS spectrum.
Maximum number of ions that can be selected for MS/MS from a single MS survey scan.	
Select (1-8)	Specifies the maximum number of individual precursors that can be selected from a single MS survey.

Stop MS/MS:

MS/MS stops when the first of these criteria are met:

Parameter	Description
Return to MS survey scan after each MSMS ion acquires	Specify the number of scans.
Return to MS Survey scan when the following criteria are met.	Select a reason or time to halt MS/MS acquisition.

Scan rate for MS/MS:

Parameter	Description
Scan Time	The dwell time of each MS/MS scan.
Data Format	The data format used to save the scan, centroid or continuum.
Tune file for MS/MS	When selected, the cone voltage defined in the tune file in the adjacent box is used. Type the file name in the box, or browse for a file.
MS/MS lock-mass	You can optionally switch the tuning parameters during the MS/MS scan and select a lock-mass that is more appropriate for the altered tuning.

[Back to top](#)

Peak Detection tab

To prevent the selection of more than one isotopic peak from a single chromatogram, a peak detection window is set. The size of this window determines which peak of an isotopic cluster is marked as a potential candidate for an MS/MS experiment. Using a small peak detection window—0.5 Da, in the case of a doubly charged species—all the peaks in a cluster can be individually selected. If you increase the peak window, only the most abundant one within the window is selected for MS/MS. Typical values for a peak detection window for peptides are 2 to 3 Da.

! **Notice:** These peak detection parameters must be used only by experienced operators with the necessary training and experience. Others should adopt the suggested settings below.

Parameter	Description
Apply no criteria other than intensity (as specified on the MS Survey tab) for peak selection	The default setting. Peak selection is by intensity only.
Charge State Peak selection	When selected, the Charge State Peak detection parameters are used to find peaks of interest. See below: Charge state peak selection parameters
Isotope Pattern selection	When selected, isotope pattern detection algorithms are used to find peaks of interest in the MS scan. See below: Isotopic Pattern
Deisotope Peak selection	When selected, a deisotope peak filter is added to the isotopic pattern algorithm. See below: Deisotope peak selection
Peak Detection Window	The window size used to detect peaks. Typical values for peptides are 2 to 3 Da.

Charge State Peak selection parameters

The Charge State Peak selection parameters are available only when the Charge State Peak selection is selected. These are typically used for multiply charged proteolytic peptide mixtures.

Parameter	Description
-----------	-------------

Parameter	Description
Select Charge States of interest	Selects the required charge states. (selected charge states are indicated by depressed buttons)
Use advanced options for	When selected, enables the advanced parameters. When cleared, the default values are used.
Number of Components	Specifies the maximum number of components you allow identified from a single MS survey spectrum. Typically a value such as 60 is used.
Tolerance Window +/-	Specifies the tolerance window for charge state determination. Typical value of 0.2 Da is used.
Peak Extraction Window	For peptide work a value of 3 Da is recommended.

Isotopic Pattern

If you specify a specific isotope pattern in the Isotopic Pattern box, only components fulfilling these criteria are selected for an MS/MS experiment.

Tip: Selecting isotopic pattern parameters is particularly useful when the targeted precursors exhibit a specific distribution, such as that caused by chlorinated or brominated molecules.

Parameter	Description
1st Mass Difference	Specifies the first mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the first mass difference.
2nd Mass Difference	Specifies the second mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the second mass difference

Detection parameters

The detection parameters allow the tolerances of the following parameters to be used to detect the presence of the specified isotope pattern and subsequently switch to perform an MS/MS experiment.

Parameter	Description
Intensity above	Considers peaks that have an intensity above the specified value.
Mass within +/-	Considers peaks that have a mass within the specified range.
Only consider peaks within a ratio of	When selected, enables the Intensity Ratio Tolerance parameter.
Intensity Ratio Tolerance	Considers peaks within the specified intensity ratio tolerance range.

Deisotope Peak Selection

When you select Deisotope Peak Selection, an area of a spectrum is extracted and passed to the Charge State Recognition routine. The extraction window needs to be wide enough to extract all the peaks within the cluster, or a false answer results. You set minimum and maximum number of charges to be considered by selecting the appropriate charge state. For proteolytic peptides, 2X-, 3X- and 4X- charged components are normally considered. The charge state algorithm first calculates theoretical isotope distributions for the identified component mass for all charge states that must be considered.

In the case of Deisotope Peak Selection the Peak Extraction window must cover all the peaks in the cluster. If the extraction window is too small, the deisotope routine would return an incorrect mass.

Parameter	Description
Use advanced options for	When selected, enables the advanced parameters. When cleared, the default values are used.
Tolerance Window +/-	Specifies the tolerance window for finding the best fit for the isotopes to determine their charge state. Usually set at half the extraction window, 3.5 Da.
Peak Extraction Window	For Deisotope Peak Selection, a typical value is 7 Da (double the overall peak cluster size).

[Back to top](#)

Exclude tab

On the Exclude tab you specify predefined masses that you want to exclude from MS/MS analysis. These are the two types of peak exclusion:

- Dynamic.
- Fixed.

Dynamic peak exclusion

Parameter	Description
Enable real time exclusion of masses from MS/MS	When selected, enables the real-time exclusion option.
Acquire once then always exclude for the rest of the acquisition	When selected, automatically excludes masses from MS/MS scans after their first acquisition. (Used for infusion type experiments.)
Acquire and then exclude for (seconds)	When selected, allows excluded masses for a specific amount of time. Using the Exclude Mass dialog box, you can specify the amount of time for which a mass is excluded. (Used for LCMS type experiments.)
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Exclude list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Exclude list.

Fixed peak exclusion

Parameter	Description
Exclude from range	When selected, you can specify mass ranges so that any masses detected within the ranges are excluded. Separate individual masses with commas. Use an underscore to denote a range of masses, for example, 510, 520, 550_600, 700. Tip: Excluding all masses in this way works for only a small number of components. For larger numbers use "Exclude from file".
Exclude from file	With this selected, the list of masses becomes a file editor. Use the Save and Save As buttons to create files or click browse and load a pre-existing file. When you add an entry, you can also optionally enter further configuration values.

To edit the Exclude list:

1. In the Fixed Peak Exclusion frame, click Add to open the Exclude Mass dialog box.

Requirement: To edit an existing entry in the Exclude list, double-click the entry.

2. Enter the values as required.

The parameters in the Exclude Mass dialog box are described below.

Exclude mass dialog box parameters:

Parameter	Description
Add/Modify	Exclude Mass (m/z) specifies the mass to exclude from the acquisition.
Exclude Time	When selected, specifies the excluded time, in minutes. If a mass is within the Exclude tab Retention Time Window value, it is excluded from the acquisition.

[Back to top](#)

Include tab

In the Include tab, you can predefine m/z values for including in the MS/MS scan. Then you can select ions from the available list or select them preferentially.

Include Mass

Parameter	Description
Enable MS/MS peak selection using an include list	Select to activate the include list.
Only select peaks that are in the include list	Include peaks in the list and exclude all others.
Preferentially select peaks that are on the include list	When selected, masses on the Include File List have priority over any other mass detected for MS/MS
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all include list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the include list.

Include File List

Parameter	Description
Include from range	When selected, you can specify mass ranges so that any masses detected within the ranges are included. Separate individual masses with commas. Use an underscore to denote a range of masses, for example, 510, 520, 550_600, 700. Tip: Including all masses in this way works for only a small number of components. For larger numbers use "Include from file".
Include from file	With this selected, the list of masses becomes a file editor. Use the Save and Save As buttons to create files or click browse and load a pre-existing file. When you add an entry, you can also optionally enter further configuration values.

To edit the Include File List:

1. Click Add to open the Include Masses dialog box.

Tip: To edit an existing entry in the Include list, double-click it.

2. Select the check box for the required option, and enter the value in the adjacent text box.

Rule: The text box becomes available only when you select an option.

The Include Masses dialog box parameters are described below.

Parameter	Description
Include Mass	Specifies the mass added to the mass list.
Retention Time	Specifies the associated retention time, in minutes. If the specified mass elutes at the specified retention time and is within the Include tab's retention time (+/-) value, the mass is considered to be on the Include File List. Masses eluting from a column are included as masses of interest.
Collision Energy	Specifies the collision energy used to fragment the detected mass for the MS/MS scan.
Cone Voltage	Specifies the cone voltage value applied during MS/MS acquisitions.
Charge State	Include the mass by the charge state.
Scan Time	The scan time, in seconds, used during MS/MS scans for the detected mass
Inter-Scan Time	The time between subsequent scans, in seconds, used during MS/MS scans for the detected mass.

[Back to top](#)

CE Control tab

Using the Collision Energy tab, you can:

- [Use Collision Energy Profile](#)
- [Use Charge State Recognition](#)
- [Use Collision Energy Ramp](#)

Use Collision Energy Profile

The collision energy profile allows a range of different collision energies to be used depending on the m/z value of the selected precursor. In addition, you can use between 1 and 5 collision energies for each precursor m/z to attempt to produce optimal fragmentation.

For ranges not specified in the profile, this option uses the collision energy taken from the Tune window unless you override it and specify an Override Collision Energy.

Use Charge State Recognition

This mode sets the collision energy according the charge state z and m/z of the precursor as determined by the software. Collision energy files for charge states from 1 to 6 may be created and selected. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so each m/z across the mass range will have different energies applied.

Example: When a mass of interest is detected in an MS experiment, its charge state is calculated. Using the mass and its charge state, you can obtain a collision energy value from the entered charge state table. The value is then used in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files in the CE Control tab, click Modify, and select the appropriate tab. When you have finished, click Save and OK.

Use Collision Energy Ramp

The collision energy (CE) ramps between start and end voltages. The voltages are specified for low-mass (LM) and high-mass (HM), with intermediate MS/MS scans using linearly interpolated values. The low and high masses are set on the MS Survey tab.

[Back to top](#)

Cone Voltage tab

Parameter	Description
Override Cone Voltage Specified in tune file	Select this override and enter your preferred cone voltage.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. The cone voltage follows a linear voltage ramp between these two values during the scan.

[Back to top](#)

Related Topics

- [Fast-DDA method](#)
- [Prepare the instrument for use](#)

Mobility DDA method

- [MassLynx > MS Method > Survey](#)

The Mobility DDA methodology utilises ion mobility in a DDA experiment to increase the sensitivity and selectivity of precursor ion selection for MS/MS. It can provide improved MS/MS quality through; 'wideband enhancement' for improved fragment ion sensitivity; or removal of interferences by precursor drift time matching.

It may be used to:

- enhance quality of MS/MS information from weak samples or
- for improved identification of compounds from mixtures of simple/medium complexity.


You assign the criteria for selecting ions for an MS/MS acquisition in this order:

1. Threshold - MS Survey tab
2. Peak Selection - Peak Detection tab
3. Exclude list - Exclude tab
4. Include list - Include tab

During the survey, once the peaks are selected, the MS/MS cycles through them until the Stop MS/MS criteria are met (MS/MS tab).

Creating a survey experiment

To create a survey experiment:

1. Click Survey  Survey .
2. Specify parameter values on these tabs:
 - Acquisition
 - MS Survey
 - MS/MS
 - Peak Detection
 - Exclude
 - Include
 - Trap CE Control
 - Transfer CE Control
 - Cone Voltage
 - Mobility

Acquisition tab

On the Acquisition tab, specify initial acquisition parameter values.

Acquisition parameters:

Parameter	Description
Acquisition Times	Enter the start and end times for the acquisition.
Source	Select the source type.
Acquisition Mode	Polarity: Select the polarity. Analyser Mode: Select the type of optics required. Dynamic range: Select either Normal or Extended.

[Back to top](#)

MS Survey tab

On the MS Survey tab, you configure these parameters:

- MS survey Da range - The mass range (m/z) for acquiring data.
- Stopping the MS survey - The conditions at which the survey scan stops and the acquisition switches to MS/MS mode.
- MS survey scanning conditions - The conditions applied when acquiring data.

MS survey Da range:

Parameter	Description
Start	The lower m/z limit.

Parameter	Description
End	The upper m/z limit.
Discard MS Survey Scans that do not contain ions selected for MS/MS analysis	Where no MS/MS scans are performed discard the MS scan data. Doing so reduces file size.

Stopping the MS survey:

Parameter	Description
TIC rises above	When selected, the MS scan stops, and the MS/MS scan starts when the total ion current (TIC) of the scan rises above the specified threshold value.
Intensity of individual ion rising above	Threshold above which the software recognizes an ion and the selection criteria specified on the Peak Detection tab are assessed.
Threshold	The value at which the switch to MS/MS acquisition occurs.

MS survey scanning conditions:

Parameter	Description
Scan Time	Specifies the duration over which TOF ion-detections accumulate in each spectrum presented to Unifi.
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

MS/MS tab

On the MS/MS tab, you can configure these parameters:

- MS/MS m/z range - The m/z range for acquiring data.
- Stopping MS/MS - The conditions at which an MS/MS acquisition stops and switches back to an MS survey acquisition.
- Scan rates and instrument conditions for MS/MS - The conditions under which MassLynx acquires data.

Tip: The parameter values you specify to stop an acquisition depend upon the type of experiment, survey or PID.

MS/MS m/z range:

Parameter	Description
Acquire MS/MS over the range:	
Start	The lower m/z limit.
End	The upper m/z limit.
Maximum number of ions selectable for MS/MS from a single MS survey scan:	
Select (1-8)	Specifies the maximum number of individual precursors selectable from a single MS survey.

Stop MS/MS:

MS/MS stops when the first of the criteria in the following table are met.

Return to MS Survey scan parameters:

Parameter	Description
Return to MS survey scan after each MSMS ion acquires	Specify the number of scans.
Return to MS Survey scan when the following criteria are met.	Select a criteria, for example elapsed time, for stopping the MS/MS acquisition.

Scan rate for MS/MS:

Parameter	Description
Scan Time	The dwell time of each MS/MS scan.
Data Format	The data format used to save the scan, centroid or continuum.

[Back to top](#)

Peak Detection tab

To prevent the selection of more than one isotopic peak from a single chromatogram, you specify a peak detection window. The window's size determines which peak of an isotopic cluster is marked as a potential candidate for an MS/MS experiment. Using a small peak-detection window, 0.5 Da in the case of a doubly charged species, you can select all the peaks in a cluster. If you increase the peak window, only the most abundant peak within the window is selected for MS/MS examination. For peptides, the typical values for a peak detection window are 2 to 3 Da.

! **Notice:** These peak detection parameters must be used only by experienced operators with the necessary training and experience. Others should adopt the suggested settings below.

Peak detection parameters:

Parameter	Description
Apply no criteria other than intensity (as specified on the MS Survey tab) for peak selection	The default setting. (Peak selection occurs by intensity only.)
Charge State Peak selection	When selected, the Charge State Peak detection parameters find peaks of interest. Note: This is the only peak selection type available for ETD fragmentation. See below: Charge state peak selection parameters
Isotope Pattern selection	When selected, isotope pattern detection algorithms find peaks of interest in the MS scan. See below: Isotopic Pattern
Deisotope Peak selection	When selected, a deisotope peak filter is added to the isotopic pattern algorithm. See below: Deisotope peak selection
Peak Detection Window	The window size used to detect peaks. Typical values for peptides are 2 to 3 Da.

Charge state peak selection parameters

The charge state peak selection parameters are available only when you select Charge State Peak selection. They are typically used for multiply charged proteolytic peptide mixtures.

Parameter	Description
Select Charge States of interest	Selects the required charge states. (Selected charge states are indicated by depressed buttons)
Use advanced options for...	When selected, enables the advanced parameters. When cleared, the default values apply. Number of Components: Specifies the maximum number of components you allow identified from a single MS survey spectrum. Typically, a value such as 60 is used. Tolerance Window +/-: Specifies the tolerance window for charge-state determination. Typically, 0.2 Da is used. Peak Extraction Window: Recommendation: Specify 3 Da for peptide work.

Isotopic pattern

If you specify an isotope pattern in the Isotopic Pattern box, only components fulfilling these criteria are selected for an MS/MS experiment.

Tip: Selecting isotopic pattern parameters is particularly useful when the targeted precursors exhibit a specific distribution, such as that caused by chlorinated or brominated molecules.

Isotopic pattern parameters:

Parameter	Description
1st Mass Difference	Specifies the first mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the first mass difference.
2nd Mass Difference	Specifies the second mass difference from the mono-isotope for isotopic selection.
Intensity Ratio	Specifies the intensity ratio for the second mass difference

Detection parameters

The detection parameters provide additional thresholds, ratios, and tolerances for detecting the presence of the specified isotope pattern and subsequently switching to perform an MS/MS experiment.

Parameter	Description
Intensity above	Considers peaks that have an intensity above the specified value.
Mass within +/-	Considers peaks that have a mass within the specified range.
Only consider peaks within a ratio of	When selected, enables the Intensity Ratio Tolerance parameter.
Intensity Ratio Tolerance	Considers peaks within the specified intensity ratio tolerance range.

Deisotope peak selection

Selecting the Deisotope Peak Selection check box extracts an area of a spectrum and subjects it to the charge state recognition function. The extraction window must be wide enough to extract all peaks in a cluster, otherwise a false answer results. You set the minimum and maximum number of charges considered by selecting the appropriate

charge state. For proteolytic peptides, 2X-, 3X- and 4X- charged components are normally considered. The charge state algorithm calculates theoretical isotope distributions for the identified component mass for all charge states considered.

The Peak Extraction window must encompass all the peaks in a cluster. When the extraction window is too small, the deisotope function returns an incorrect mass.

Deisotope peak selection parameters:

Parameter	Description
Use advanced options for	When selected, enables the advanced parameters. When cleared, the default values are used.
Tolerance Window +/-	Specifies the tolerance window for finding the best fit for the isotopes to determine their charge state. The typical value is half the extraction window, 3.5 Da.
Peak Extraction Window	For deisotope peak selection, a typical value is 7 Da (double the overall peak cluster size).

[Back to top](#)

Exclude tab

On the Exclude tab, you specify predefined masses that you want to exclude from MS/MS analysis. These are the two types of peak exclusion:

- Dynamic
- Fixed

Dynamic peak exclusion:

Parameter	Description
Enable real time exclusion of masses from MS/MS	When selected, enables the real-time exclusion option.
Acquire once then always exclude for the rest of the acquisition	When selected, automatically excludes masses from MS/MS scans after their first acquisition. (Used for infusion-type experiments.)
Acquire and then exclude for (seconds)	When selected, allows excluded masses for a specific duration. Using the Exclude Mass dialog box, you can specify the duration for which a mass is excluded. (Used for LCMS type experiments.)
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Exclude list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Exclude list.

Fixed peak exclusion:

Parameter	Description
Exclude from range	When selected, you can specify mass ranges so that any masses detected within the ranges are excluded. Separate individual masses with commas. Use an underscore to denote a range of masses: for example, 510, 520, 550_600, 700. Tip: Excluding all masses in this way works for only a small number of components. For larger numbers use "Exclude from file".

Parameter	Description
Exclude from file	<p>When selected, specifies the file name of an excluded masses file.</p> <p>Tip: This parameter provides a useful way of excluding large numbers of components. For example, you can use it to exclude components or impurities that are not permanently present.</p> <p>The list is also a file editor. Use the Save and Save As buttons to create files of excluded components.</p>

To edit the Exclude list:

1. In the Fixed Peak Exclusion frame, click Add.

Tip: To edit an existing entry in the Exclude list, double-click the entry.

2. In the Exclude Mass dialog box, specify settings, as required.

The parameters in the Exclude Mass dialog box are described below.


Parameter	Description
Add/Modify	Exclude Mass (m/z) specifies the mass to exclude from the acquisition.
Exclude Time	When selected, specifies the excluded time, in minutes. If a mass is within the retention time window, it is excluded from the acquisition.

[Back to top](#)

Include tab

On the Include tab, you can predefine m/z values for including in the MS/MS scan. Then you can select ions from the available list or select them preferentially.

Include File List frame, with file loaded:

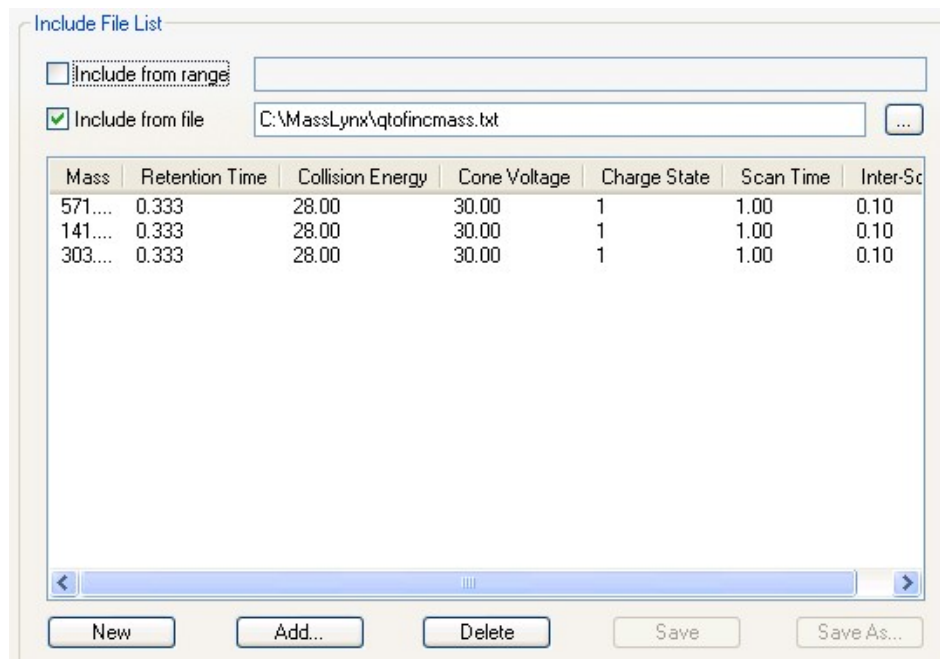

Include Mass parameters:

Parameter	Description
Enable MS/MS peak selection using an include list	Select to activate the Include list.
Only select peaks that are in the include list	Include peaks in the list, and exclude all others.
Preferentially select peaks that are on the include list	When selected, masses on the Include File list have priority over any other mass detected for MS/MS
Exclude peaks with +/- ...mDa of entries on the exclude list	Specifies a mass tolerance value that applies to all Include list masses.
Retention Time window +/- (seconds)	Specifies a tolerance to apply to the retention times in the Include list.

Include File list parameters:

Parameter	Description
Include from range	When selected, you can specify mass ranges so that any masses detected within the ranges are included. Separate individual masses with commas. Use an underscore to denote a range of masses, for example, 510, 520, 550_600, 700. Tip: Including all masses in this way works for only a small number of components. For larger numbers use Include from file.
Include from file	When selected, specifies the name of a file containing included masses file. Tip: This parameter provides a useful way of including large numbers of components. The list is also a file editor. Use the Save and Save As buttons to create files containing included mass. When you add an entry, you can also enter additional configuration values.

To edit the Include File list:

1. Click Add.

Tip: To edit an existing entry in the Include list, double-click the entry.

2. In the Include Masses dialog box, select the check box for the required option, and enter a value in the adjacent text box.

Rule: The text box becomes available only when you select an option.

Include Masses dialog box:

Parameter	Description
Include Mass	Specifies the mass added to the mass list.
Fragmentation mode	Choose the fragmentation mode used for the mass.
Retention Time	Specifies the associated retention time, in minutes. If the specified mass elutes at the specified retention time and is within the Include tab's retention time (+/-) value, the mass is considered to be on the Include File list. Masses eluting from a column are included as masses of interest.
Collision Energy	Specifies the collision energy used to fragment the detected mass for the MS/MS scan when you choose the CID fragmentation mode.
Cone Voltage	Specifies the cone voltage value applied during MS/MS acquisitions.

Parameter	Description
Charge State	Include the mass by the charge state.
Scan Time	The scan time, in seconds, used during MS/MS scans for the detected mass.
Inter-Scan Time	The interscan time, in seconds, used during MS/MS scans for the detected mass.

[Back to top](#)

Trap CE Control tab

Using the Trap CE Control tab, you can:

- [Use Collision Energy Profile](#)
- [Use fixed CE](#)
- [Use Collision Energy Ramp](#)

Use Collision Energy Profile

The collision energy profile allows a range of different collision energies to be used depending on the m/z value of the selected precursor. In addition, you can use between 1 and 5 collision energies for each precursor m/z to attempt to produce optimal fragmentation.

For ranges not specified in the profile, this option uses the collision energy taken from the Tune window unless you override it and specify an Override Collision Energy.

Use fixed CE

This mode sets the collision energy according to the charge state z and m/z of the precursor as determined by the software. Collision energy files for charge states from 1 to 6 may be created and selected. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so each m/z across the mass range will have different energies applied.

Example: When a mass of interest is detected in an MS experiment, its charge state is calculated. Using the mass and its charge state, you can obtain a collision energy value from the entered charge state table. The value is then used in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files in the CE Control tab, click Modify, and select the appropriate tab. When you have finished, click Save and OK.

Use Collision Energy Ramp

The collision energy (CE) ramps between start and end voltages. The voltages are specified for low-mass (LM) and high-mass (HM), with intermediate MS/MS scans using linearly interpolated values. The low and high masses are set on the MS Survey tab.

[Back to top](#)

Transfer CE Control tab

On the Transfer CE Control tab, you can select these options:

- [Use Fixed Collision Energy](#)
- [Use Collision Energy Ramp](#)
- [Use Charge State Recognition](#)

Rule: If ETD fragmentation mode is set on the Acquisition tab, the settings on this tab are applied only to MS/MS acquisition functions using CID fragmentation. For acquisition functions using ETD fragmentation, collision energy is not used.

Use Fixed Collision Energy

Select this option if you want to use a different collision energy value from the one specified in the Tune window.

Use Collision Energy Ramp

For MS/MS scans, enter the start and end energy values.

Use Charge State Recognition

Rule: This option is available only if ETD fragmentation mode is set on the Acquisition tab.

The Use Charge State Recognition option sets the collision energy according the charge state z and m/z of the precursor, as determined by the software. Collision energy files for charge states 1 through 6 can be created and selected. The files specify a list of m/z values and associated collision energies. The applied energy is interpolated linearly between the specified masses within the file, so different energies are associated with each m/z across the mass range.

Example: When a mass of interest is detected in an MS experiment, its charge state is calculated. Using the mass and its charge state, you can obtain a collision energy value from the charge state specified in the table. The value is then used in the MS/MS experiment.

CS files are text files with a list of CS, comma-separated, mass and collision-energy value pairs. To edit the files in the CE Control tab, click Modify, and select the appropriate tab. When you finish, click Save and OK.

Cone Voltage tab

Cone Voltage parameters:

Parameter	Description
Override Cone Voltage Specified in tune file	Enter the preferred cone voltage to override that specified in the .ipr tune file.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. The cone voltage follows a linear voltage ramp between these two values during the scan.

MS/MS Cone Voltage parameters:

Parameter	Description
Override Cone Voltage Specified in tune file	Enter the preferred cone voltage to override that specified in the .ipr tune file.
Ramp the Cone Voltage during the scan	Select this option, and enter the start and end values. The cone voltage follows a linear voltage ramp between the two values during the MS/MS scan.

[Back to top](#)

Mobility

Data Collection

Select the Add Drift Time Function to save collapsed retention time data to a separate function in the Raw file.

MS Survey:

Parameter	Description
-----------	-------------

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave settings will be used.
Apply Rule File for Charge State / Drift Time Stripping	Select this option to filter the acquired data using a rule file (.rul).

MS/MS:

Parameter	Description
Maintain Mobility Separation	Select this option to maintain mobility separation. If enabled, Tune Page Transfer T-Wave settings will be used.
Fragmentation	<p>Use Values Specified on the Trap CE and Transfer CE Control Pages:</p> <ul style="list-style-type: none"> Select to use Collision Energy settings as specified on the Trap CE and Transfer CE Control tabs. <p>Fragment before Ion Mobility Separation:</p> <ul style="list-style-type: none"> Select Trap CE Control to fragment the precursor ion before mobility separation. Settings on the Trap CE Control tab and Tune page settings for the Transfer CE are be used. <p>Fragment after Ion Mobility Separation:</p> <ul style="list-style-type: none"> Select Transfer CE Control to fragment the resulting ions after mobility separation. Leave it de-selected to use Tune page settings. This option is selected automatically if you have the Use Collision Energy From Lookup Tables option selected on the Transfer CE Control tab. Selecting Use Precursor Drift Times Only will automatically apply a rule to MSMS data that filters out drift times at which the precursor ion is not present.
Sensitivity	<p>Normal:</p> <ul style="list-style-type: none"> Select Normal for normal sensitivity. <p>EDC:</p> <ul style="list-style-type: none"> Select EDC (Enhanced Duty Cycle) to boost the sensitivity for the selected masses. For this option you have only 5000.0 or 100000.0 end mass selection options. <p>HDC:</p> <ul style="list-style-type: none"> Select to use settings in the HDC lookup file on the MSMS data.

[Back to top](#)

Related Topics

- [Adding lock mass correction to methods](#)



Product ion Discovery (PID) method

- [MassLynx > MS Method > PID Product | PID Neutral](#)

The Product Ion Discovery method provides the opportunity to target compounds of interest for MS/MS analysis based on 'exact neutral loss' (detection of specific 'neutral loss' events with exact mass accuracy) or 'product ions'. The diagnostic neutral loss or product ion events can be specified by the user in the Product Ion or Neutral Loss tabs

respectively, to help drive MS/MS acquisitions for a particular class of molecules or modification (e.g. hydroxylated small molecules, glycosylated peptides etc).

To create an PID experiment:

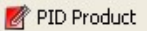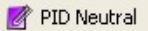1. From the main MassLynx window, click MS Method.
2. Click PID Product  or PID Neutral .
3. To avoid acquiring data when the product ion or neutral-loss is not found, on the MS/MS tab select the option to Stop the MS/MS acquisition if the expected product ion or neutral loss is not found.
4. Use the instructions for setting up MS Survey experiments with the additional selection criteria on the tabs described below.
 - Creating data directed MS and MS/MS survey experiments (see Related Topics)
 - PID Product (below)
 - PID Neutral Loss (below)

PID Product tab

Define your list of product ions on the Product Ion tab according to the descriptions in the following table.

Product ions criteria parameters:

Parameter	Description
High Energy (volts)	Specifies the collision energy used during the high-energy survey scan.
Low Energy (volts)	Specifies the collision energy used during the low-energy survey scan.
Retention Time Window +/- (seconds)	Specifies a tolerance to apply to the retention times in the product ion list.
Mass Tolerance Window +/- (mDa)	Specifies the mass tolerance to be placed around each mass on the product ion list.
Peak Threshold	The minimum peak intensity of the product ion to consider the detection a positive identification.
Ramp High Collision Energy	Ramps the collision energy in the high-energy survey scan between two specified values within each scan.

Product ions file list:

Parameter	Description
All product ion entries must be present in a single MS Survey scan for a valid mass selection	Specifies that all masses in the product list must appear in a single spectrum for a switch to MS/MS to occur.
Masses from File	Enter the file name.

PID Neutral Loss tab

Using the Neutral Loss tab, you can configure a list of neutral loss masses to search for during the parent switch experiment. The Neutral Loss tab parameters are described in the following tables.

Neutral loss criteria:

Parameter	Description
-----------	-------------

Parameter	Description
Use Neutral Loss peak selection	When selected, allows product ions to be used during the acquisition.
High Energy	Specifies the collision energy used during the high-energy survey scan.
Low Energy	Specifies the collision energy used during the low-energy survey scan.
Mass Tolerance Window +/- (mDa)	Specifies a tolerance to apply to the retention times in the Neutral Loss list.
Peak Threshold	The minimum peak intensity in counts/second for which the neutral loss is considered found.
Ramp High Collision Energy	Ramps the collision energy in the high-energy survey scan between two specified values.

Neutral Loss File list:

Parameter	Description
All Neutral Loss entries must be present in a single MS Survey scan for a valid mass selection	Specifies that all masses in the Neutral Loss list must appear in a single spectrum for a switch to MS/MS to occur.
Masses from File	Specifies the Include Mass file name to use during the acquisition.

To edit the Neutral Loss list

1. On the Neutral Loss tab, click Add to open the Neutral Loss dialog box.

Tip: To edit an existing entry in the Exclude list, double-click the entry.

2. Enter the values required.

Related Topics

- [Fast-DDA method](#)
- [Prepare the instrument for use](#)

Time Aligned Parallel (TAP) fragmentation method

- [MassLynx > MS Method > TAP](#)

Tap fragmentation provides a distinct advantage for structural elucidation of small molecules compared to traditional MS or MS/MS techniques.


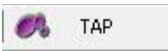
The TAP fragmentation data acquisition methodology enables first and second generation product ions to be acquired for up to 15 target compounds from a single injection, by specifying target m/z and rt time values for each component of interest. The method includes LockMass to enable exact mass measurement, allowing elemental composition analysis and fragment spectra analysis with MassFragment.

Multiple components of interest can be individually selected for TAP fragmentation in a UPLC gradient and then subjected to two stages of CID which provides extensive fragmentation with high resolution and exact mass to aid unambiguous structural elucidation.

Using an associated application, MS^E data viewer software, Ion mobility separation plays a key enabling role, separating first generation fragments and second generation fragments and then, through 'drift time' values, aiding the automated association of fragments.

Note: Data Viewer is a separate MassLynx application accessed as a Windows PC tool, it is not part of the MassLynx application.

To define target compounds and set collision energies:

1. In the MassLynx home page, click  .
2. Click  .
3. In the Function Mobility dialog, set Acquisition and TOF MS/MS requirements and Set Mass (fixed) to m/z of interest.
4. Select the Mobility tab and make sure the Maintain Mobility Separation box is checked.
5. Set CE values as table.
Note that the CE provided are an example only and may need to be appropriately optimized for ions of interest.

Example:

Parameter	Function 1 Value	Function 2 Value
Trap CE ramp	15 - 35 V	Off
Transfer CE ramp	15 - 35 V	15 - 45 V

6. Click OK.

To start MS^E Data Viewer:

At the workstation Windows desktop, click Start > All Programs > MS^E Data Viewer.

Result: The application opens with a set of empty data grids.

Setting parameters for TAP data in MS^E Data Viewer:

1. Click Tools > Peak detection parameters and answer Yes to override settings query.
2. In the Editor dialog, under Peak Width injection, make sure auto, boxes are selected.
3. If data is to be acquired by LockMass, enable LockMass and specify accurate mass of LockMass compound in charge state 2 box.
4. Set filtering as shown in table



Filter	Counts
Low energy intensity threshold	250.0
High energy intensity threshold	100.0

5. Click Update Parameters.

Tip: An application BatchApex PLG is supplied, which enables you to build up a batch of data and thus speed up loading into MS^E Data Viewer.

To open and investigate TAP data files in MS^E Data Viewer:

1. In Data Viewer, click File > Open and browse to available TAP raw files then select RT + DT.
2. Select the box you want de-isotoped. and click OK.
3. Click the HDMS low energy button to set DT **off** for middle (1st generation fragments) spectrum.

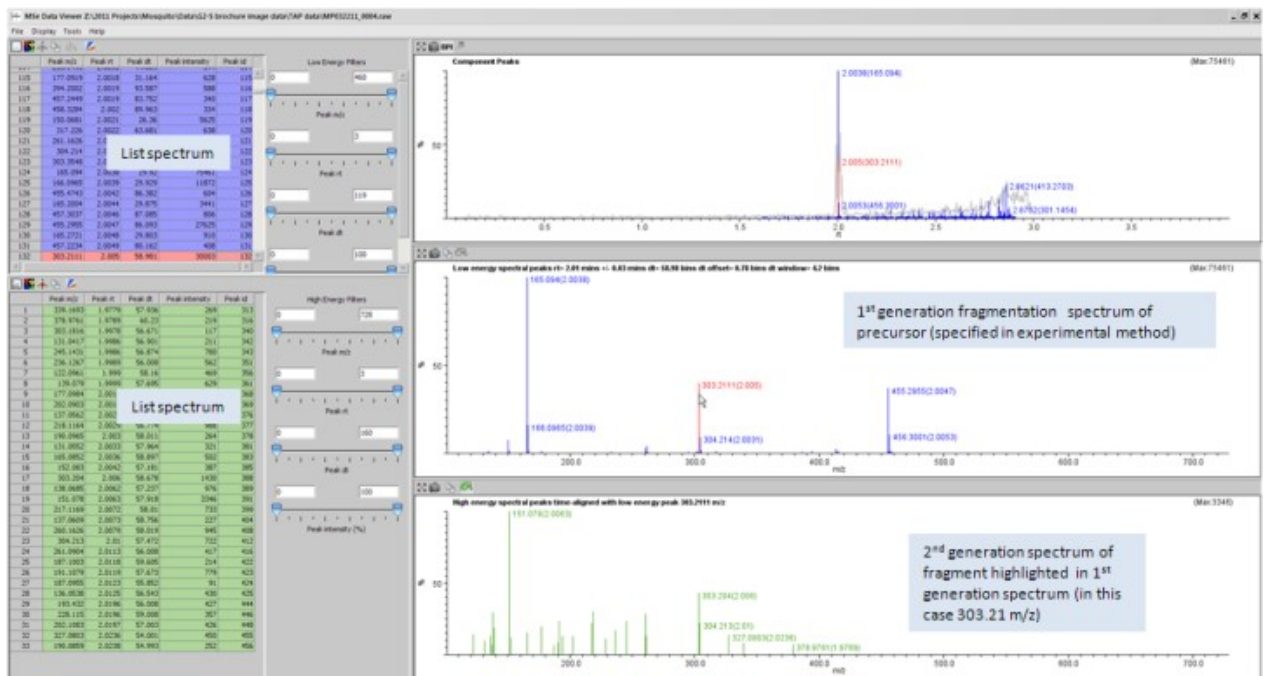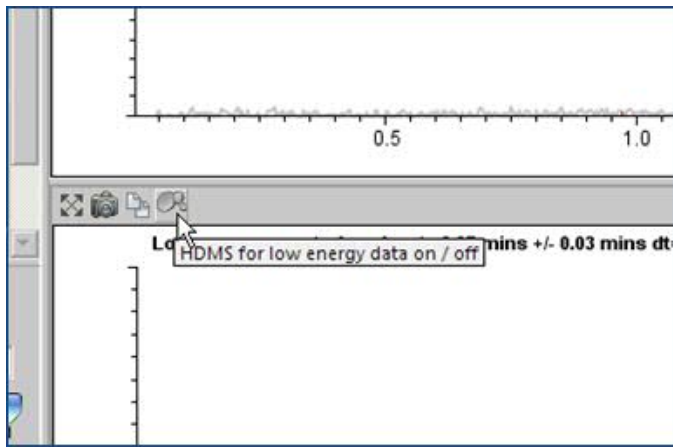

Result: Processing commences and upon completion, data is displayed in each of the panels.

4. Select the chromatic peak of interest by either clicking its mass from the low energy table in the top left panel, or by clicking the peak in the chromatogram display.

Result: Peak turns red

Result: 1st generation fragment appears in middle display panel.

5. Select a 1st generation fragment of interest in the middle panel.

Result: corresponding 2nd generation fragments are displayed in bottom panel for selected 1st generation fragment.

To exit MS^E Data Viewer:

Click File > Exit.

Related Topics

- Creating Mobility MS and MS/MS methods

Create Mobility mode lookup tables

- MassLynx > MS Tune > Tune window > Setup > Look-up Tables

Use lookup tables to finely control the Triwave ramping of parameters. To do so, create tables in the table editor:

1. Transfer CE – Applies ramps to the transfer collision energy. Reference from a function created in the method editor.
2. IMS wave height – Applies a ramp to the wave height.
3. IMS wave velocity – Applies a ramp to the wave velocity.

By default, the gradient created by the table editor is divided equally among the 200 slots of the drift time.

Automatically generating tables

To create a linear lookup table:

1. Select Linear Population.
2. Enter a start and end value.
3. Select or clear the Full Table check box.

Tip: If you select full table, the gradient is created across the full range. If not, you can apply, different gradients across different parts of the range. See "To create several gradients across a range".

4. Click Populate.

Result: The table fills with data, and you can view the created ramp in the graph.

To create several gradients across a range:

1. Clear the Full Table check box.

Result: Start item and End item parameters are enabled.


2. Enter a start value and an end value in the Start Item and End Item boxes.
3. Click Populate.

Result: The table is populated up to the value entered in the End item box.

4. Enter in the Start item field a value 1 greater than was initially entered in the End item box.
5. Click Populate.

Result: The table fills with data up to the value entered in the End item box.

6. Continue with the procedure until the table fills.

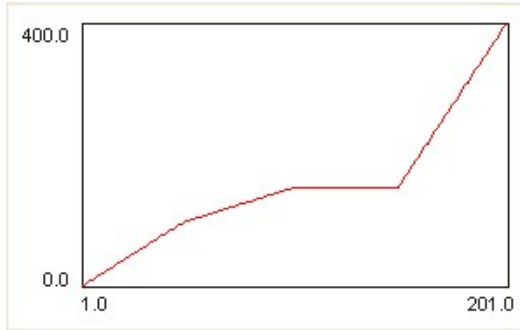Table editor graph with four gradients applied:**Manually creating and editing tables**

You can create or edit tables manually by using the add, remove, and insert functions.

To add a value:

Enter an item value, and then click Add.

Result: The value appends to the end of the table.

To insert a value:

1. Highlight a position in the table.
2. Enter an item value, and then click Insert.

Result: The value is inserted above the highlighted position.

To remove a value:

Select the value in the table, and then click Remove.